



INSTRUCTIONS MANUAL

TYL-TB5011

CORDLESS WATER FLOSSER

For easy, effective flossing



www.tylrhome.com

Table of Contents

IMPORTANT SAFETY INSTRUCTIONS.....	2
KNOW YOUR CORDLESS WATER FLOSSER.....	5
CHARGING YOUR CORDLESS WATER FLOSSER	6
USING YOUR CORDLESS WATER FLOSSER.....	7
Filling the Water Reservoir	7
Inserting and Removing Nozzles	7
Using the Water Flosser.....	7
CLEANING YOUR CORDLESS WATER FLOSSER	9
TROUBLESHOOTING	10
WARRANTY TERMS AND CONDITIONS	11

IMPORTANT SAFETY INSTRUCTIONS

When using electrical appliances, basic safety precautions should be followed, including the following:

DANGER! To reduce the risk of electrocution:

- Do not handle the charger with wet hands.
- Keep the charger away from water. Do not place or store it over or near water contained in a bathtub, wash basin or sink.
- Do not immerse the charger in water or any other liquid.

WARNING: To reduce the risk of burns, electrocution, fire, or injury to persons:

- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- Do not charge this unit with any charging devices other than the one provided.
- Do not operate this unit with a damaged charger, cord or plug or after the unit malfunctions, or has been damaged in any manner or dropped into water. Return unit to an authorized service center for examination and repair.
- If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons


in order to avoid a hazard.

- Keep the cord away from heated surfaces.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Avoid applying excessive nozzle tip pressure on the teeth and gums.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Do not drop or insert any foreign object into any opening on the product.
- Do not use the appliance with mouthwash which may cause corrosion.
- Do not fill the water tank with water exceeding 40°C to avoid scalding.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Consult your doctor if you have had any

dental surgical operation in past two months or are under dental treatment.

**SAVE THESE INSTRUCTIONS. FOR HOUSEHOLD USE ONLY.
NOT FOR COMMERCIAL USE.**

KNOW YOUR CORDLESS WATER FLOSSER

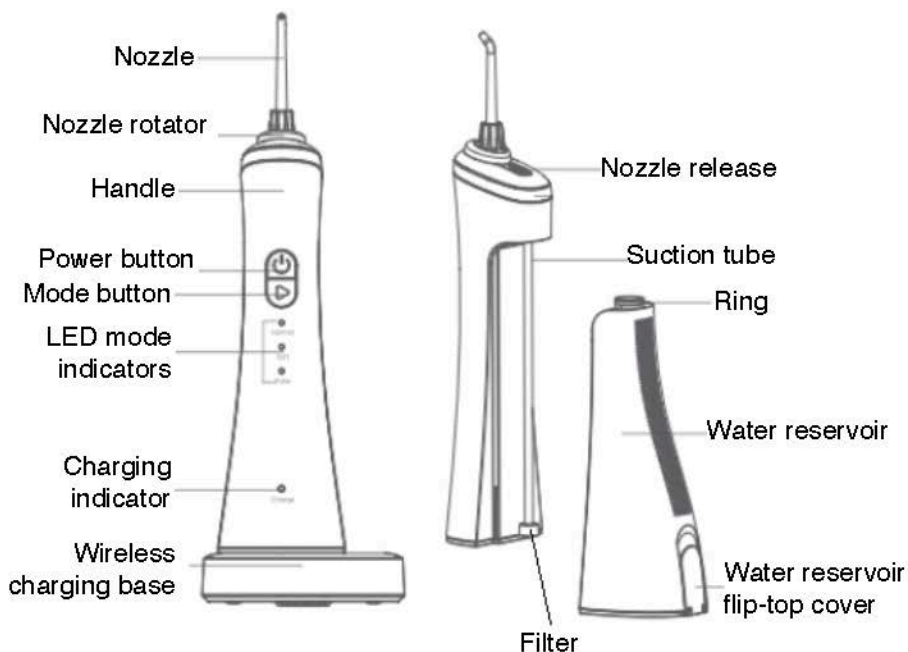

CHARGING YOUR CORDLESS WATER FLOSSER


Charge for 12 hours before first use.

1. Place the charging base on a level, non-skid surface, then plug into a power outlet.



2. Place the handle on the charging base. The charging indicator on the handle will slowly blink when charging and will stay steady when the unit is fully charged. A full charge will last 30 minutes of use.

When the Water Flosser needs charging, the charging indicator light will blink before and after use to indicate that you should recharge the unit.

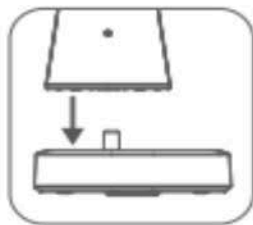

USING YOUR CORDLESS WATER FLOSSER

Filling the Water Reservoir

Lift water reservoir flip-top cover and fill reservoir with water. Close the cover.

For ease of adding water, you may remove the reservoir and add water from the opening the top opening of the reservoir. The water reservoir may be removed from the handle by sliding it towards the base of the handle. Turn on the unit and let it run to drain the residual water from nozzle.

Inserting and Removing Nozzles

Insert nozzle into the ring on the handle. Press nozzle firmly until the nozzle clicks into place. To remove the nozzle from the handle, press the nozzle release button and pull nozzle from the handle.

IMPORTANT: *Do not press the nozzle release button while unit is running.*



Using the Water Flosser

1. Press the Mode button to select preferred mode: **Normal**, **Soft** and **Pulse**.

NOTE: *Start with **Soft** setting when using the Water Flosser for the*

first time, then switch to other modes once you have adapted to its use.

2. Lean over the sink and place nozzle in the mouth. Aim nozzle towards the gumline and press the Power button turn the unit on. Start with the back teeth and work towards the front. Glide nozzle along the gumline and pause briefly between teeth, keeping the unit upright as much as possible. Direct the jet stream at a 90°

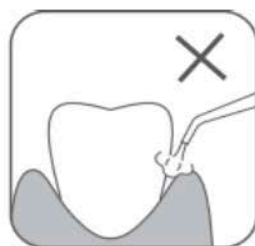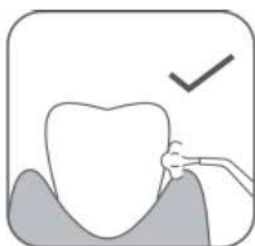

angle. Slight close lips to avoid splashing but allow water to flow freely from the mouth to the sink.

3. Press Power button to turn the unit off.

NOTE: *The unit will automatically turn off after running for 2 minutes. Press the Power button again to turn the unit on.*

4. Remove residual water by emptying the water reservoir, and by turning the unit on and letting it run for a few seconds to drain the residual water from nozzle.
5. Wipe the unit with a dry cloth or towel.

CLEANING YOUR CORDLESS WATER FLOSSER

Never immerse unit in water. Before cleaning, unplug from the electrical outlet.

- Clean the handle by using a soft cloth and mild non-abrasive cleanser.
- Detach nozzle from handle and rinse with water. Wipe with a soft cloth.
- Rinse water reservoir with water. Wipe with a soft cloth. Empty the water reservoir completely if unit will not be used for an extended period. Turn the unit and let run until no water comes out of the nozzle.
- Unplug charger from outlet. Clean charger with a dry cloth.

TROUBLESHOOTING

Problem	Possible Cause	Solution
The unit does not work.	<ul style="list-style-type: none">• The unit is used for the first time or has not been used for over 3 months.	<ul style="list-style-type: none">• Charge the unit for 24 hours before use.
The water pressure from the nozzle is low.	<ul style="list-style-type: none">• There is a blockage in the nozzle.• The nozzle is deformed.• The filter is blocked.	<ul style="list-style-type: none">• Replace the nozzle.• Clean the filter.
The unit runs for a short time even when fully charged.	<ul style="list-style-type: none">• Battery is at the end of life.	<ul style="list-style-type: none">• Contact your authorized service center.
The charging indicator does come on while charging.	<ul style="list-style-type: none">• Unit is not properly plug to the outlet.• Circuit has malfunctioned.	<ul style="list-style-type: none">• Check that the plug is properly connected to the outlet.• Contact your authorized service center.

WARRANTY TERMS AND CONDITIONS

ONE-YEAR LIMITED WARRANTY ON PARTS AND LABOR

TYLR provides a warranty to the original purchaser of this product against defects in materials and workmanship for one year. If a product covered by this warranty is determined to be defective within the warranty period, **TYLR** will repair or replace the defective part with a new or remanufactured part free of charge subject to terms and conditions stated herein.

To obtain warranty service, proof of purchase in the form of a sales invoice or copy thereof is required to show that a product is within the warranty period. If any one of the documents cannot be produced, corresponding fees for labor and replacement parts will be charged.

Parts and Labor

There will be no charge for parts and labor during the warranty period. Replacement parts or units may be new or recertified and are upon **TYLR's** option and sole discretion. Replacement parts and units are warranted for the remaining portion of the original warranty. Our decision on all questions relating to complaints as a result of defects, either workmanship or materials, shall be conclusive and you shall agree to abide by such decision. Any unit or defective part that have been replaced shall become **TYLR's** property. After the lapse of this warranty, a charge will be made on all labor and replacement of parts.

Types of Services

Defective Products must be sent to any **TYLR** service center to obtain warranty service. **TYLR** is neither responsible for transportation costs to the service center nor **TYLR** will cover return shipping to the customer. Product returns to **TYLR** service centers must utilize either the original carton box and shipping material or packaging that affords an equal degree of protection.

For home service, the transportation fee of **TYLR** service personnel will be charged to the customer, depending on the location. Our service personnel reserve all rights to refuse to attend to any equipment installed in an inaccessible location or any other hazardous situation.

LIMITATIONS AND EXCLUSIONS

TYLR one-year limited warranty only covers defects in materials and workmanship; however, this warranty does not cover the following:

1. Damage resulting from accidents, transport, misuse, abuse, alteration, tampering, or failure of the purchaser to follow normal operating procedures outlined in the user's manual.
2. Damage caused by negligence and poor maintenance.
3. Defects or damage due to spillage of food/liquids or improper usage of electrical supply or electrical circuit, major fluctuations in voltage, exposure to rain or moisture, wrong fuel or oil used.
4. Damage, losses, and defects as a result of fire, flood, or other Acts of God.
5. Normal wear-and-tear, corrosion, rusting, or stains.

6. Defects & damage arising from improper testing, operation, usage of the wrong component, demonstration, maintenance, installation, adjustment, or any alteration or modification of any kind.
7. Scratches & damage to the outer surface areas and externally exposed parts that are due to normal customer use.
8. General maintenance & routine servicing.
9. Claims for damaged/missing parts after 7 days from the date of purchase.
10. If any part or parts of the unit are replaced with a part or parts not supplied or approved by us, or if the unit has been dismantled or repaired by any person other than a **TYLR** authorized technician.
11. Any equipment/product which has its serial number removed or made illegible/tampered with.
12. If the warranty card is altered, defaced, or erased in any manner whatsoever.
13. If the unit is used for purposes other than household use, such as commercial use.

The six (6) months warranty period covers the rechargeable batteries and the one (1) year warranty period covers the mainboard and control PCB, motor and power supply. However, this warranty does not cover parts (i.e. accessories, charging base, bristles, nozzles, fuse, gears, interior and exterior casing) that are subject to wear and tear, rust, stains, or corrosion.

This one-year limited warranty does not cover Products sold **“AS IS”**, **“FACTORY RECTIFIED”**, or by a non-authorized reseller.

THERE ARE NO EXPRESS WARRANTIES OTHER THAN THOSE LISTED OR DESCRIBED ABOVE. ANY IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE PERIOD OF TIME SET FORTH ABOVE. **TYLR'S** TOTAL LIABILITY FOR ANY AND ALL LOSSES AND DAMAGES RESULTING FROM ANY CAUSE WHATSOEVER INCLUDING **TYLR'S** NEGLIGENCE, ALLEGED DAMAGE, OR DEFECTIVE GOODS, WHETHER SUCH DEFECTS ARE DISCOVERABLE OR LATENT, SHALL IN NO EVENT EXCEED THE PURCHASE PRICE OF THE PRODUCT. **TYLR** SHALL NOT BE RESPONSIBLE FOR LOSS OF USE, COMMERCIAL LOSS, LOST REVENUE OR LOST PROFITS, OR OTHER INCIDENTAL OR CONSEQUENTIAL DAMAGES. THIS WARRANTY IS SUBJECT TO CHANGE WITHOUT NOTICE.

IMPORTANT: Please present sales invoice as proof of purchase whenever you require our service during the warranty period.

www.tylrhome.com