

PROLON

DENTAL STERILISER



MODEL: PDS001

PROTECT YOUR DENTAL HEALTH

Thank you for purchasing our dental steriliser PDS001. Should you have any questions or concerns about using your new product, feel free to reach out to our helpful customer support team. We hope you enjoy your new product.

Customer Service Contact

NZ Customers: 0800100345
International: +6463501617

IMPORTANT


Please read and comply with all of the instructions and warnings provided in this manual before using the product. Failure to comply with the instructions and warnings provided here may result in inaccurate results and/or damage to the product itself.

WARNING

DO NOT stare at the UV light or expose skin directly to UV light.
DO NOT disassemble or repair the product.
DO NOT operate or place the product on or near hot surfaces such as ovens, heaters, or fireplaces.
Keep this product out of reach from children.
Keep the product away from water or other liquid.
Store this product away from direct sunlight or chemicals.
For the safest and fastest charge, use original or certified cords.
Avoid dropping.
Essential oils should not come into contact with denture appliance.

1

PRODUCT DIAGRAM

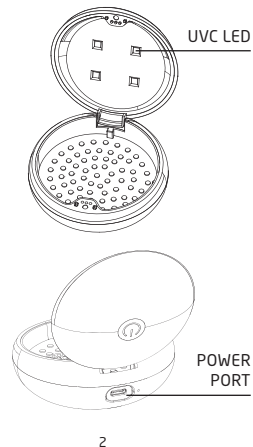

SPECIFICATIONS


NAME:	PROLON UV DENTAL STERILISER
MODEL:	PDS001
DIMENSION:	Ø84.7*44MM
RATED VOLTAGE:	5V
RATED POWER:	5W
BATTERY CAPACITY:	500MAH
CHARGING TIME:	2H




Innova Products Limited
322 Neilson Street, Onehunga,
Auckland, New Zealand.
Phone: 0800 100 345
Email: nzinfo@innovaproducts ltd.co.nz
www.prolon.co.nz

STERILISING METHOD

1 Open the lid.

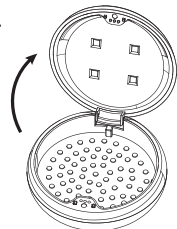

2 Place mouthguard or dentures in the machine.

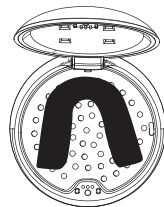


4

3 Press the button to start and light will illuminate.



4 The UV lamp will start to work for 3min, and the top indicator light will be on when it is working. Light will go off at end of UV cycle.

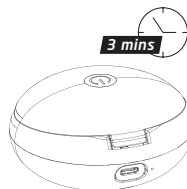


5

5 After cycle is complete, remove mouthguard or dentures and wear with confidence.



6 The indicator light will flash when the battery level is lower than 10%, indicating charging is required. At this time, plug in the power charging cord to charge.

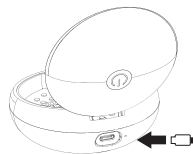

6

TERMS AND POLICY

We warrant all products to be of the highest quality in material, craftsmanship and service for 2 years, effective from the date of purchase.

Under the limited warranty, we will replace or refund any product found to be defective due to manufacturer.

It is recommended you retain your receipt of purchase for this product and use the product as per the instruction booklet provided.

Any warranty claims in the first instance must be to the place of purchase but should this not resolve your issue please contact our customer service team on our free phone within New Zealand or email.

Refer to back of booklet for customer contact details.

7