



Sonic / 3000

Operating Instructions

WELCOME TO THE KYOUI FAMILY!
Thank you for choosing us!

Our goal is to improve people's overall health by providing more effective oral care tools.

ATTENTION

Don't forget to charge your toothbrush for 24h before using it for the first time!

BATTERY TYPE:

Ni-MH Battery

CHARGER INFORMATION:

International Voltage (100 - 240V)



Kyoui USA
info@kyoui.com

Before operating this product, please read these

STEP- BY- STEP GUIDE FOR BETTER ORAL HYGIENE

1. **ONCE A WEEK** **PICK HEAD - PLAQUE/TARTAR REMOVER**
Use the PICK HEAD to clean your gum line and the back of your teeth.
2. **MORNING TIME OR TWICE A DAY** **CLEANING HEAD**
Use the BIG HEAD (wavy bristles) to clean your teeth in the morning or twice a day.
3. **NIGHT TIME** **PERIO/POLISHING HEAD**
Use the SMALL HEAD if you have sensitive gums and to polish your teeth once a day .



* Need replacement heads? Please visit our website: www.Kyoui.com

Dental Professional Endorsements:

“The Sonic 3000 is the most versatile and innovative toothbrush I have ever used”
Sondra Schachtmeister, RDH, USA

“The Angle Neck makes all the difference when it comes to hard to reach areas, it’s more effective” Jean Pierre Squaglia, DDS, France

“After using the Kyoui brush any other brush doesn’t seem to clean well enough”
Timothy Bradbury DDS, USA

“The Kyoui Sonic 3000 is the most effective tool for oral care available in the market today.” Guy Ananou DDS, France

“The Angle Neck and Dupont Bristles are the best combination for an outstanding cleaning and safe gums.” Cyrille Tremblay DDS, Canada

Important

If you experience discomfort or pain in your mouth, discontinue use and contact a dentist.

- If you have had oral surgery within the past 2 months, please consult your dentist before using this product.

The toothbrush has been tested and is compliant with safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer prior to use.

Do not operate this product if: It has a damaged cord or plug, it is not working properly, it has been dropped or damaged, or the charger has been submerged in water.

Keep the cord away from heated surfaces.

Do not block the air openings of the product or place on a soft surface, such as a bed or couch which may potentially block air openings. Keep air openings free of lint, hair, and any other foreign debris.

Do not use product while sleepy or drowsy.

Use of product contrary to what has been described in this manual is expressly prohibited.

Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

Configuration

a. Patented Angle Neck PERIO brush head

b. Identification ring

c. Toothbrush Handle

d. On/Off button

e. Interdental brush head (pick)

f. Patented Angled Neck
CLEANING brush head #02

g. Charger


h. Charging indicator ring

i. Mode button


j. Mode indicator

k. Charger cover

l. Power cord

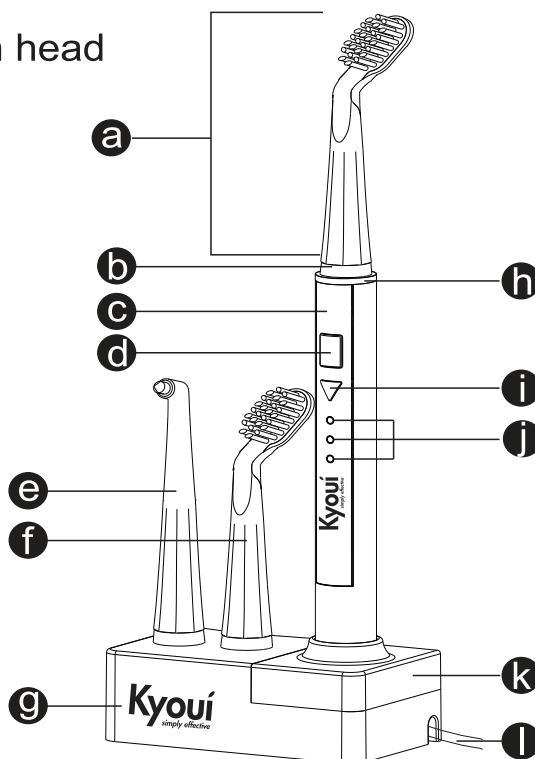


Getting Started

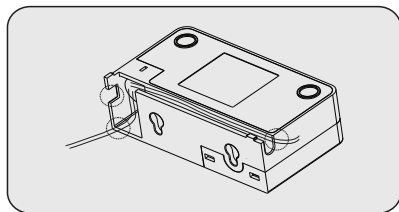


Installing the charger Power Cord

Place cord in the right, left or center cord guide, depending on where you'll place it.

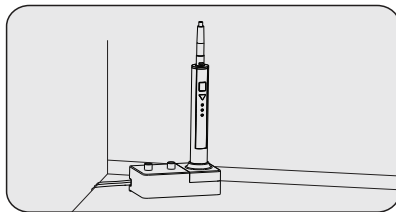


Installing near a sink or on a shelf

Place the charger on a secure surface. Insert the plug into an electrical outlet.

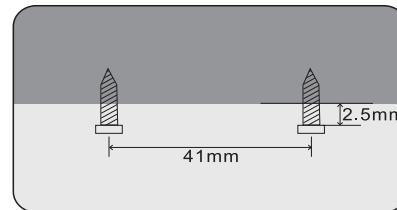


Mounting on the wall

Purchase two round wood screws (diameter: 2.8~3.2mm), fasten the screw on the wall, so that they extend about 2.5mm from the wall and are positioned 41mm apart.

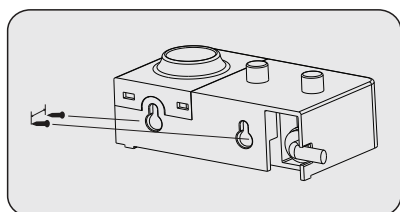

Align the holes on the back side of the charger with screws and secure in place.

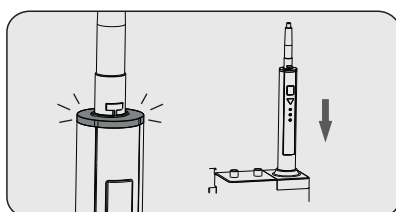



Charging the toothbrush: be sure to turn off the toothbrush before charging. Place handle in the charger. The charging indicator ring will light up red.

Notes:




- ① It takes at least 24 hours to charge the battery fully. Indicator light will change into blue after fully charged.
- ② The toothbrush must be charged for 24 hours, if:
 - a. this is the first time using the toothbrush;
 - b. the toothbrush has not been used for more than three months; or
 - c. if the charging indicator blinks red during use.
- ③ To keep the battery fully charged, we recommend keeping toothbrush in the charger when not in use.
- ④ During charging, the toothbrush will get warm, but this is normal.
- ⑤ The operating duration will become shorter as the battery gets older. If toothbrush operates only 1 short time after being fully charged, this means that the battery is nearing the end of its life cycle.

How To Use

- ① Align the slot of the brush head to the protusion of the handle and slide down and then turn the brush head until it clicks into place.
- ② Wet the bristles and apply toothpaste
- ③ Place the bristles against the teeth at a slight angle (45 degrees) towards the gums line.
- ④ Press the mode button to select mode.
 - a. Normal: Standard mode for superior teeth cleaning.
 - b. Soft: Gently cleans your teeth or any other sensitive areas.
 - c. Pulse: Pulsates for gentle gum stimulation.

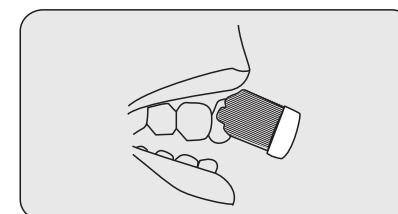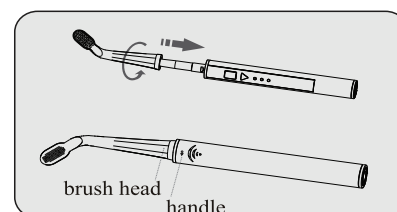

- ⑤ Press the on/off button to turn on the toothbrush
- ⑥ Apply light pressure and gently move the brush head slowly across teeth with bristles in a 45 degree to the gum line
- ⑦ In order to maximize cleaning effectiveness, divide your mouth in 4 sections. Start with Section 1 and brush each section for 30 seconds. The brushing action will pause between each 30 second interval, indicating when to move on to the next section. The toothbrush will automatically stop after 2 minutes of use.
- ⑧ Turn the brush head in direction as shown, and then pull it to detach from the handle.
- ⑨ Rinse both the brush head and handle under running water wipe them dry.
- ⑩ Put the handle back into the charger

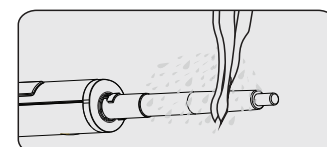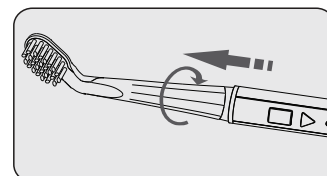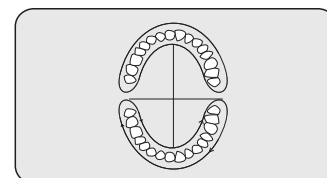


Maintenance


Be sure to remove the plug from the electrical outlet before cleaning the product.

Charger

Pull the charger cover to remove. Using a soft cloth, wipe off any dirt on the charger and the charger cover. Clean the plug periodically to ensure it remains in good working order. The accumulation of dust or dirt on the plug can cause electrical shock, overheat or fire.

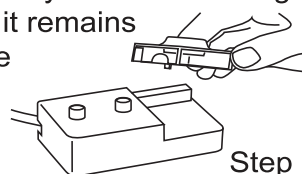

Step 1

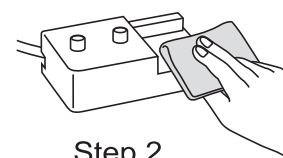

Step 2

Toothbrush

Remove the brush head. Wipe surface areas with a damp cloth and place the handle back into the charger. Rinse the brush head and bristles after each use. Do not clean in the dishwasher.

Kyoui Sonic 3000 System Warranty Disclaimer

Thank you for your interest in the products and services of Kyoui Oral Care.

This Limited Warranty applies to Kyoui Sonic 3000 system (here referred as "Physical Goods"), and only for Kyoui Sonic 3000 system, purchased from Kyoui Oral Care.

This Limited Warranty covers any defects in material or workmanship under normal use during the Warranty Period of 2 years from the date of purchase. During the Warranty Period, will repair or replace, at no charge, products or parts of your Kyoui Sonic 3000 system that proves defective because of improper material or workmanship, under normal use and maintenance.

Our production department will repair the Physical Goods at no charge, using new replacement parts. The Warranty Period for Physical Goods is 2 Years from the date of purchase.

A replacement Physical Good or part assumes the remaining warranty of the original Physical Good or 2 Years from the date of replacement or repair, whichever is longer.

This Limited Warranty does not cover any problem that is caused by:

- Conditions, malfunctions or damage not resulting from defects in material or workmanship
- Conditions, malfunctions or damage resulting from negligence, improper maintenance or modification
- Damaged by natural causes including but not limited to lightning, flood, or other natural disaster.
- Theft or loss of the Physical Goods

What do you have to do?

To obtain warranty service, you must first fill out a Warranty Claim to determine the problem and the most appropriate solution for you.

Deliver the Physical Goods, in either its original packaging or packaging providing an equal degree of protection, to the address specified by Kyoui Oral Care.

You can contact us by any of the following contact methods:

- By email: warranty@kyoui.com
- By visiting this page on our website:

www.Kyoui.com --> Warranty Claim (locate on the top right side of the home page.

Implied Warranties

EXCEPT TO THE EXTENT PROHIBITED BY APPLICABLE LAW, ALL IMPLIED WARRANTIES (INCLUDING WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE) SHALL BE LIMITED IN DURATION TO THE DURATION OF THIS LIMITED WARRANTY.

Limitation of Damages

EXCEPT TO THE EXTENT PROHIBITED BY APPLICABLE LAW, SHALL NOT BE LIABLE FOR ANY INCIDENTAL, INDIRECT, SPECIAL, OR CONSEQUENTIAL DAMAGES RESULTING FROM ANY BREACH OF EXPRESS OR IMPLIED WARRANTY OR CONDITION, OR UNDER ANY OTHER LEGAL THEORY.

Troubleshooting

| Problem | Possible cause | Solution |
|--|--|---|
| The toothbrush doesn't work. | You have just purchased the product or it has not been used for 3+ months. | Charge the main unit for at least 24 hours. |
| Weak or diminished vibration | The toothbrush is pressed hard against your teeth. | Place the brush head gently against your teeth. |
| The toothbrush operates only for a short period of time even immediately after being charged | The battery has reached the end of its life cycle. | Discard battery (see instructions). |
| | Charging time was insufficient. | Charge the handle for at least 24h. |

Kyoui Oral Care

6787 W Tropicana Ave ste 274
Las Vegas, NV, 89103

Kyoui Customer Care:
info@kyoui.com
www.Kyoui.com

For Warranty Claim, please visit:

www.Kyoui.com

Fill out the Warranty Claim on
the top right side of the
home page.