


## In-Use Screen

After a Massage program is selected, In-use screen displays all settings to that program.

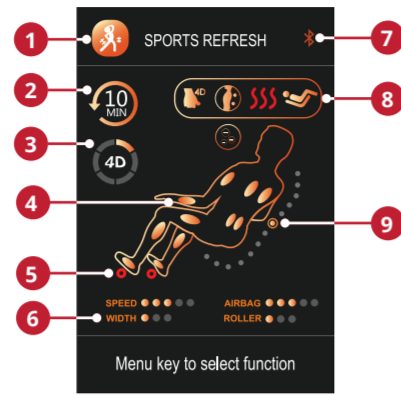

- |                  |                                    |
|------------------|------------------------------------|
| 1 Program        | 6 Indicators                       |
| 2 Time left      | 7 Bluetooth connected              |
| 3 4D intensity   | 8 Status Icons                     |
| 4 Active airbags | 9 Back-roller position & direction |
| 5 Active rollers |                                    |

## Chromotherapy

From the **HOME** screen, select the **LIGHTS** menu. Navigate to **THEMES** and select your desired light setting.

**Solid:** Single color light.

**Breathe:** Single color that dims and brightens on a loop.

**Wake Up:** Blue enriched lights to boost energy levels - gradually increasing brightness.

**Wind Down:** Colors of sunset - gradually decreasing brightness.

**Multi-wave:** All the colors of the rainbow morphing into each other.

**Multi-rotate:** Colors rotate from one to another.

## Intelligent Voice Commands

Awaken Intelligent Voice Commands by saying: **'Switch to Voice Mode'**.

You may now give your chair any of the commands below **within 5 seconds**. Make sure to speak clearly and slowly. After 5 seconds the Intelligent Voice Command function will turn off automatically and must be awakened again.

### Commands

Power On	Deep Shiatsu
Switch Off	Healthy Breath
Activate Massage	Massage Extend
Sports Refresh	Relaxation
Extension	Air Bag Massage
Rest & Sleep	Zero Gravity
Working Relief	Start Heat
Neck & Shoulder	Heat Off
Waist & Spine	

## Sounds

From the **HOME** screen, select the **SOUNDS** menu. Toggle the following settings to play Nature sounds from your chair.

**Sound:** Select from Nature sounds.

**Volume:** Set volume level.

**BEEP:** Set to OFF to stop sounds from remote control.

To play Nature sounds, Bluetooth must be **ON** prior to selection.

## Quick Access Controller

### Button Operation

**On/Off:**


Press and hold.

### Change Auto Program:

Short press; new program will be announced and start.

### Dial Operation


Turn left or right to recline or incline.

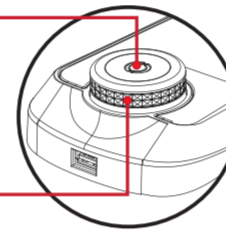

**KYOTA** | The World's Most Advanced Massage Chairs®

# KYOTA®

## Quick Start Guide

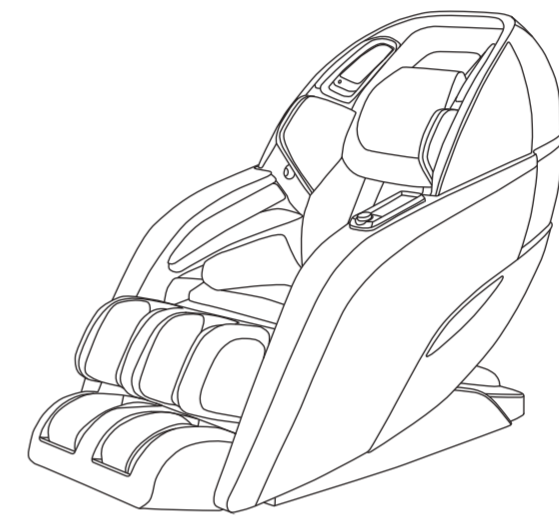

## Yutaka M898 4D Massage Chair

# KYOTA®

**Power** – Powers up the remote or returns the chair to the exit position and powers down remote.

**Home** – Shortcut to the home page of all menu screens.

**OK** – After navigating to the command desired on the display, push OK to select.

**Zero Gravity** – Press for zero gravity 1, press again for zero gravity 2, press again to restore recline.

**Foot Roller** – Press to turn off foot roller, press again to resume foot roller.

**Footrest** – Press and hold up arrow to bring footrest up, press and hold down arrow to bring footrest down.

**Leg Length** – Press and hold right arrow to extend the leg length, press and hold left arrow to retract the leg length.



**Multi-Functional Screen** – Use the directional buttons below to navigate these menu screens. Then press OK to select.

**Pause** – Press pause to immediately stop all massage chair activity, press again to resume.

**Auto** – Shortcut to the auto programs menu screen.

**Directional Button** – Press up, down, left or right to navigate the display screen.

**Back Adjustment** – Back mechanism position adjustment; press and hold up to raise, press and hold down to lower. Use during initial body scan fine tuning or in manual point or partial zone mode.

**4D** – Press to increase or decrease 4D intensity of back mechanism. In-use screen displays intensity dial.

**Seat Back Position** – Press and hold down arrow to recline the seat back, press and hold up arrow to upright the seat back.