


EXTRA WIDE TALKING SCALE

In order to get the utmost accuracy and lifespan of our product, please read the following instruction before use.

Weight memory records (from 01 to 30)

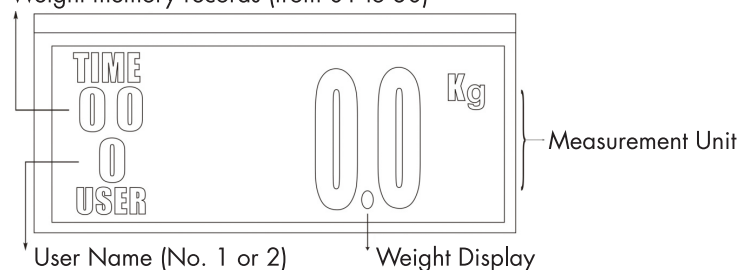

Function and characteristics

- Voice function (English or French)
- Memory for 2 users. (No. 1 or 2)
- Personal data memory: 30 records for each user
- Auto ON

Definitions of the keys

MODE:

- Turns the scale ON
- Displays user's most recent measured weight
- "TIME" record goes back to "00" when viewing recorded weights
- ▲ Select user's number (No. 1 or 2)
- Checks the user's memory records by pressing once, the historical measurement record on the screen will go forward once (from 01 to 30).

NOTE: 01 is the most recent measured weight

- ▼/U – Changes measuring unit: kg-lb. You can only operate this key when you turn on the scale and the screen displays "0.0".
- By pressing ▼/U once, the historical measurement record on the screen will go backward once.

NOTES

1. When the display is flashing, no function can be operated.
2. Make sure to place the scale on flat ground.
3. Press "MODE" key or step onto the scale to turn it on.

Weighing without memory function

Place the scale on flat ground and press "MODE" to turn on the scale. The scale then will display "0.0". Step on the scale and your weight will be displayed; and then read out loud twice. Stepping onto the platform steadily will also turn on the scale, and your weight will be displayed and then read out loud twice. The scale will turn off automatically after reading your weight out loud.

Weighing with memory function

Turn on the scale and it will display "0.0". Press the ▲ key and the user number will begin to flash. Press the ▲ key to select user (No. 1 or 2). Wait 4 seconds, when scale shows "0.0" step onto the platform and your weight will be displayed and read out loud twice. This weighing data will be saved as the user's most recent weight. Each user can save up to 30 recordings.

Changing the language

Press and hold the "MODE" key for few seconds to turn on the scale. "L-01" or "L-02" will be displayed then you can press "MODE" key again to select the language :

- "L-01" for English
- "L-02" for French

Viewing of the historical measurement results

Turn on the scale and it will display "0.0", and then press the ▲ key to select user (No.1 or 2). When the scale shows "0.0", press the "MODE" key which will display the user's latest weight. The memory time will flash; press the ▲ key to view recorded weights from your 1st to your 30th recording. When you are viewing the recorded weights, the scale will not read them out loud.

Deleting memory results

Turn on the scale and it will display "0.0", then press the ▲ key to select user (No.1 or 2). When the scale displays "0.0" press the "MODE" key displaying the users latest weight, then press both "MODE" and "▼/U" for a few seconds. Then all records on that user will be deleted. If you want to delete the other user's records, re-start the scale and repeat the procedure above.

Battery installation

Slide off the battery cover on the back of the scale and connect one 9V battery (not included) onto the proper terminal heads. Slide the battery cover back onto the unit.

Troubleshooting

1. "Lo" : weak battery, replace with a new battery (not included).
2. "OL" : scale overload.
3. "Err": scale incorrect operation.

Technical specification

Max. capacity: 250kg (550 lbs)
Min. capacity: 2kg (4.4 lbs)
Conversion: 0.1kg (0.22 lb)
Weight mode: kg, lb
User memory: 2 users
Operation condition: T: 0-40C RH<85%
Power: 1x9V battery (not included)

NOTES

1. Place scale on flat ground.
2. Keep your scale in a dry place to protect the electronics. Avoid extremely high or low temperature environments.
3. The scale is a high precision electronic product; please avoid excessive impact or vibration to the unit.
4. To prevent scratches on the glass panel, do not wear footwear.
5. Clean with a cloth, do not use solvents or immerse scale in water.
6. The product is designed for household use and should not be used for commercial purposes.
7. Remove the battery if the scale will not be used again for a long period of time.

WARNING: Do NOT use the Memory Storage Talking Scale on carpet, that will reduce the accuracy. Use on dry and flat ground only.

ONE YEAR WARRANTY

If this product proves to be defective in material or workmanship within one year of purchase, please return it to the address below. It will be repaired or replaced without charge upon receipt of the unit prepaid with \$20.00 to cover handling, packaging and return postage. Please include proof of purchase, your full name, address, daytime phone number or email address.

MANUFACTURED BY / FABRIQUÉ PAR:

BIOS | Medical

16975 LESLIE STREET
NEWMARKET, ON L3Y 9A1
MADE IN CHINA / FABRIQUÉ EN CHINE