

Introduction

Thank you for purchasing the BIOS Living Jumbo Digital Display Bath Scale. This scale was developed with state-of-the-art electronics. This precision device has a weight capacity of 550 lbs / 250 kg and is equipped with an extra large blue backlight screen.

Features

- Auto ON / OFF
- Auto zero reset
- Large blue backlit LCD (4" x 3½" / 10.16 x 8.89 cm)
- LB/KG/ST switchable
- Low battery / over-load indication
- Weight capacity: 550 lbs / 250 kg

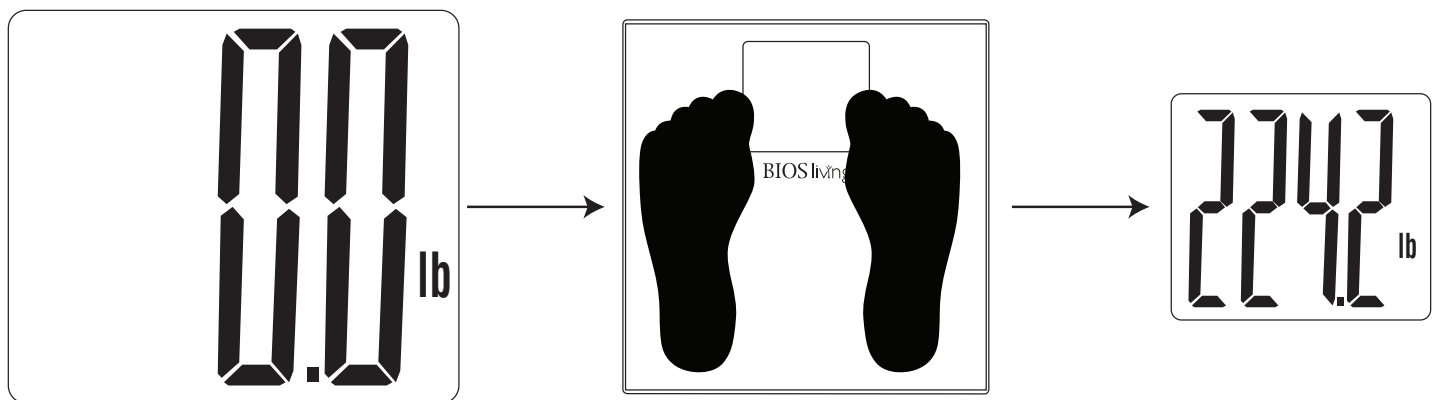
Battery Installation (3 x AAA batteries included)

Open the compartment on the back of the scale and insert the batteries according to their polarity. Please be aware that inserting the batteries incorrectly may cause permanent damage.

NOTE: If you do not intend to use this scale for a long period of time, it is advisable to remove the batteries before storing it away.


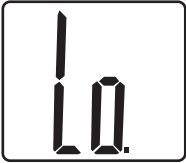
Using Your Scale

- Place scale on a hard flat surface
 - Step onto the scale gently, placing your feet evenly apart, and stand still until the reading blinks 3 times.
- NOTE:** If the scale is placed on a soft surface, such as carpet or an uneven floor, the results may be inaccurate or an error message may occur.



- The scale will shut off automatically after 15 seconds or once weight is removed.

Troubleshooting

Error	Problem	Solution
	An error has occurred when taking a measurement.	Step off the scale and wait a few minutes before trying again. Note: Do not use if you exceed the maximum weight capacity.
	Battery power is running low.	Remove batteries, wait 10 seconds and reinsert them into the scale. If error occurs again, replace the batteries.

Care and Maintenance

- Always keep the scale in a horizontal position and do not put in wet or hot conditions
- Avoid hitting or dropping the scale
- Do not overload the scale, as this will damage the sensor
- Do not immerse this product in water
- Clean the scale with a soft, damp cloth. Do not use chemicals or harsh abrasives
- Keep the scale surface dry and do not use on a slippery floor

Five Year Warranty

If this product proves to be defective in material or workmanship within five years of purchase, please return it to the address below. It will be repaired or replaced without charge upon receipt of the unit prepaid with \$5.00 to cover handling, packaging and return postage. Please include proof of purchase, your full name, address, daytime phone number or email address.

BIOS Medical
16975 Leslie Street
Newmarket, ON L3Y 9A1
www.biosmedical.com