Acupressure Mat

Thank you for purchasing this Bios Living Acupressure Mat Set to aid in your home therapy treatments.

Acupressure treatment has been used for more than a 10,000 years in China. It is said to aid in relaxation, pain relief, the release of toxins, relieve tension in muscles and stimulate the nerve endings to boost circulation and healing.

ATTENTION

***This product is not suitable for everyone. Certain health conditions and their restrictions do not recommend the use of an acupressure therapy mat.**

Please consult a healthcare professional before using in order to determine suitability.

- This product is not recommended for users who suffer from thin skin disorders, diabetes or who are taking blood thinners.
- Do not use this mat if you have a rash, broken skin or are susceptible to infection.
- Do not use for children, Infants and pregnant women

Please note before using your mat set:

The points are sharp, in order to give you the ultimate acupressure experience.

You must use caution when you are lying down.

Lie down slowly allowing your body to adjust to the pressure.

It is recommended to start by using a light shirt or sheet separating your skin from the mat until you become accustomed to the pressure.


The recommended use for the mat is 20 minutes maximum per exercise.

Some people may feel sore after using the mat or experience lightheadedness. Pay attention to your body, move slowly and adjust your usage accordingly.

There are many exercises to use with the mat set - Lying down, sitting, or standing all enable you to target different areas of stress and tightness. We have shown a few examples to get you started:

Lower Back / Lumbar / Hips

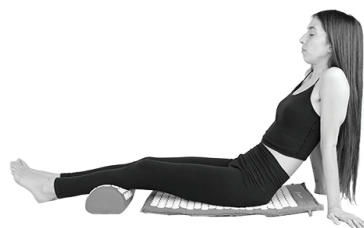Calves / Glutes

Back / Neck / Shoulders

Glutes / Lower Back / Feet

Please consult our website for more options

Caring for your mat and pillow set

Do not machine wash or dry as that may damage the acupressure points.

Remove the cover and hand wash gently with the cover inside out. Air dry.

1 Year Limited Warranty


The Acupressure Mat and Pillow Set has a Limited 1 Year Warranty to be free of manufacturing defects in materials and workmanship under normal applications, for the original owner. If this product becomes inoperable due to defect, return the product with all components and proof of purchase to the address listed below. This warranty does not cover any shipping/transport costs. This warranty does not apply to normal wear and tear or if the product is subject to misuse, neglect, rough handling or damage.

Ship the unit prepaid and insured (at owners option) to:

Thermor Ltd.

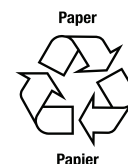Attn: Returns Department
16975 Leslie Street
Newmarket, ON L3Y 9A1

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THERMOR LTD.