iComfort® ic4500

**MASSAGE CHAIR
FAUTEUIL DE MASSAGE**

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Thank you for purchasing this product. Please read this instruction manual carefully before operating this massage chair.

Please pay attention to the important safety warnings and keep this manual for future use.

Note: EBI reserves the right to revise the design and description of this product without notice.

CONTENT

SAFETY WARNING
RESTRICTIONS
MAINTENANCE
TROUBLESHOOTING
FEATURES AND BENEFITS
INSTALLATION
COMPONENTS
OPERATION
MOVING THE CHAIR
FLOOR PROTECTION
PRODUCT SPECIFICATIONS
WARRANTY
FRANÇAIS

1. SAFETY WARNING

- This product is not recommended for children.
- Disconnect the chair before cleaning.
- Please pull out the AC power cord after use or before cleaning to avoid injury or damage to the product.
- Please operate this product according to the instructions in this manual.
- Use only accessories recommended by iComfort®.
- This product is designed for indoor use only.
- Recommended use time is 20 minutes or less.
- Please avoid use if the chair is damaged and call your iComfort® technician or dealer.
- Make sure behind the chair is well ventilated and not obstructed.
- Please make sure small objects are not dropped into the chair mechanism.
- Use of this chair under the influence of drugs or alcohol is not recommended.
- Do not use wet hands to pull the AC cord.
- Do not spill liquids onto this product to avoid electric shock or cause damage to this product.
- Maintain your chair and avoid damaging the wires.
- Do not use a wet cloth to clean the electric parts.
- In the event of a power failure disconnect the chair.
- If the chair works abnormally, stop using it, disconnect the chair and call your iComfort® service technician or customer service representative at 1-866-967-7333 ext 227.
- This chair should be used only in a well ventilated, normally dry, and dust-free environment.

2. RESTRICTIONS

People with osteoporosis, heart disease and heart pacemakers, pregnancies, injuries, or any other health-related issues should consult their doctor before using this massage chair.

Stop immediately the use of this massage chair if you feel unwell and consult your doctor.

3. MAINTENANCE

- This product should only be maintained by iComfort® certified technicians. Otherwise, EBI reserves the right to void the warranty.
- Disconnect the power cord after if you do not plan to use the chair for a long period.
- Make sure the wall outlet is in good working order before use.
- If this chair will not be used for an extended period, store in a dry and dustless environment.
- Avoid chair exposure to high temperatures and direct sunlight.
- Clean this product with a dry cloth. Do not use thinner, benzene or alcohol.
- When moving the chair on a wood floor, rough ground or in a narrow space two people should lift the chair together.
- We recommend using the massage chair intermittently and not for long periods.

4. TROUBLESHOOTING

- It is normal for the motor to make noise while in use.
- If the controller doesn't normally work, please check that the plug and socket are properly connected, and the switch is turned on.
- When the programmed working time is over, the chair will power off automatically.
- After several uses, the chair might power off because of overheating. It can be used again after half an hour of cooling time.

5. FEATURES AND BENEFITS

Enjoy the new iComfort® ic4500 massage chair. Its vertical to horizontal movable quad-rollers will massage from your neck all the way to your thighs. Very versatile, the iComfort® ic4500 massage chair reproduces massage actions with the sensation of a professional massage in the comfort of your very own home! It even includes shoulder, seat, forearms, and foot massages.

6. INSTALLATION

The iComfort® ic4500 massage chair comes partly assembled.

- Assemble this massage chair according to the installation instructions provided in this guide.
- Check the parts list carefully to make sure nothing is missing.
- Sitting in the massage chair before it is completely assembled could result in permanent damage to the chair and personal injury.
- **This massage chair is equipped with a cord having a grounding conductor. The cord must be connected into an appropriate and grounded AC wall outlet.**
- If the massage chair fuse is burned out, please use the chair's main On/Off button to turn it off and unplug the AC power cord. **Can call iComfort® at 1-866-967-7333 ext 227 if you need help.**



CAUTION: Never substitute a missing part with one of your own. If a part is missing, contact customer service at: 1-866-967-7333 ext 227.

Unpacking


- 1) Remove the massage chair, cushions, and accessories from the main box.
- 2) Remove all plastic wraps.

Controller bracket installation instructions

- 1) After the screw thread on the support knob of the controller rest assembly aligns to the screw thread position of the support mounting seat on the armrest of the massage chair, just turn clockwise and screw tightly.

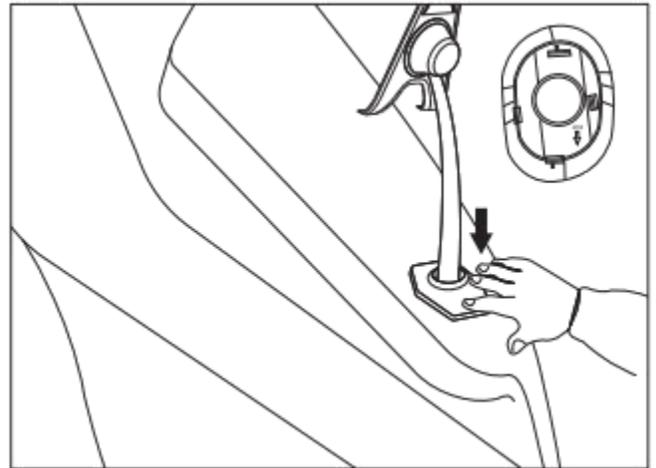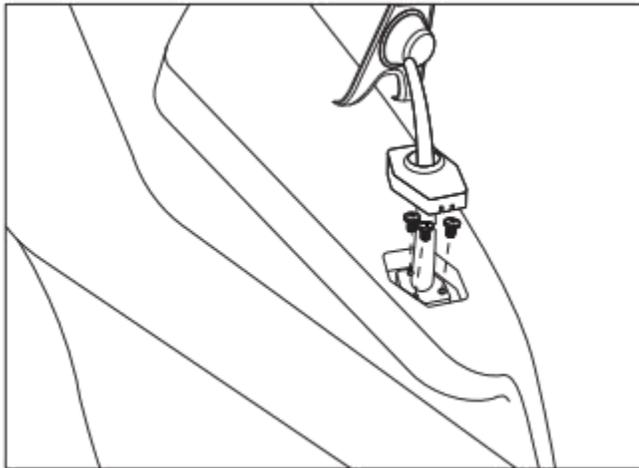
2) Connect the main remote-control cable to the connector under the seat.

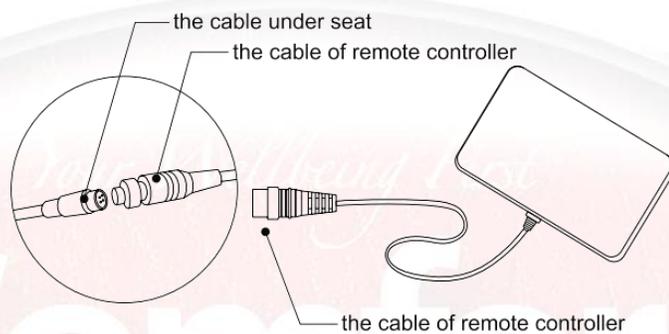Positioning

Make sure there is at least 50 cm in front and behind the massage chair before use.

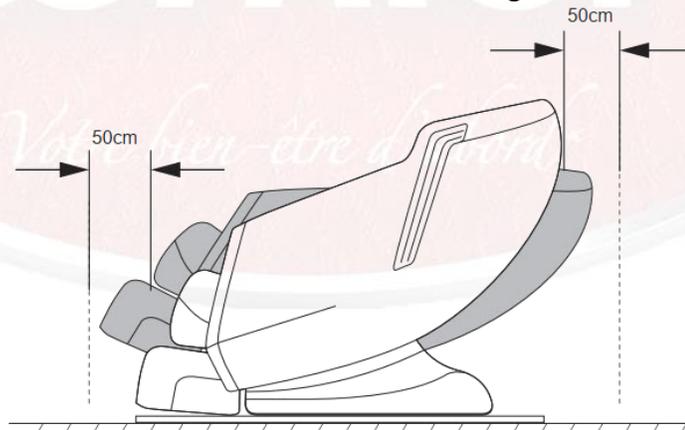Lifting

At least two people are recommended to lift and move the chair to avoid injury.


Back and head cushions

Install the back and head cushions using the cushion's zippers.


Final Installation


- 1) Insert the AC power cord into the massage chair power socket and the other end into the AC wall power outlet.

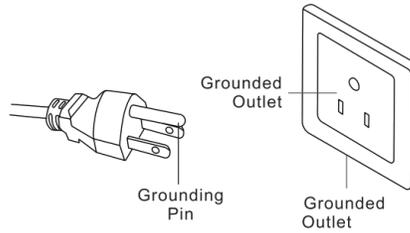
- 2) Turn on the main power switch at the bottom, behind the massage chair.

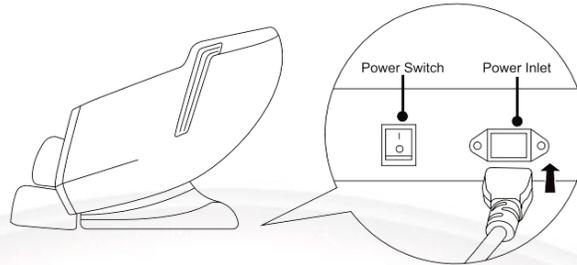7. COMPONENTS

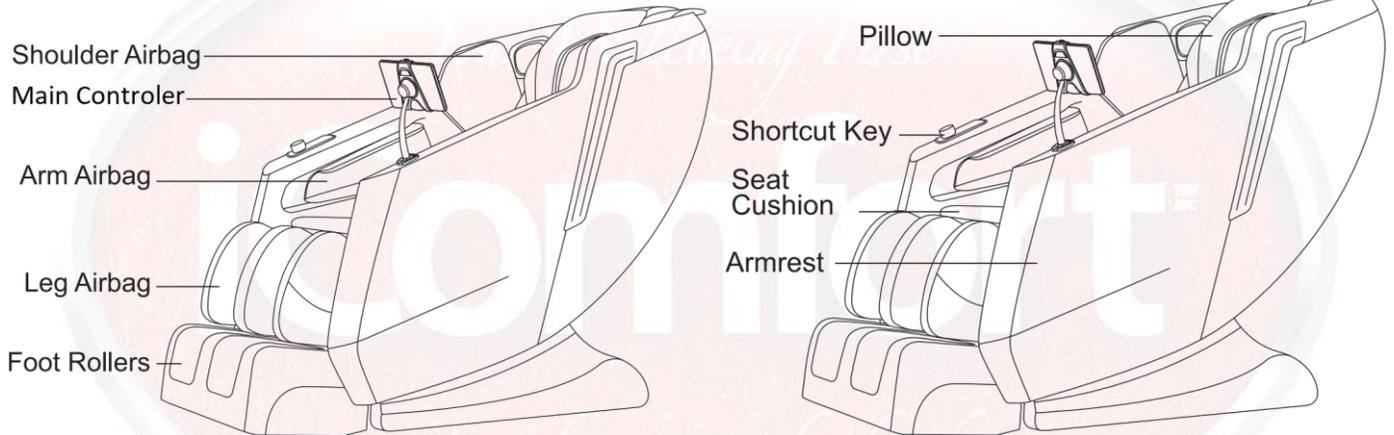Main Remote Control

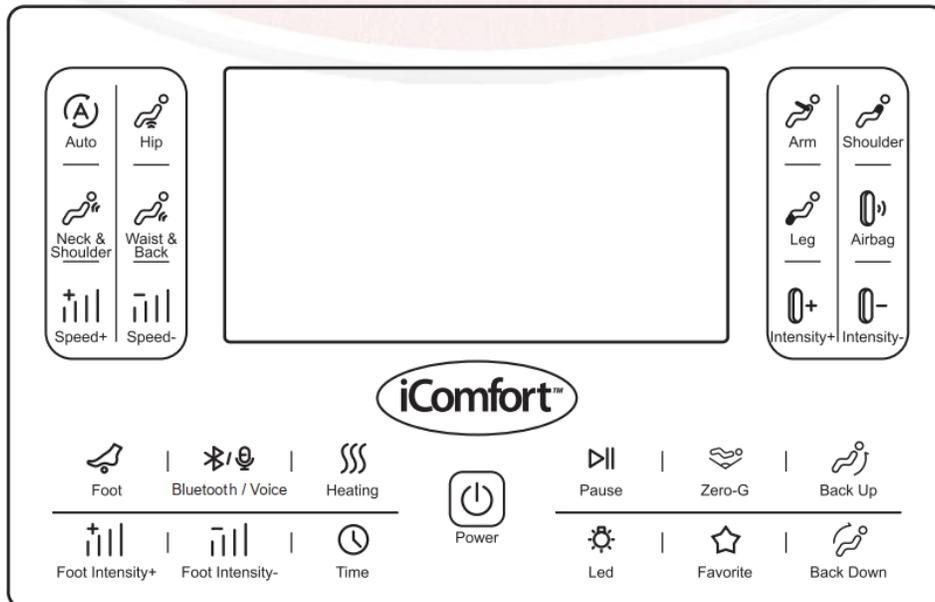8. OPERATION

Once your massage chair is connected to an AC wall outlet, you can turn on the chair by power on/off key.

Press the power key  on the main remote control to start the massage.

Automatic massage programs

Press repeatedly on the automatic massage button to activate one of your 6 preferred automatic massage program.

Refresh	Full body massage, the intensity of the roller and airbags is strong.
Relax	Full body massage, the massage intensity of the roller and airbags is medium.
Yoga	Full body massage, this massage simulates a Yoga exercise.
Neck & Shoulder	Partial massage: the focus is on improving neck pain and shoulder muscle strain. The massage is more focused on the neck and shoulder.
Waist & Back	Partial massage: the focus is on the waist and back to improve back strain.
Dream	Full body massage, the intensity of the rollers is light.

Press this button to turn on/off massage rollers on the hip.

Press this button to turn on/off massage rollers on the shoulder and neck.

Press this button to turn on/off massage rollers on the waist and back.

Press this button (+) to increase the massage speed.

Press this button {-} to decrease the massage speed.

Press this button to turn on/off the arm airbags massage.

Press this button to turn on/off the shoulder airbags massage.

Press this button to turn on/off the leg airbags massage.

Press this button to turn on or off the airbag massage function.

Press this button {+} to increase the message airbag.

Press this button {-j to decrease the massage airbag.

Press this button to turn admit the roller massage function.

Press this button to increase the massage speed of the roller massage function.

Press this button to decrease the massage speed of the roller massage function.

Bluetooth is enabled by default. A long press on the button will turn on/off the Bluetooth. Voice commands are enabled by default, press to turn on/off it and adjust the volume of the voice on the massage chair speakers. Open "Bluetooth" on your smart device and Bluetooth will automatically search for the corresponding Bluetooth model (ic4500). Click the "pair" button. Turn on "Music Player" in the smart device and select the song to play.

Press this button to tum on/off the heating function.

Press this button repeatedly to set the massage time (15/20/25/30min).

Press this button to pause all massage functions. Press again to resume the massage.

Press this button to turn on/off the LED light.

Press this button repeatedly to activate one of the three zero gravity positions and deactivate it.

While a massage program is being performed, hold this button for 3 seconds to record this program to your favorite program, then you can short press this button to start this program. (Notes: only one program can be saved). The previously recorded program will be cleared when a new one is added.


Press this button continuously to raise the backrest and lower the leg rest. Release it to lock the preferred angle.

Long press this button to lower the backrest and raise the leg rest. Release it to lock the current angle.

Armrest Remote Control

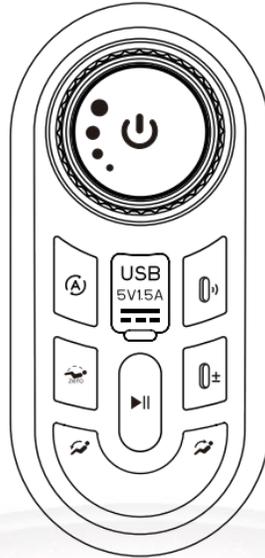Press this button to activate or deactivate the massage chair.

Rotate clockwise to change the massage chair angle to a more horizontal position. Rotate counterclockwise to change the massage chair angle to a more vertical position. There are three possible adjustments.

Press this button repeatedly to activate one of the automatic program functions. There are 6 in total. The display on the remote will show the program you have chosen.

Press this button to pause all massage functions. Press again to resume the massage.

Press this button to activate or deactivate the massage airbags.

Press this button to adjust the airbag intensity from levels 1 to 3.

Press repeatedly this button to activate one of the three zero gravity modes.

Press the backrest down button to lower the backrest. The leg rest will go up at the same time.

Press the backrest up button to lower the backrest. The leg rest will go down at the same time.

USB

You can remove the protective cover tab to recharge your phone using a USB cable (*cable not provided*).

Voice Control (English Only)

When the massage chair is turned on or in standby mode; say 'Hello massage chair' to enable the voice control feature. Use the following commands to activate or deactivate the feature of your choice:

COMMAND WORDS	RESPONSE BY CHAIR
Hello massage chair	Hello Master
Turn on massage chair / Power On	Please select your favorite program to start your.
Turn off massage chair / Power off	During shutdown, please leave or restart the massage chair after the chair is completely reset. Thank you and hope you enjoy the massage.
Pause massage	Massage paused
Continue to massage	Massage resumed
Turn on Bluetooth	OK, Bluetooth on
Turn off Bluetooth	OK, Bluetooth off
Turn off foot roller	OK, Foot roller off
Turn on foot roller	OK, Foot roller on
Turn on heating	OK, Heating on
Turn off Heating	OK, Heating off
Foot roller off	OK, Foot roller off
Foot roller on	OK, Foot roller on
Increase foot roller intensity	OK, foot roller intensity up
Decrease foot roller intensity	OK, foot roller intensity down
LED light off	LED light off
LED light on	LED light on
OK, Back heating on	OK, Back heating on
Back heating off	OK, Back heating off
Zero gravity	OK, Zero gravity
Back up	OK, Backrest up

COMMAND WORDS	RESPONSE BY CHAIR
Back down	OK, Backrest down
Start airbag massage	OK, Start airbag massage
Stop airbag massage	OK, Stop airbag massage
Massage for five minutes	OK, Massage for five minutes
Massage for ten minutes	OK, Massage for ten minutes
Massage for fifteen minutes	OK, Massage for fifteen minutes
Massage for twenty minutes	OK, Massage for twenty minutes
Massage for twenty five minutes	OK, Massage for twenty five minutes
Massage for thirty minutes	OK, Massage for thirty minutes
Relaxing massage	Start Relaxing massage
Yoga massage	Start Yoga massage
Dream massage	Start Dream massage
Refreshing massage	Start Refreshing massage
Waist&Back massage	Start Waist & Back massage
Neck&Shoulder massage	Start Neck & Shoulder massage
Lower body massage	Start lower body massage
Upper body massage	Start upper body massage
Volume up	Volume up
Volume down	Volume down
Maximum volume	Maximum volume
Minimum volume	Minimum volume
Increase strength	Increase back roller intensity
Decrease strength	Decrease back roller intensity
Increase air intensity	Air intensity up


Note: When using the voice control function, the speech speed should not be too fast, and the voice commands should be issued after 1 to 2 seconds.

Massage end

During a massage, press the remote-control power button  to turn the equipment off. The massage chair will immediately turn off all massage functions, reset the backrest unit and the leg unit, or when the massage timing time arrives, all massage functions will automatically be turned off.

MOVING THE CHAIR

Make sure that the AC power cord is pulled out of the chair. Turn the back backward to a proper angle so the centre gravity falls onto the casters. Pull the massage chair forward or backward as you wish.

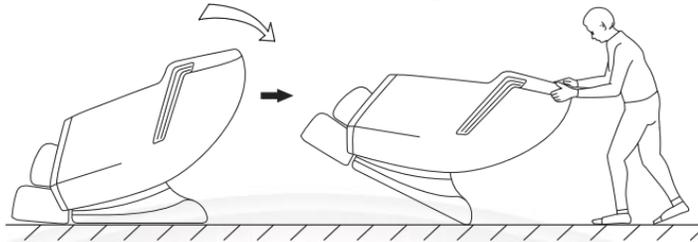FLOOR PROTECTION

Prolonged use of this massage chair can damage wood floors. To prevent permanent damage, place a carpet or mat under the chair.

PRODUCT SPECIFICATIONS

Rated voltage: 110-120V~ 60Hz.

Rated power: 100W.

WARRANTY

LIMITED ONE (1) YEAR WARRANTY

This product is guaranteed against all manufactured faults and defects for a period of one (1) year from the date of purchase.

For repairs or product defects, please call our customer service department at: 1-866-967-7333. Please have your proof of purchase on hand. If we are not able to replace your defective part, we will send you a service technician for repairs.

Important notes: Any damage caused by an accident, or an abusive usage are NOT covered under warranty. We also decline all responsibility concerning indirect or consecutive damages. In certain Canadian provinces, it is not permitted to limit indirect or consecutive damages. Therefore, it is possible these limitations are not applicable to you. This warranty gives you specific legal rights. You may also have other legal rights which vary from one province to another.

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