Setting Alarm Tone


Your new pill dispenser features 3 different alarm tones as well as a silent (light alert only) option to accommodate users with different hearing levels. At the top of the display, symbols "A", "B", and "C" indicate the selected tone. When set to silent mode, none of the three symbols will be displayed.

- To change alarm tone, begin by pressing and holding Button 3 (10).
- While holding down Button 3 (10), press and hold Button 1 (8) as well.
- Alarm tones will sound for a few seconds each, allowing you to choose.
- To select a particular tone, simply release all buttons when that tone sounds.


Testing Your Dispenser

- Press and hold Button 2 (9) for 5 seconds. "CH" will be displayed. Device is now in testing mode.
- Press Button 1 (8), and device will automatically proceed through following tests:
 - "CH 1" will be displayed, and all display characters will then light up.
- "CH 2" will be displayed, and tray will then rotate one compartment.

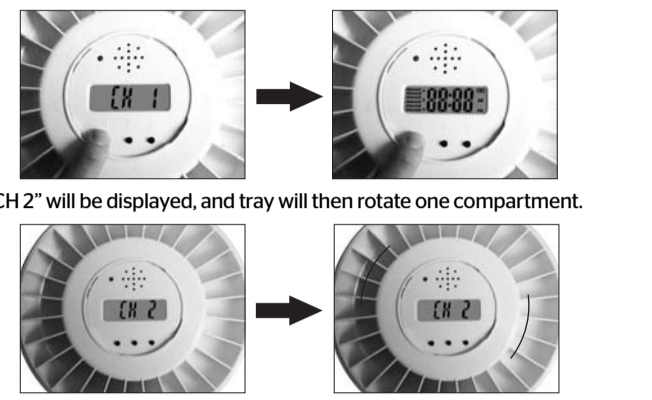c. "CH 3" will be displayed, and alarm will then sound and red light will flash.

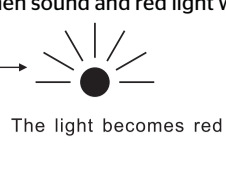
The light becomes red

- Lift dispenser and turn it over as if pouring medication into palm. Alarm and red light will shut.
- To exit test mode and return to clock mode, press Button 3 (10). Note! If dispenser is not turned over within 10 seconds, "CH 4" will be displayed. If tray does not rotate properly, motor will stop and "CH 6" error message will be displayed. In event of test malfunction, reset by removing and replacing batteries.

Locking Your Dispenser

- Locate Cover Lock (12) on bottom of unit.
- Insert included key and turn counterclockwise.
- To unlock, insert key and turn clockwise.

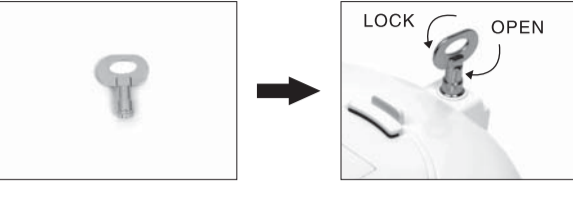Important Notes

Your pill dispenser is equipped with a special sensor. To stop the alarm and dispense the pills, lift the dispenser and turn it over to pour medication into your palm or a suitable container. The alarm and red light will then shut automatically.

A battery symbol with an "X" over it indicates that the batteries are low and should be replaced to ensure proper function. To save power, pill dispenser will stand by and shut down if no operation in several minutes. To wake it up, tilt up or press the Button 3 (10).

Do not place your pill dispenser in direct sunlight or strong ultraviolet light. It may cause operational malfunction.


Customer Service:

support@camarketing.com
1-888-390-7622

Distributed by
C&A Marketing, Inc.
114 Tived Lane East
Edison, NJ 08837 USA

C&A Marketing (UK) LTD.
167 Hermitage Road
Unit 1 Crusader Industrial Estate
London, UK, N4 1LZ

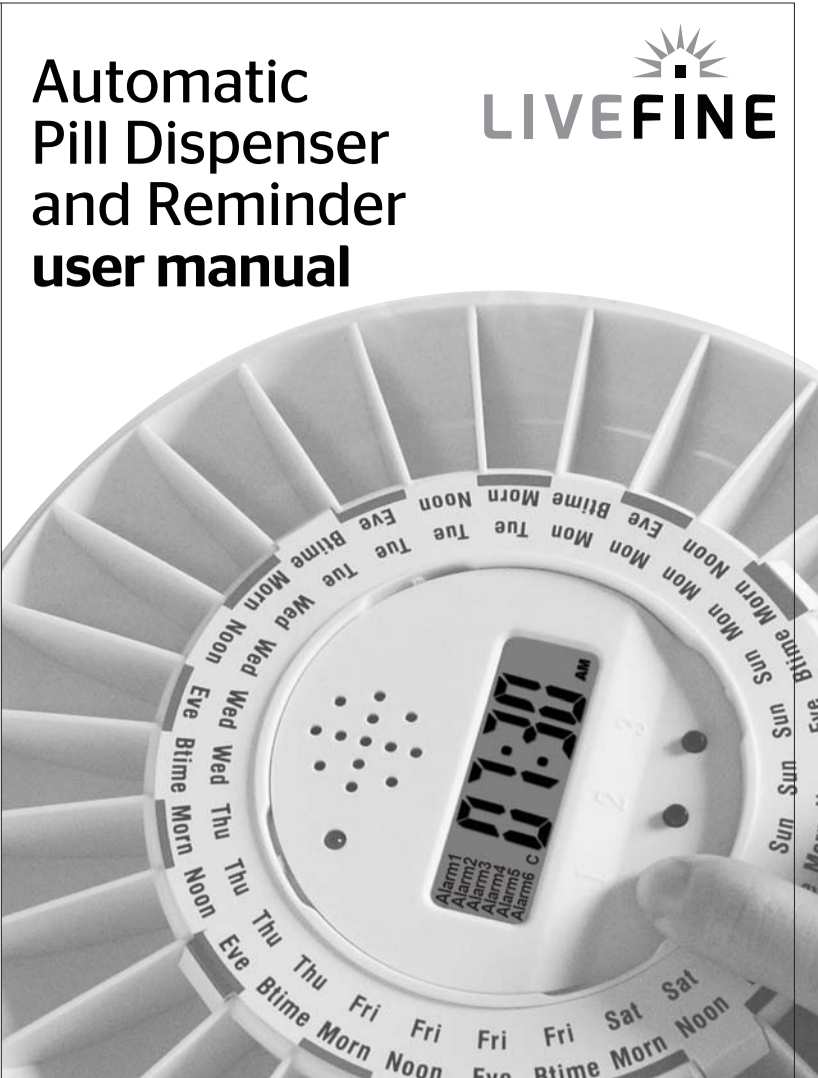Made in China
LIVEFINE is a trademark
of C&A IP Holdings, LLC
© 2016 All right reserved

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de C&A IP Holdings, LLC