Instant Ear Thermometer

Model 09-340

THESE BRIEF INSTRUCTIONS ALLOW YOU TO QUICKLY BEGIN USING YOUR INSTANT EAR THERMOMETER. PLEASE REFER TO THE INSTRUCTION MANUAL FOR MORE DETAILED INSTRUCTIONS COVERING THE USE AND FEATURES OF THE THERMOMETER.

User Guide:

1. Clean the lens at the tip of the thermometer, using alcohol and a soft cotton swab. Cleaning the lens before each use ensures the most accurate reading and prevents cross-contamination.
2. Press and release the red "SCAN" button. The thermometer will perform a brief self-test, during which all indicators on the display will be shown, *Fig. 1*, followed by the last reading taken.
3. The thermometer will beep, indicating that the unit is ready for measurement. Gently insert the probe into the ear canal until a slight resistance is felt, *Fig. 2*.
4. Press and release the SCAN button again, while keeping the probe in place. A beep will sound when the measurement is complete.
5. The temperature reading will be shown on the display, *Fig. 3*.
6. After about one minute after use, the thermometer will automatically beep and shut off.
7. Clean the thermometer lens again after use.

Fig. 1

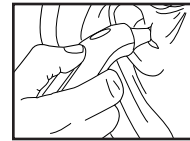
Fig. 2

Fig. 3

Please Note:

- Ear temperature readings are equivalent to oral readings.
- For accurate results, wait at least 10 seconds between measurements.
- If the temperature measurement is below 100.4° F, the beep indicating the reading is complete will be a single, short beep. A happy face will be displayed next to the reading.
- If the reading is 100.4° F or above, the unit will sound a long beep, followed by two short beeps. A sad face will be displayed.
- If, after pressing the SCAN button to activate the thermometer, a flashing "°F" does not appear, corrective action may need to be taken before use. See the Instruction Manual for interpretation of error messages and troubleshooting tips.
- Illness may be present even without an elevated temperature, particularly in the elderly, in infants under 3 months of age, and those with compromised immune systems. If illness is apparent without a fever, consult with a healthcare professional.
- Certain ear conditions, such as ear canal drainage, ear infection, or prior ear surgery may make ear thermometer readings impossible or inaccurate. If special ear conditions are present, consult with your healthcare professional.
- Never leave a child unattended with a thermometer.

Please read the Guidebook for this thermometer prior to taking your first reading. If you have any questions or are experiencing difficulties, please call the toll-free Customer Service Help Line for further assistance at 1-866-326-1313.

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