


CÉLAGE

B I O S C I E N C E

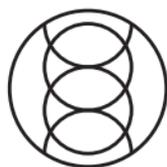

THE CÉLAGE TREATMENT SYSTEM
INSTRUCTION MANUAL



TABLE OF CONTENTS

How to use The Célage Treatment System	2
Three-Day Treatment Cycle.....	2
Treatment Days.....	4
Non-Treatment Days	7



WARNING:

Products for external use only. Keep out of reach of children. If you experience irritation, rinse from skin and discontinue use. Avoid direct contact with eyes. If contact occurs, rinse thoroughly with water.

HOW TO USE THE CÉLAGE TREATMENT SYSTEM

The Célage Treatment System is a scientifically-advanced skincare regimen designed for at-home use. This instruction booklet provides a step-by-step guide to help achieve optimal results. Like any treatment designed for long-term benefit, it is important to maintain a regular schedule, adhering to the following instructions as closely as possible.

Everyone who uses The Célage Treatment System begins with their own expectations and is likely to experience different results. When performed consistently and according to these directions, our findings suggest a '30/60/90' distribution of results. This means that approximately 30% of Célage Treatment System users experience improvement to skin within 30 days, 60% of users experience improvement within 60 days, and 90% of users experience improvement within 90 days.

THREE-DAY TREATMENT CYCLE

The Célage Treatment System is based on a 3-day cycle; 1 Treatment Day followed by 2 Non-Treatment Days. If you happen to miss a Treatment Day, simply add a Non-Treatment day and perform a Treatment on the next day.


On Treatment Days: You'll exfoliate and perform a MediStamping™ Treatment, using products A-D as shown in the illustrated diagrams of products below (*shown on page 3-5*).

On Non-Treatment Days: You'll perform daily maintenance, using products D & E (*shown on page 6*).

Each shipment contains a supply of approximately 10-12 treatments spanning five to six weeks, when used as recommended.

The Célage Treatment System includes the following components:

- (A) The PreTreatment Exfoliant
- (B) MediStamp Disinfecting Pads
- (C) The MediStamp 25™ Microneedling Tool
- (D) The Treatment Serum
- (E) The Daily Active Complex

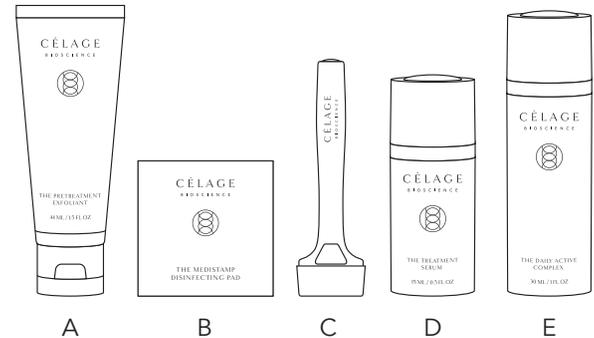

TREATMENT DAYS

MediStamping takes about five minutes to perform and is designed for use every third day. It consists of three steps: skin preparation, MediStamping and serum application. The Treatment may be performed at any time of day, but may be better suited to evenings as no other skincare item other than The Treatment Serum should be applied for at least one hour after MediStamping.

Step 1: Prepare with The PreTreatment Exfoliant

When: On Treatment days, immediately before the MediStamping Procedure.

What: This PreTreatment step prepares skin for MediStamping by gently sweeping away dead skin cells and surface oils.

How: After removing makeup or sunscreen, wet the entire area to be treated, (i.e., face only or face and neck). Place a pea-sized amount of The PreTreatment Exfoliant in the palm of your hand and mix with desired amount of cleanser or water. Apply to skin with gentle, circular motions, then rinse thoroughly. The exfoliant formula is concentrated and designed to be mixed with cleanser or water, so that you may adjust the amount to your desired level of abrasion.

Next, apply one pump of The Treatment Serum to face. Skin is now conditioned and ready for a MediStamping Treatment.


Step 2: The MediStamping Treatment

When: On Treatment days

Before the treatment: Disinfect The MediStamp by firmly pressing a freshly-opened Disinfecting Pad on the MediStamp surface. Hold for 5 seconds, then allow the MediStamp surface to air dry.

How to use the MediStamp correctly:

Holding the MediStamp perpendicular to the skin, press the stamp surface directly onto skin and lift away. Stamp with gentle but firm pressure, placing even pressure on all 4 corners of the MediStamp surface. Stamp in an even pattern so that all desired areas receive the MediStamping treatment. For optimal results, rotate the stamp head to create a random orientation of needles while adding a second or third layer of stamping to each desired area. You may choose to make more passes over specific areas of your face, like the crow's feet, forehead, cheeks or chin. Do not stamp the inner soft tissue of the lips or the eyelid skin. After use of the MediStamp, please replace the protective cap.

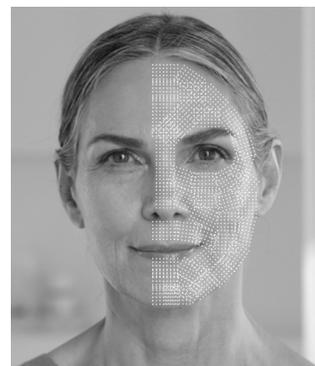

Frequency and Continuity: For best results, perform the MediStamping Treatment every third day, with each procedure separated by two Non-Treatment days. You may also wish to adjust the frequency based on skin sensitivity and personal preference. For example, if you prefer to Medistamp on the same two days of the week (rather than every third day), simply make sure that each treatment is separated by a minimum of two days (such as Sunday and Wednesday).

Please note: After a MediStamping treatment followed by The Treatment Serum, your skin may feel a tightening sensation. This is normal.

Step 3: Apply The Treatment Serum

Immediately after MediStamping, apply one pump of The Treatment Serum to treated areas and gently massage into skin. Wait one hour after MediStamping before using additional skincare or makeup on treated areas.

NON-TREATMENT DAYS

On Non-Treatment days, daily maintenance takes about 30 seconds in both morning and evening. Consistent use of these products in between MediStamping treatments helps maintain the look and feel of skin health throughout the week.

Step 4: Daily Maintenance with The Daily Active Complex

When: Morning and evening on Non-Treatment days

How: On Non-Treatment days, gently massage one pump of The Daily Active Complex to clean, dry skin - both morning and evening - with a special focus on areas of treatment. Using gentle, upward and outward circular motions not only ensures better application, but also helps massage the Complex into the skin. Immediately follow with your regular moisturizer.

The Daily Active Complex is a highly concentrated and active serum designed to maintain the progress of The Célage Treatment System throughout the week. *It is not a moisturizer, so please continue use of your favorite moisturizer after applying the Daily Active Complex.*

Step 5: Nightly Eye Maintenance with The Treatment Serum (Optional)


When: Evenings on Non-Treatment days

How: Using the soft pads of your fingertips, gently massage one pump of The Treatment Serum into the delicate skin surrounding both eyes. Take special care to keep product from entering the eye. If contact occurs with the eye, rinse thoroughly with water.



Visit our website for instructional videos on how
to best use The Célage Treatment System

WWW.CELAGE.COM

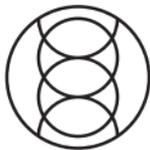

QUESTIONS?

Give us a call at 800-589-5536
Monday - Friday between 9 AM - 6 PM PST,
and our Concierge will be happy to help.

Célage Bioscience, Inc.
Portland, OR 97209

WWW.CELAGE.COM