

Changing the Replaceable Toothbrush Head


1. Firmly hold the DUO SONIC ADVANCED Toothbrush base section with one hand and the base of the brush section with the other hand.
2. Turn the brush section counter-clockwise and when the brush disengages, pull the head off the base.

Note: If you do not rotate the brush head you risk damaging the brush head or the main body.

3. Place the new brush head on the stem of the base and turn clockwise until the head is secure on the base.

Ordering extra brush heads

Replacement DUO SONIC ADVANCED brush heads are available. Please check with your place of purchase for more information or visit www.violife.com.

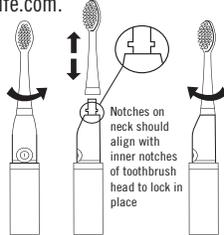


Replacing the Battery

Alkaline batteries are recommended

1. Grip the DUO SONIC ADVANCED Toothbrush at the base of the brush section with one hand and the bottom of the base section with the other.
2. Turn the base section counter-clockwise and pull down, exposing the battery box.
3. Remove the old battery. Insert a new AAA battery.
4. Replace the base cover and turn clockwise until the base section is securely in place.

Note: When battery is low, vibration gets weaker.

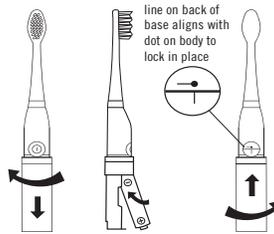

duo
sonic™
advanced
toothbrush



Congratulations on the purchase of your new Violife DUO SONIC ADVANCED Toothbrush!
The DUO SONIC ADVANCED Toothbrush operates at 2 speeds including 22,000 & 31,000 brush strokes per minute for cleaner, healthier teeth and gums.

Note: Features fading bristles (turn from blue to white) that indicate when to replace the brush head (approximately every 3 months).

How to Use

1. Put toothpaste on the brush.
2. Press button once for high speed, twice for low speed and third time to turn off DUO SONIC ADVANCED Toothbrush. 2 minute Quad timer consists of 30 second intervals with a short pause in between to remind you to move to the next quadrant of your mouth. Toothbrush automatically turns off after the 2 minute timer is complete.
3. Brush teeth as recommended by your dental professional.

How to Care for Your DUO SONIC ADVANCED Toothbrush

1. After every use, rinse the brush head thoroughly with water.
2. Shake excess water from bristles and dry your brush head before covering and storing.
3. Remove batteries if you won't be using the item for an extended period of time.
4. To clean your DUO SONIC ADVANCED Toothbrush, wipe surface with a damp cloth. Make sure the battery box is closed securely to prevent water from getting into the motor or battery compartment.

Do not immerse your DUO SONIC ADVANCED Toothbrush in water.

Keep out of reach of children and pets as this product may be dangerous for them.

Warning: Choking hazard - this product contains small parts. Seek immediate medical attention if swallowed or inhaled.

Questions?

email customerservice@violife.com or visit www.violife.com

ACTIVATE YOUR WARRANTY!

Visit www.violife.com to register your DUO SONIC ADVANCED Toothbrush and receive discounts and special offers on Violife products!

Limited 1 Year Warranty

Violife LLC warrants to the original purchaser of this product that it is free from defects in materials and workmanship for a period of 1 year from date of purchase.

If the unit becomes defective within the 1 year period, email us at customerservice@violife.com for replacement or troubleshooting assistance.

Within 1 year of purchase, Violife will repair or replace any part of the product we deem defective at no charge, provided the product has not been abused, misused, or damaged and was used according to the instructions.

This warranty gives you specific legal rights and you may also have other rights which vary from state to state.