



rockee the toothbrush that rocks


EN

Congratulations on the purchase of your new Violife Rockee toothbrush!
Rockee's unique design prevents it from ever laying on a surface, keeping the toothbrush and counter clean.

⚠ WARNING: CHOKING HAZARD - this product contains small parts.
Not for children under 3 years of age. Adult Supervision Required.

BEFORE EACH USE, inspect your Rockee toothbrush for any damage or loose bristles. Check that the brush head fits securely in place. Test connection outside your mouth. If the brush head is loose or detaches easily from handle, DO NOT USE.

- Seek immediate medical attention if swallowed or inhaled.
- If an unexpected event happens when operating, please report it to www.violife.com/customerservice.

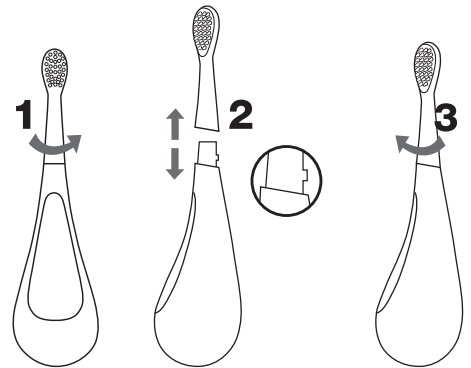


Changing the Replaceable Brush Head

- Twist brush head counterclockwise to disengage and pull up to remove the existing brush head.
- To place new brush head line up open slot on bottom of brush head to the notch on the base and rotate until secure.
- Make sure brush head is locked in place before use.

Ordering Extra Brush Heads

Please check with your place of purchase or visit www.violife.com for replacement brush heads.

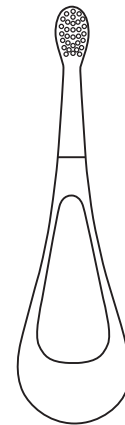

1. Rinse brush head after each use.
2. If needed, clean base with a damp cloth.
3. Replace heads every 3 months.