# Acupuncture + Air Compression Foot Massager with Heat

Thank you for

- choosing this product to complement your healthcare needs.
- Before using please read this user manual carefully for the correct usage and operation of this product.
- Please keep the user manual on hand for easy reference.

## Important Safety Instruction

When using the product, all instructions and basic precautions should always be followed, including the following:

### DANGER

To reduce the risk of electric shock:

Always unplug this product from the electrical outlet immediately after use and before cleaning.


### WARNING

To reduce the risk of burns, fire, electric shock, or injury to persons:

1. The product should never be left unattended when plugged in. Unplug from outlet when it is not in use and before putting on or taking off parts.
2. Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock or injury to persons.
3. Do not carry this product by it's supply cord or use cord as a handle.
4. Never operate the product with the air openings blocked. Keep the air openings free of lint, hair and the like.
5. Do not drop or insert any object into any opening.
6. Do not use outdoors.
7. Do not operate where aerosol (spray) products are being used or where oxygen is administered.
8. To disconnect, switch off all controls and remove plug from outlet.
9. Close supervision is necessary when this appliance is used by, or near children, invalids or disabled persons to ensure they are in no danger of the product.
10. Keep the cord away from heated surfaces. May cause serious burns. Do not use over sensitive skin areas or in the presence of poor circulation.
11. The appliance has a heated surface. Persons sensitive to heat must be careful using this appliance.
12. The appliance must only be supplied at safety low voltage corresponding to the marking of the appliance.
13. The appliance must only to be used with the power supply unit provided with the appliance.
14. This appliance is not intended for use by persons, including children, with reduced physical, sensory, or mental capabilities unless they have been given supervision concerning the use of the appliance by a person responsibly for their safety.
15. Special attention to components most susceptible to wear and tear to replace defective components immediately and/or keep the equipment out of use until repair.

# Name and Function of The Machine Parts

## Configuration Figure of Functions
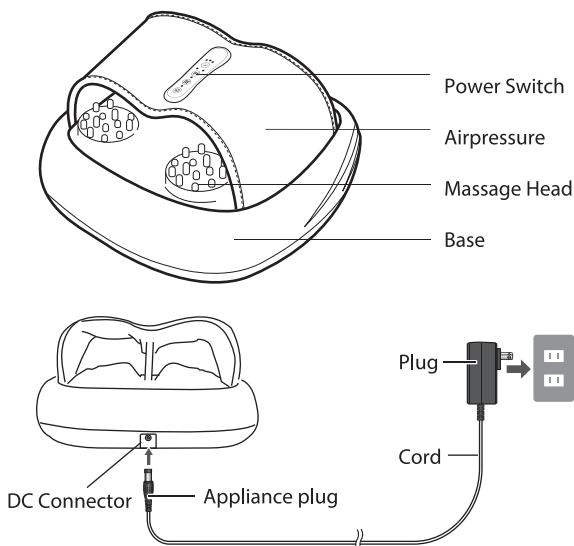
## How to Operate


- (Power ON/OFF):** Press this button, machine power-on. All keys are effective. Automatically setting time is 15 mins. The Acupuncture and Air pressure function turn on when power on, and the air-pressure of the feet is opened by default in weak intensity, Acupuncture speed defaults to slow speed. Press this button again, the whole machine is shut down, and all massages are stopped and the system is in standby mode. When time's up, machine stop massaging, enter standby state.
- Speed:** It is defaulted to be in Slow speed when power on, press the button to adjust the speed as you need. Speed switches every time when press button, shows by 2 LED indicator lights, respectively indicate the Slow and Fast speed.
- Heating:** It is defaulted to open heating function in startup state. Press button to close Led light and heating function. Followed by recycling.
- Air pressure:** It is defaulted to be in Weak Intensity on power-on state, Press button to switch intensity level to be in Medium, corresponding indicator light on.

1. Power button ON/OFF

2. Heating button: ON/OFF

3. Massage speed (Low / High)

4. Air Compression (Low / Medium / High)

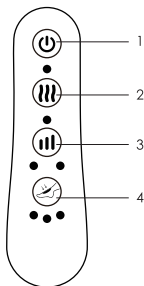This device is for household use ONLY and must not be used for commercial or medical purposes. If you have health concerns, consult your doctor before using this massage unit.

The massage unit must only be used in enclosed spaces. (Indoor use)

Do NOT use this product in damp rooms. (e.g. when bathing or showering)

- Do NOT wash!
- Do NOT dry clean!

## Operations for Use


Symptom	Cause	Solution
The product does not work	The plug is not inserted properly	Insert the plug properly and safely
	The power key is not started	Start the power key
	The product will be turned off automatically after it finishes a 15 minutes' massage	If you want massage for another time, press power key again
The product is hot or cannot be used suddenly in midway	You have performed several groups of 15 minutes' massage already.	Take off the power plug and let it cool. Always take off the power plug after using it.

If you still have any problem after using the methods above, please contact Redfern Ent.

## To clean and maintain the product

### 1. Cleaning

Before you clean, please make sure that the switch is in the off location, and the power plug has been pulled out from the socket. Please use dry cloth to clean the product.
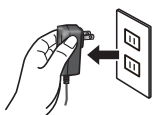
### 2. Maintenance

Keep the product clean and far from the wet environment.

### Attention

Please do not use benzene or thinner and other aggressive detergents to clean the product or it may fade.

## Technology parameter

Name: Acupuncture + Air Compression Foot Massager with Heat

Model No.: RFE-100

Rated Voltage: DC12V  $\overline{=}$  2.5A

Rated Power: 30W

DC 24V is not allowed.

### 1 YEAR LIMITED WARRANTY

Redfern Enterprises warrants this product to be free from defects in material and workmanship for a period of one (1) year from the date of original purchase, when utilized for normal household use, subject to the following conditions, exclusions and exceptions. Should any defect covered by the terms of the limited warranty be discovered within one (1) year, Redfern Enterprises will repair or replace (at Redfern's discretion) any defective part provided the unit is returned by the original purchaser, freight prepaid to:

Canada

Redfern Enterprises, 101-127 Commercial Drive, Calgary, AB, T3Z 2A7

Toll Free: 1-877-387-7770 FREE ext.221

Calgary Local: 403-735-5130 ext. 1

Email: info@redferment.com

USA

Redfern Enterprises, Unit 102, 4114 North Pecos Road, Las Vegas, NV, 89115

Toll Free: 1-877-462-5757

Las Vegas Local: 702-795-3314

Email: infous@redferment.com

Proof of purchase and \$29.95 for handling and return shipping charges must be pre-paid.

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