

3 in 1 Pro Massager **Owners Manual**



Redfern

Smart Products for Home + Health

Canada www.redferment.com

United States www.redfermentus.com

This manual should be kept on hand for easy reference.
It is your guide to a safe and efficient operation.

IMPORTANT SAFETY INSTRUCTIONS

When using the product, all instructions and basic precautions should always be followed, including the following:

DANGER

– To reduce the risk of electric shock:
Always unplug this product from the electrical outlet immediately after use and before cleaning.

WARNING

- To reduce the risk of burns, fire, electric shock, or injury to persons:
 - 1) The product should never be left unattended when plugged in. Unplug from outlet when it is not in use and before putting on or taking off parts.
 - 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
 - 3) Close supervision is necessary when this appliance is used by, or near children, invalids, or disabled persons.
 - 4) Use this product only for its intended use as described in this manual. Do not use attachments which are not recommended by the manufacturer.
 - 5) Never operate this product if it is dropped, has a damaged cord or plug, a malfunction, or if it is dropped into water. Return the product to service center for examination and repair.
 - 6) Do not carry this product by its supply cord or use cord as a handle.
 - 7) Keep the cord away from heated surfaces.
 - 8) Never operate the product with the air openings blocked. Keep the air openings free of lint, hair, and the like.
 - 9) Never drop or insert any object into any opening.

- 10) Do not use outdoors.
- 11) Do not operate where aerosol (spray) products are being used or where oxygen is administered.
- 12) To disconnect, switch off all controls and remove plug from outlet.
- 13) This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Should not be used unsupervised by children or people with sensory disorders
- 14) Children should be supervised to ensure that they do not misuse the product.
- 15) The appliance has a heated surface. Persons sensitive to heat must be careful when using the appliance as may cause serious burns. Do not use over sensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous. The appliance has a heated surface. Persons sensitive to heat must be careful when using the appliance.
- 16) The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- 17) The appliance is only to be used with the power supply unit provided with the appliance.


SAVE THESE INSTRUCTIONS

This product is intended for household use only


HOW TO OPERATE

1. Connect the power adapter to the electrical outlet.
2. Press the power button to begin the massage. Heating function is fully activated.
3. To deactivate the heating function, press the heating button and vice versa.
4. Press the "DR" key, the massage heads go anticlockwise and vice versa.
5. Press the power button to stop the massage at any time.

Massage position

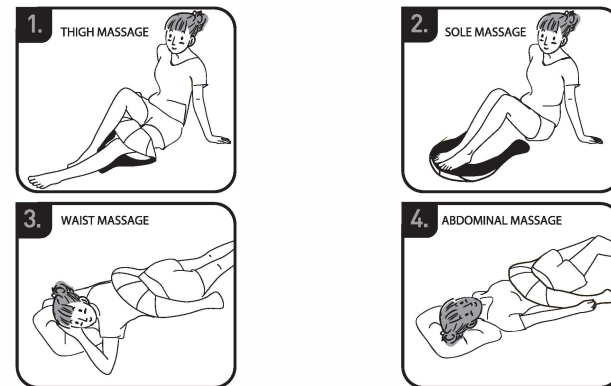

Name and function of the machine parts

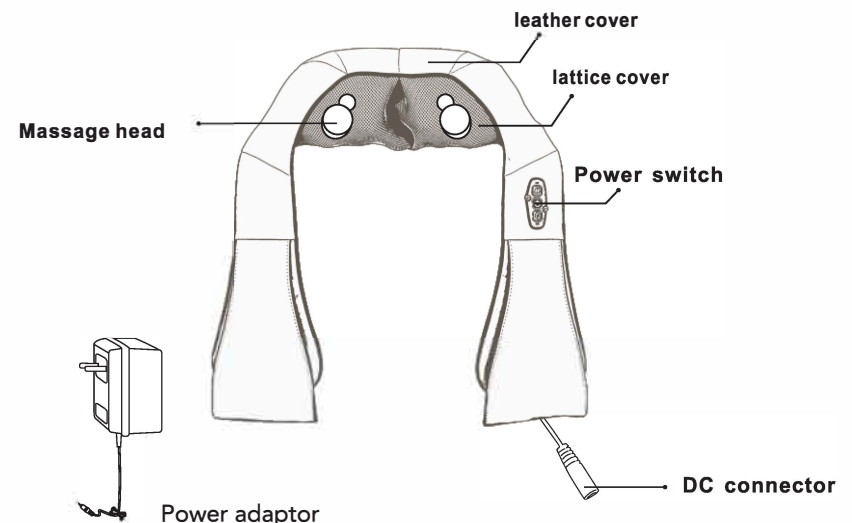

To clean and maintain the product

1. Cleaning

Before you clean, please make sure that the switch is in the off location. And the power plug has been pulled out from the socket. Please use dry cloth to clean the product.



Attention

- Please make sure that do not use benzene or thinner and other aggressive detergents to clean the product. Or the product may fade.

2. Maintenance

Make sure that the switch is off, and the power plug has been pulled out from the socket. Keep the product clean and far from the wet environment.

Technology parameter

Message Rated Voltage: 12V \Rightarrow 2.5A
Message Rated Power: 30W

Warranty

1 Year Limited Warranty

This product is warranted for 1 year from the original date of purchase, against defects in materials and/or workmanship when utilized for normal household use. During this period if the unit is proved defective it will be repaired or replaced at the discretion of the service technician.

The unit shall be returned, by original purchaser, freight prepaid with proof of original purchase date and \$39.99 return shipping/handling charges included. Redfern Enterprises does not accept C.O.D.'s and will not reimburse for any fees associated with the return or shipping of the item.

This warranty does not apply to any defect arising from a buyers misuse, negligence, failure to follow instructions, use on current or voltage other than that stamped on the unit or repairs not authorized by a Redfern Ent. technician will null and void the warranty.

If you believe that your product is defective, contact Redfern Enterprises Customer Services Department:

CANADA
101-127 Commercial Drive,
Calgary, AB T3Z 2A7
TF: 1-877-387-7770 (Toll Free)
info@redferment.com www.redferment.com

UNITED STATES
4168 N North Pecos Road, Ste 102
Las Vegas, Nevada, 89115
1-877-462-5757 ext. 1 (Toll Free)
infous@redferment.com www.redfermentus.com

Obtain a repair authorization number (RA Number) and bring in or send prepaid to the Service Center with proof of original purchase date, and prepaid return shipping of \$39.99. See contact info above for RA Number.