



Compression Massager



Redfern

Smart Products for Home + Health


Caution:

1. New users should limit the usage time to 5 min/section. You may increase the time usage by 5 minutes per section when you feel your body is well adapted.
2. Always unplug the unit from the electrical outlet when not in use.
3. Never insert finger or object inside the cushion.
4. Never use the product in the bathroom or in wet places.
5. Store in a cool dry place.
6. Do not hold the product by the wire as this may result in malfunctioning of the product.
7. Do not sit on this unit to prevent any break down.
8. It is advisable to use the device at least one hour after meals.
9. An adult should supervise use of the device by a child or the physically handicapped.
10. Do not try to repair the device yourself or through any unauthorized persons or agent, or it will void the warranty.
11. Should the device malfunction, disconnect it from the power source and contact the authorized agent immediately.

Description of each part



Controller

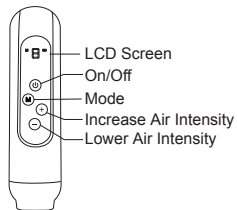

Product Specifications

Model No: BP-CM2022
Product Type: Electrical
Massager For use on Calfs, Feet
and Arms Lithium Battery:
900mah Voltage: DC5V 1A
Charging time: 1H
Using time: 1H Timer:
15 mins Gross Weight:
1.1KG USB line: 1m

USB Cable



Storage Bag



Usage



1. Fully charge the remote control for 1hr prior to using the product.
2. Unplug USB cord and store for next use.
3. Connect the rubber tube from the remote to the correct side massager cuffs (i.e LEFT or RIGHT).
4. Comfortably (not tight) wrap around the chosen body part with the Air Bag Cuffs.
5. Choose desired level or mode on remote

Controller Operation

1. Plug the rubber tubes into the corresponding hole in the remote control.
2. Turn unit ON to use or OFF when finished.
3. Choose your mode with the M
4. Adjust your intensity level + up or - down



Cleaning and maintenance



WARNING

- Switch off the device each time before cleaning.
- **Clean the device with a slightly damp cloth only.**
- Under no circumstances may liquid enter the device or accessories.
- Do not use the device again until it is completely dry.
- Do not clean the device in a dishwasher.

- Do not use any abrasive cleaning products or hard brushes.

- Whenever the product is not in use, please store it in a cool, dry place in the original packaging. Do not place other items on top as the pressure could damage the product.

Disposal



IMPORTANT

- Batteries should not be disposed of in household waste. Check with your local recycling to find out where to properly dispose of batteries.

Important Safety Instructions

1. Use item only for its intended use.
2. This item is intended for household and residential use only. This product is not approved for commercial applications.
3. If cord is damaged it must be replaced by manufacturer, its service agent or similarly qualified technician in order to avoid a hazard.
4. Not intended for use by children. Close supervision is necessary for any item being used near children. Unsupervised use can result in fire or personal injury. Do not leave item unattended while plugged in.

DISCLAIMER: Content is for reference purposes and is not intended to substitute any advice given by a physician, pharmacist, or other licensed health-care professional. You should not use this information as self-diagnosis or for treating a health problem or disease. Contact your health-care provider immediately if you suspect that you have a medical problem. Redfern Enterprises assumes no liability for inaccuracies or misstatements about products.

1 Year Limited Warranty

This product is warranted for 1 year from the original date of purchase, against defects in materials and/or workmanship when utilized for normal household use. During this period if the unit is proven defective it will be repaired or replaced at the discretion of the authorized service technician.

The unit shall be returned, by original purchaser, freight prepaid with proof of original purchase date and \$39.99 return shipping/handling charges included. Redfern Enterprises does not accept C.O.D.'s and will not reimburse for any fees associated with the return or shipping of the item.

This warranty does not apply to any defect arising from a buyers misuse, negligence, failure to follow instructions, use on current or voltage other than that stamped on the unit or repairs not authorized by a Redfern Ent. technician will null and void the warranty.

If you believe that your product is defective, contact Redfern Enterprises Customer Services Department:

CANADA

101-127 Commercial Drive, Calgary, AB
T3Z 2A7

TF: 1-877-387-7770 X 0 (Toll Free) Local 403.735.5130 X 0
info@redferment.com www.redferment.com

UNITED STATES

4168 N North Pecos Road, Ste 102
Las Vegas, Nevada, 89115

1-877-462-5757 X 1 (Toll Free) Local 703.795.3314 X1
infous@redferment.com www.redfermentus.com

For Service you must first obtain a Return Authorization Number (RA Number) Original Proof of purchase must be provided to customer service for Warranty. Product must be shipped prepaid, with RA# on it and return shipping of \$39.99. Return Shipping is subject to change as the rates change.