
BACKplus  **SHIATSU FOLDABLE
MASSAGE CHAIR**

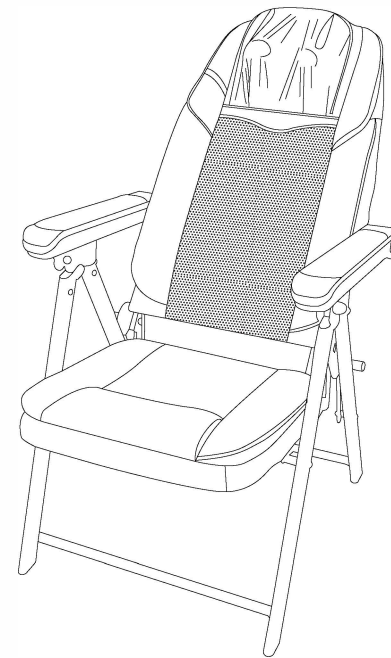MANUAL

General Warnings

1. Keep cord away from heated surfaces.
2. Never block the air openings of the product or place it on a surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint and hair.
3. Never use while sleeping.
4. Never drop or insert objects into any opening.
5. Do not operate where aerosol products are being used or administered.
6. Do not use on or near eyes or other highly sensitive areas.
7. Do not operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury.
8. Do not carry this product by its supply, or power cord, or use the cord as a handle.
9. Do not use directly on one concentrated area of the body. It is meant to be used in a continuous motion, constantly changing the massage area.
10. To disconnect, turn all controls to the off position then remove the plug from the outlet.
11. Do not perform any servicing on the product.
12. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.



To reduce the risk of burns, fire, electric shock, or injury:

1. Always unplug the product from the electrical outlet immediately after using and before cleaning.
2. Do not use while bathing or in a shower.
3. Do not place or store product where it can fall or be pulled into a tub, sink, or other body of water.
4. Do not place in or drop into water or other liquid.
5. Do not reach for an electrical appliance that has fallen into water. Unplug it immediately.


1

In the Box

Massage Chair

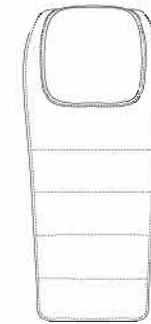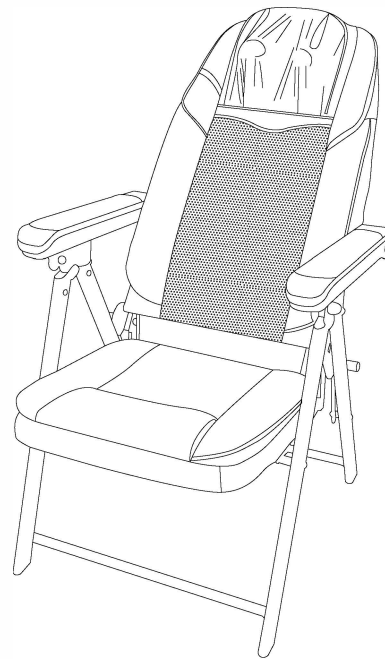
Mantle

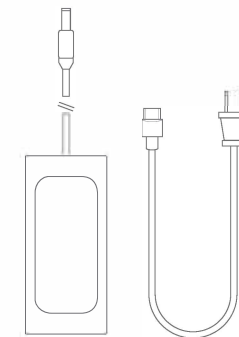12V $\overline{\text{---}}$ 5A Power Adapter and Cord

2


Location of Controls (Massage Chair)

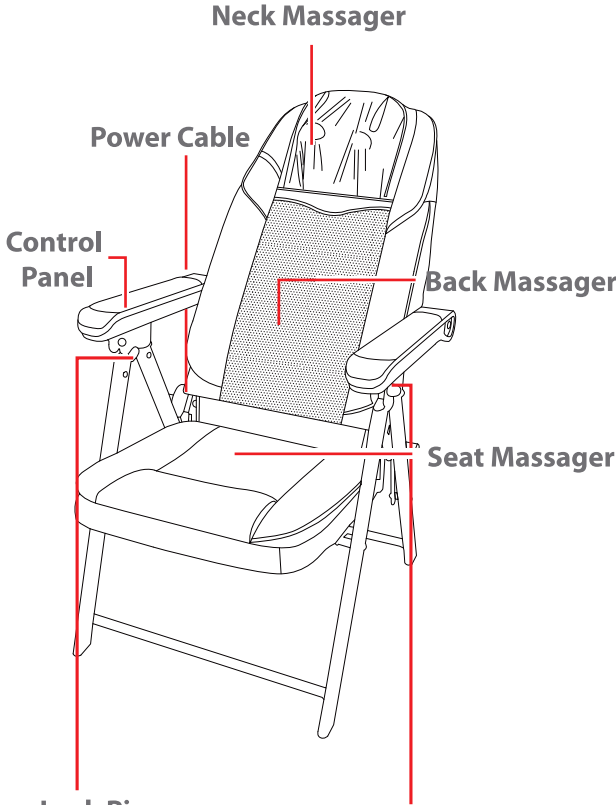Lock Pin
Adjust the back of the chair, turn the pin upward and lock into place.

Location of Controls (Control Panel)

1. **Neck Massage Adjustment**
Adjust the height of neck massage nodes
2. **Neck Massage Mode**
Switch between manual (left button) and automatic (right button) modes.
3. **Back Massage Mode**
Cycle through *FULL*, *UPPER*, or *LOWER* back massage.
4. **Vibration**
Cycle through 3 seat and massage intensities.
5. **Target**
Massage specific target area.
6. **Target Adjustment**
Adjust targeted massage area.
7. **Heat**
Adjust heat level.
8. **Power**
Power ON/OFF.

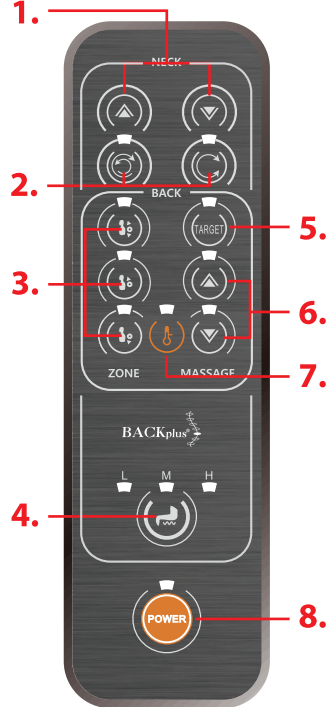
Setting Up the Massage Chair

Setup

Remove the massage chair from its packaging and fully unfold the chair, standing it upright. Make sure that you have room to unfold the chair and that it is located near a wall outlet.

Connect the pin of the **supplied 12V power adapter** into the **Power Cable** on the massager.

Plug the power cable into the power adapter and plug the other end into a wall outlet.

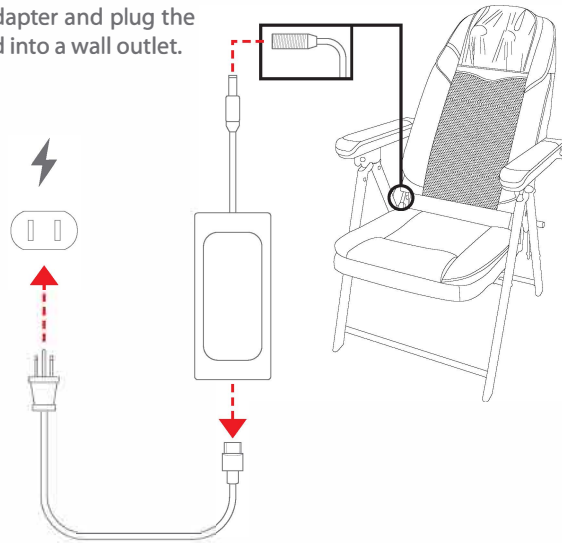Adjusting the Seat Back

To adjust the chair, undo the lock pin, adjust the back of the seat to your liking, and return the lock pin to place.

Using the Massage Chair

Using the Massage Chair

NOTE: When using the massage chair, always sit fully upright in the chair. Do NOT slouch or slide forward, do NOT sit on the edge of the chair, and do NOT tip the chair forward or backward.

After setting the chair up, the control panel on the right arm will allow you to power the chair ON/OFF, set your massage modes and intensity, activate heat, etc.

Powering ON/OFF

- Press the **Power** button on the control panel to turn the massage chair ON/OFF.
- If not in use, the chair will **automatically power off** after **15 minutes**.

Adjusting Heat

- Press the **Heat** button to adjust the heat function.

Adjusting Vibration

- Press the **Vibration** button on the control panel to cycle through the 3 seat massage intensities.

Neck Massage

- Press the **Neck Massage Adjustment** and **Mode** buttons to adjust the intensity and area of the neck massager nodes and to switch between **manual** (left button, below) and **automatic** massage (right button, below).

Back Massage


- Press the **Back Massage** buttons to adjust the different areas of the back to be massaged.
- Press the **Target** button to massage a specific area on the back and use the arrow buttons to adjust the targeted massage area.

Additional Safety Warnings

Storing the Massager

The massager should be stored in a dry, cool location, away from sources of water or heat. It is best stored in a secure box or an airtight bag to keep it free of dust.

To fold the massager for storage:

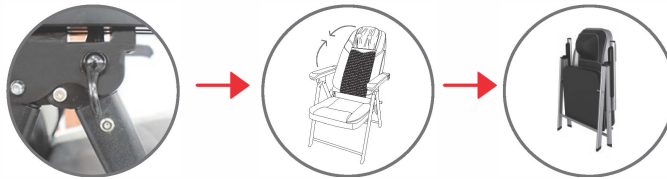Lock the pin on both sides

Fold the back and arm rest inward

Folded for storage

Cleaning the Massager

To clean the massager, unplug the massager and let the pad cool down if the heating function was just used.

Wipe the pad down with a soft, slightly damp sponge. Avoid using concentrated soaps, such as dish soap, as it will be very difficult to clean from the pads once applied.

- Do NOT allow any liquid or cleaning solutions to come into direct contact with the massager.
- Do NOT immerse massager in water.
- Do NOT use any abrasive cleaners or chemicals on the pad.
- Do NOT attempt to repair any damage to the massager.

Additional Safety Warnings

1. Consult your doctor before using this product.
2. Individuals with pacemakers should not use this item.
3. Never leave the product operating or unattended, especially if children are present.
4. Never use the product on open skin eruptions, swollen, or inflamed areas.
5. Do not operate after taking painkillers, sedatives, or alcohol. Consult your doctor if you are under these medications.
6. Do not continuously operate for 15 minutes or longer. It is recommended to turn the product off and allow it to rest (cool) for 30 minutes before re-using. This will prolong the life of the product.
7. The product is designed as a non-professional massager to soothe tired and aching muscles. It is not intended to be a substitute for medical attention. If your ailments and symptoms persist or worsen after use, contact your doctor immediately.

WARNING: Any individual who is pregnant, has a pacemaker, suffers from diabetes, phlebitis, and/or thrombosis, or is at an increased risk of developing blood clots due to recent surgery, should consult their physician before using a massaging device designed for home use.

One Year Limited Warranty

Caution: Any changes or modifications not expressly approved could void the user's authority to operate the equipment.

This product is warranted for 1 year from the original date of purchase, against defects in materials and/or workmanship when utilized for normal household use. During this period if the unit is proved defective it will be repaired or replaced at the discretion of the service technician.

The unit shall be returned, by original purchaser, freight prepaid with proof of original purchase date and \$39.99 return shipping/handling charges included. Redfern Enterprises does not accept C.O.D.'s and will not reimburse for any fees associated with the return or shipping of the item.

This warranty does not apply to any defect arising from a buyers misuse, negligence, failure to follow instructions, use on current or voltage other than that stamped on the unit or repairs not authorized by a Redfern Ent. technician will null and void the warranty.

If you believe that your product is defective, contact Redfern Enterprises Customer Services Department:


CANADA
101-127 Commercial Drive,
Calgary, AB T3Z 2A7
TF: 1-877-387-7770 (Toll Free)
info@redferment.com www.redferment.com

UNITED STATES
4168 N North Pecos Road, Ste 102
Las Vegas, Nevada, 89115
1-877-462-5757 ext. 1 (Toll Free)
infous@redferment.com www.redfermentus.com

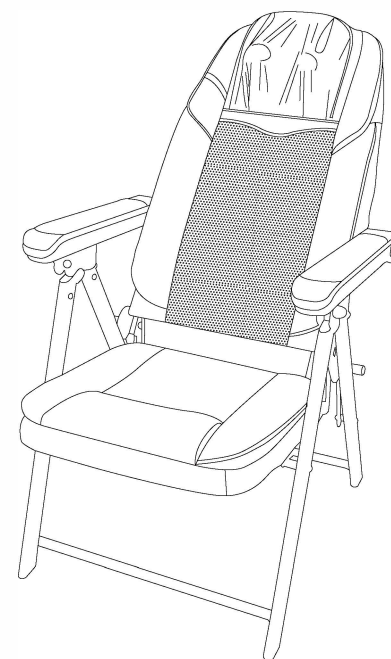Obtain a repair authorization number (RA Number) and bring in or send prepaid to the Service Center with proof of original purchase date, and prepaid return shipping of \$39.99. See contact info above for RA Number.

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