BACKplus® Percussive Massager

INSTRUCTIONS

BACKplus® Percussive Massager


A Cordless state of the are percussive massage that aids in relaxing sore stiff muscles. It helps to increase circulation and aids in range of motion to help improve the overall health of the body soft tissue.

- CORDLESS/RECHARGEABLE USB 5V**
Lithium-Ion Battery
- LIGHTWEIGHT DESIGN**

SPECIFICATIONS

Massage Heads:

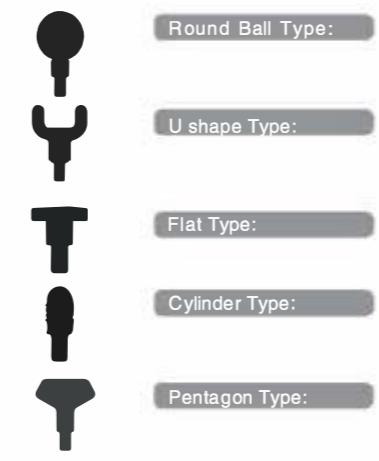LCD Display type:
Voltage:7.4V
Frequency:50Hz
Power:24W
Charger:uSB 5V

MAINTENANCE, CLEANING STORAGE

- Use a slightly damp towel to clean device and a dry with a soft cloth.
- For storage or travel remove battery and store in supplied carton box.

INSTRUCTIONS


PLEASE READ PRIOR TO USE

CHARGING

- Fully charge battery for up to two hours before first use either attached or detached from the device.
- To charge, connect DC end of the supplied USB cable to the battery's charging port and plug the USB end into a standard USB 5v supply.

OPERATING

- Press OA/OFF button to switch on or switch off the machine; Press the "+" button to intensify the power and press M button to change modes. Chargeable Type: Plug the charger port into the DC connector on the bottom of the machine and the USB port into the socket. When the indicator light shows LL, the charging process is normal and when the light shows LL, the charging process is over. The charging process will take about 2 hours.

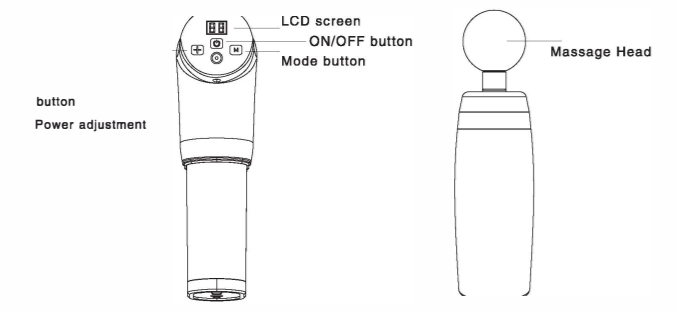- Press the button to turn on the machine
- Press it again to change speed. The speed change will show on the indicator light. You can choose any speed you want.

DISCLAIMER

WARNINGS & CAUTION

- Do not immerse in water. Keep away from liquids or heat sources. Keep ventilation Ports free from dust and debris.
- Do not remove screws or attempt to disassemble.
- Do not operate continuously for more than one hour. Allow device to rest for 30 minutes before re-using.
- Unplug the charger from device battery after charging or prior to use.
- For adult use only. Do not use if injured. Consult your doctor before using this product.

SAFETY INSTRUCTIONS

DANGER: TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

- For adult use only
- Use only on a dry clean surface of the body over clothing by lightly pressing and moving across the body for approximately 60 seconds per region.
- Use the it only on the body's soft tissue as desired without producing pain or discomfort. Do not use on head or any hard or bony area of the body.
- Use only the applicators that provide the best desired effect.
- Brusing can occur regardless of control setting or pressure applied. Check treated areas frequently and immediately stop at first sign of pain or discomfort.
- Keep fingers, hair or any other body part away from the shaft and back side of the applicator as pinching may occur.
- Do not place any objects into the ventilation holes of the unit.
- Do not immerse in water or allow water to enter the ventilation holes of the unit.
- Do not drop or misuse the unit.
- Carefully examine device before each use.
- Do not tamper or alter unit in any way.
- Never leave the unit operating or charging unattended.

PLEASE DO NOT USE THE THIS UNIT OR ANY VIBRATION DEVICE WITHOUT FIRST OBTAINING APPROVAL FROM YOUR DOCTOR IF ANY OF THE FOLLOWING APPLY
Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pacemakers, recent surgery epilepsy or migraines herniated disks. Spondylolisthesis, spondylol-ysis, or spondylosis, recent joint replacements or IUD's, metal pins or plates or any concerns about your physical health.

ONE YEAR WARRANTY

1 YEAR LIMITED WARRANTY

This product to be free from defects in material and workmanship for a period of one (1) year from the date of original purchase, when utilized for normal household use, subject to the following conditions, exclusions and exceptions. Should any defect covered by the terms of the limited warranty be discovered within one (1) year, manufacturer will repair or replace any defective part provided the unit is returned by the original purchaser, Purchaser to package item safely and prepay shipping. A prepaid return shipping charge of \$19.99 must be included with the product if repair is needed.

This warranty does not cover normal wear of parts and does not apply to any unit that has been tampered with or used for other than stated purposes. This warranty does not apply to any defect arising from misuse of the product, negligence, failure to follow manual, use on current or voltage other than that stamped on the product, alteration or repair not authorized by manufacturer, or damage due to faulty packaging or mishandling in transit (by a common carrier). Repair or disassembly by anyone other than a manufacturer authorized service center will void the warranty.

This warranty is extended to the original purchaser of the unit and excludes all other legal and/or conventional warranties. The responsibility of manufacturer, if any, is limited to the specific obligations expressly assumed by it under the terms of the limited warranty. In no event is manufacturer liable for incidental or consequential damages of any nature whatsoever.

For more information please contact (In Canada) Redfern Ent. at Toll Free 1-877-387-7770 ext 222 email info@redferment.com or (in the United States) at Toll Free 1-877-462-5757 ext 1 or infous@redferment.com

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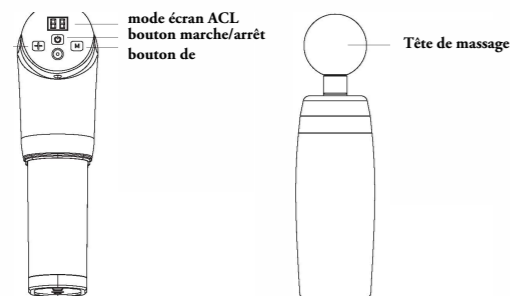
SANS FIL/RECHARGEABLE
Lithium-Ion Battery

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