



| ACCUMAST® *Plus*



USER GUIDE



ABOUT US:

FERA Diagnostics and Biologicals combines innovation and science to bring practical products to the animal health field that will enhance the health, productivity and welfare of pets and production animals.

We are a young and highly energetic research-driven company. Our growing, scientifically-proven product portfolio is designed to improve animal health, create value and contribute to the success of our customers.

FERA Animal Health, LLC was founded in 2014 by Dr. Rodrigo Bicalho, and named in honor of his two sons. The company's first product, an early version of Accu-Mast®, was developed in response to a large demand from Central New York dairy farms for a new mastitis testing methodology.

In February 2018, FERA Animal Health, LLC was restructured and rename FERA Diagnostics and Biologicals to reflect our growth and future plans.

Accumast plus is NEW SYSTEM for diagnosis of all types of Mastitis with the ability to differentiate Streptococcus Agalactiae.

The same precision and convenience of Accumast®, plus the capability of distinguishing S. Agalactiae from other Streptococcus.

LABORATORY SET UP

On-Farm Laboratory Area.

A counter or horizontal work surface that can be easily disinfected and kept free of clutter is needed to set up your incubator and culture plating space in a dedicated area for your on-farm laboratory. Avoid places that exhibit large fluctuations in room temperature and are not free of drafts. This area should always be kept clean in order to minimize the risk of contamination of the culture plates, which may result in false positive or influence the accuracy of your results.

Workers should always wear new, clean disposable gloves when working with mastitis samples or culture plates and wash their hands when the work is completed.

Incubator Setup and Operation.

The environment inside an incubator is controlled for warmth and moisture.

The goal is to provide a controlled environment that allows bacteria to grow and multiply, therefore your incubator must maintain its temperature at 37°C (98.6°F). Please keep an easy-to-read thermometer inside your incubator to allow you to monitor that your incubator is working properly. It is critical to verify if the door of the incubator is well closed and latched at all times.

The humidity inside the incubator should be about 50%. As a rule of thumb, keep 1/2 cup of water in the middle reservoir in the bottom of the incubator. Add water as needed to keep level 1/2 full. Add warm tap water. Without the proper temperature and humidity in your incubator, disease causing bacteria may fail to grow in your samples, leading to a false negative result of "No growth". Temperature and water level in the incubator should be checked daily.

Note: if you are not culturing plates everyday do not turn the incubator off. Leave the incubator on at all times so the temperature inside the incubator maintains fairly constant.

RECORD KEEPING:

Complete and accurate records are an essential component for any laboratory. We recommend a written log book be kept at the incubator. At a minimum you should record the Sample Date, Cow ID, Affected Quarter, and the final Culture Result. You may wish to capture additional information that will help you with quality control in your on-farm laboratory.

Culture results may also be entered in your computer records system (such as DairyComp 305) to help you monitor what organisms are causing mastitis on your dairy, and analyze treatment success or failure. Your veterinarian can help you design a consistent scheme for recording culture results.

SUPPLIES:

Disposable gloves. Sterile 2-ounces leakproof tube or vial for general collection and transportation of samples. Milk sample vials are sterile until opened, and should remain closed until the milk sample is collected.

Alcohol pads or cotton ball soaked in 70% alcohol. Racks for holding sample vials. Waterproof marking pen.

Disinfectant for cleaning teats. Paper or clean cloth towels. Surface disinfectant. Sterile disposable cotton tipped swabs. Sterile disposable cotton tipped swabs should be in a clean dry place. Their packing should not be opened until you are ready to plate a sample. Open the packaging at the wooden stick end so that when you remove a swab from the packing you do not touch the cotton-tipped end. Unused swabs should be kept in the original packaging.

Unused media plates should be stored upside down (lid side down) in the refrigerator. Media should not be used if the expiration date has passed. Do not freeze the media plates.

Used plates, swabs, milk samples and vials should be disposed of properly. Please watch our Demo video (Diagnosis of bovine mastitis- from laboratory to farm) at feraah.com for more details.


INSTRUCTIONS:

The unique and proprietary media formulated and manufacturing process are one of a kind. The media is formulated for rapid growth (fastest in the market) and also for high accuracy and easy interpretability.

Our proprietary chromogenic system is linked to both biochemical tests and to markers specific to the footprint of each organism to provide a multi-layered, high accuracy system. Accumast plus takes the subjectivity as well as the time out of interpreting microbiology, making it possible for the first time to bring these diagnostics back in-house without sacrificing accuracy or simplicity.

ACCUMAST PLUS MEDIA

The Accumast plus tests use a four-section plate system, enabling veterinarians, farmers and herd managers to accurately diagnose almost all commonly encountered infections faster.

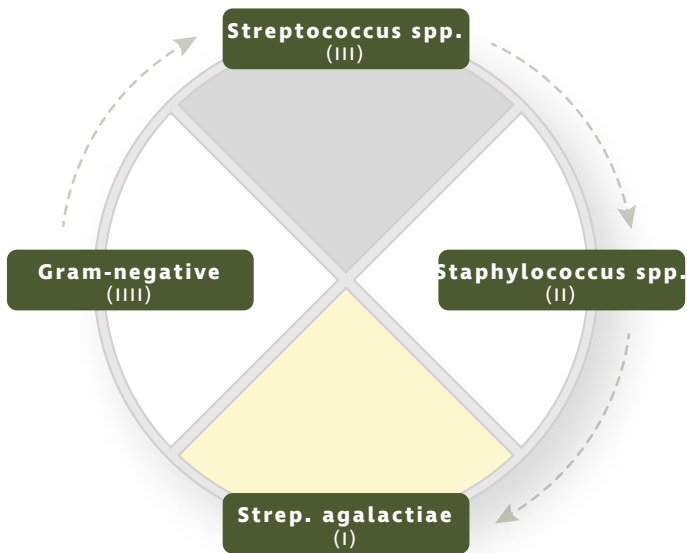

Gram-negative:

E.coli; Klebsiella; Pseudomonas

Streptococcus spp:

Streptococcus spp. Enterococcus spp.
Lactococcus spp.

Staphylococcus spp:

Coagulase negative Staphylococcus;
Staphylococcus aureus

Streptococcus plus:

Streptococcus Agalactiae

STORAGE:

Accumast plus may be shipped and stored for brief periods at room temperature.

Upon receipt store at 2-8°C. Away from direct light. Media should not be used if there are any signs of deterioration (shrinking, cracking, or discoloration), contamination, or if the expiration date has passed. Product is light and temperature sensitive; protect from light, excessive heat, moisture and freezing.

It is recommended Accumast plus to be refrigerated once the pack is open. Expiration is 6 months from the date of manufacturer. It can last up to 8 months if kept protected from light, contamination and dehydration

MAKING THE CULTURE:

Dip a swab cotton tipped swab or loop into your milk sample. If needed, tap excess moisture/debris off on the side of your vial. Use a side to side motion to lightly streak onto the surface of first section of the plate, covering the entire media. Re-dip the swab into the sample and repeat for each of the next sections, dipping into your sample between each. Label the plates with the cow info on the bottom and place upside down (on lid) into incubator at 37°C.

All bacterial organisms will grow within 16 hours, although in many cases the causative organism can be identified in as little as 8-10 hours.

HOW TO READ

When reading the plate, always place the white media at the top and the yellow media at the bottom. Hold the plate facing the surface of the agar – lid side up.

IMPORTANT:


If you cannot distinguish the solid medium by its color, remember that the Gram-negative medium is in the compartment III (written on the back of the plate), and the Staphylococcus medium is in compartment II.

A dark-colored cow with a yellow ear tag is shown in profile, looking down at the ground. The background is a blurred barn interior with other cows. The overall lighting is dim and moody.

INTERPRETING



THE RESULTS

GRAM-NEGATIVE

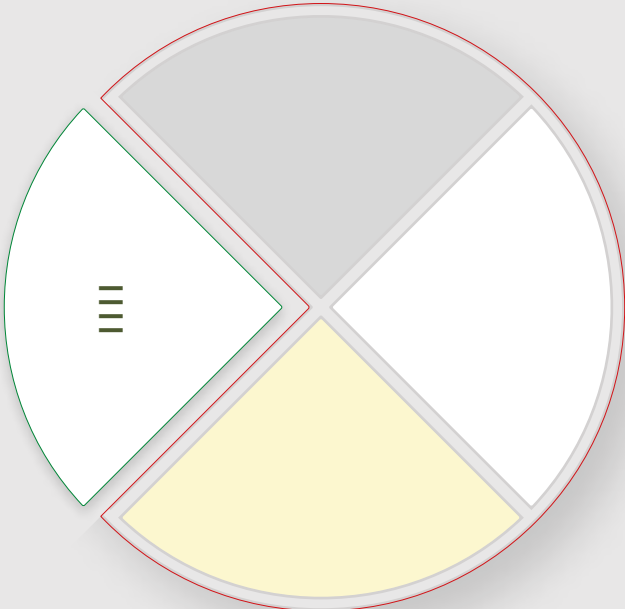



E. COLI

Pink

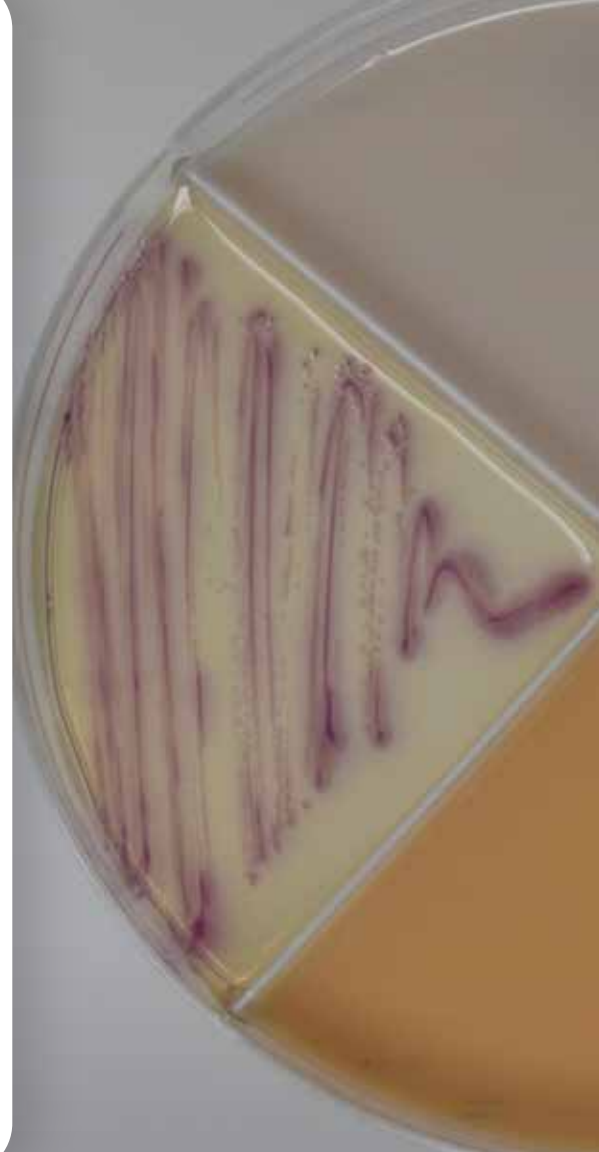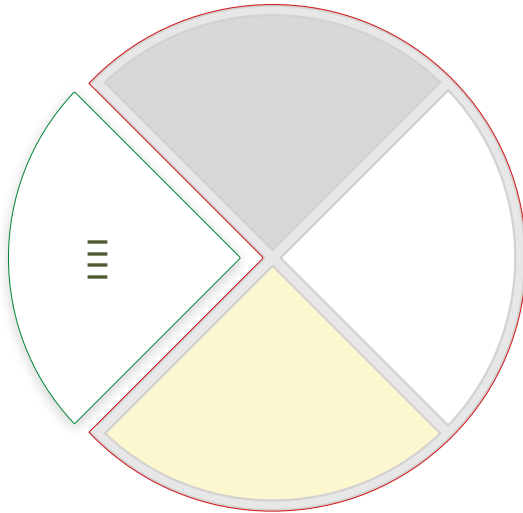


PSEUDOMONAS AERUGINOSA

White colonies with yellow media

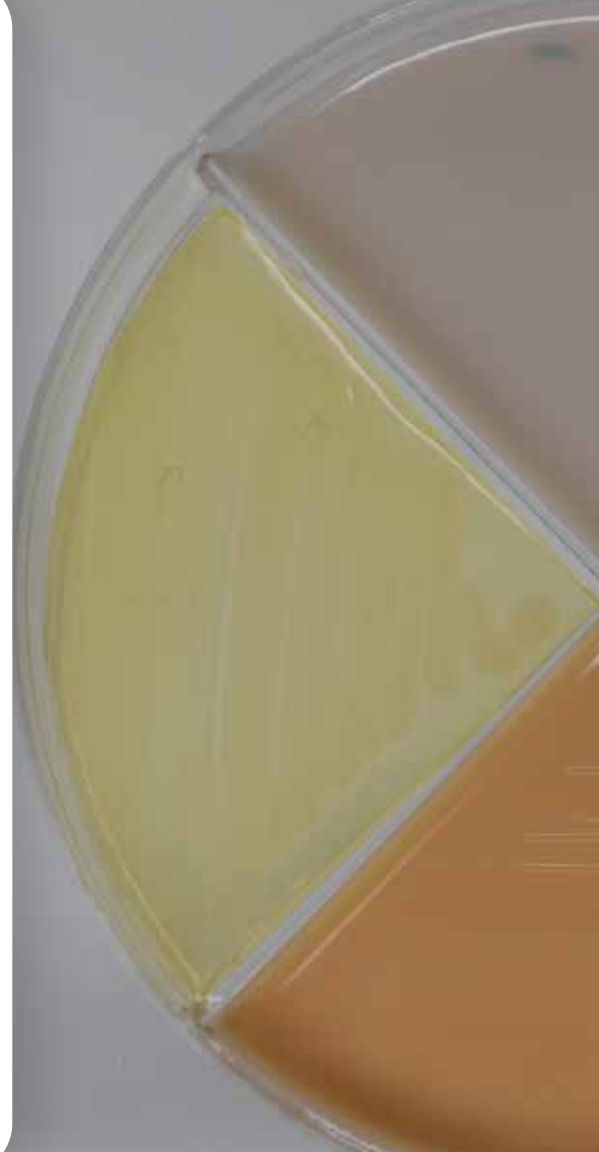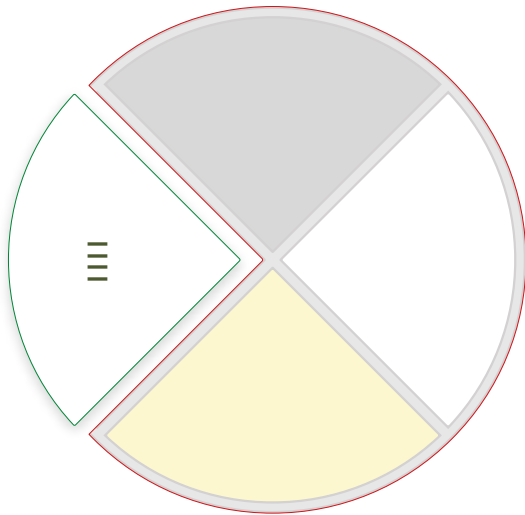


KLEBSIELLA SPP.

Blue


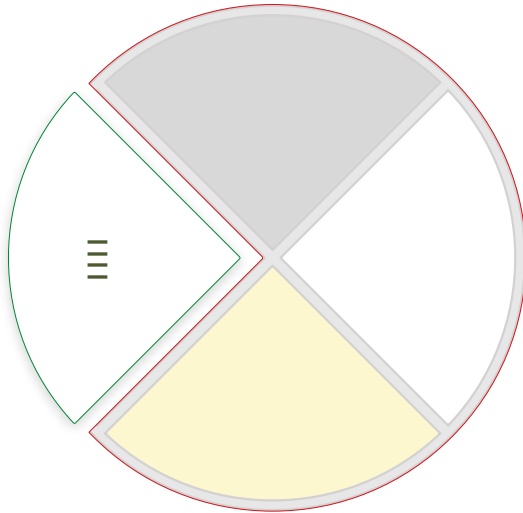


STREPTOCOCCUS SPP.

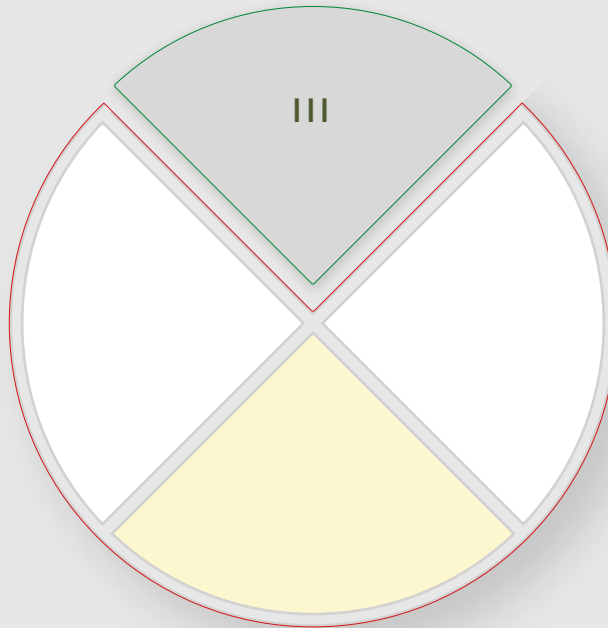


STREPTOCOCCUS AGALACTIAE & DYSGALACTIAE

Light blue

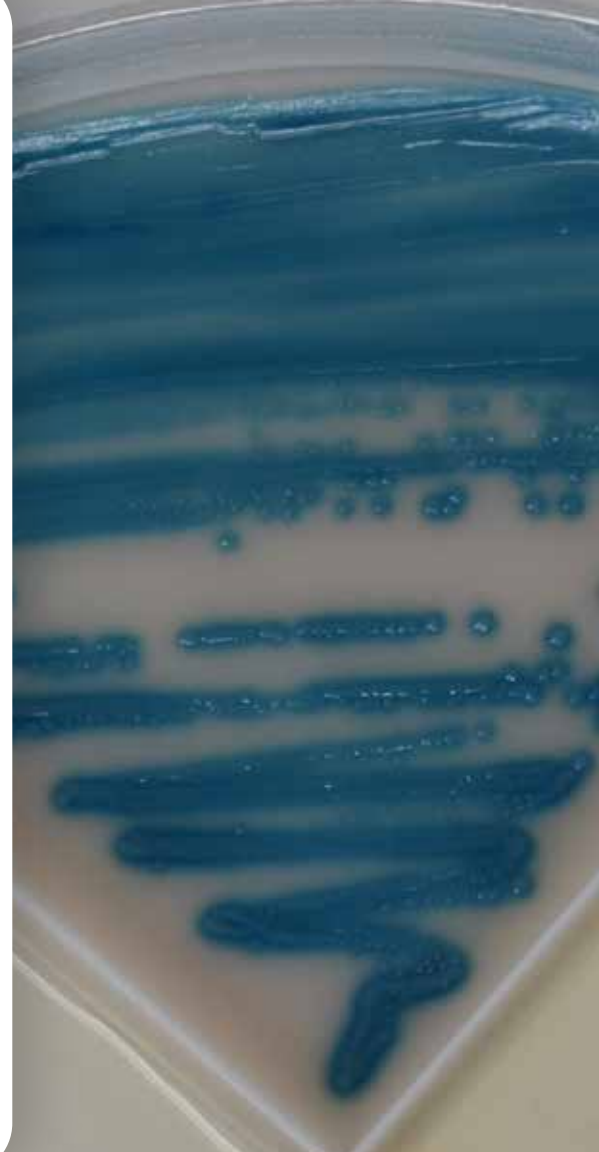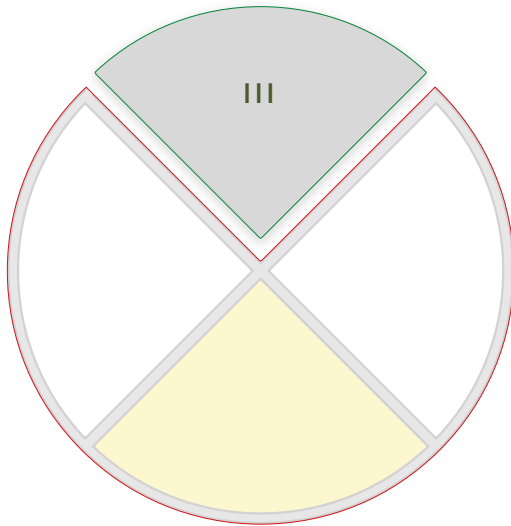


STREPTOCOCCUS UBERIS

Dark blue

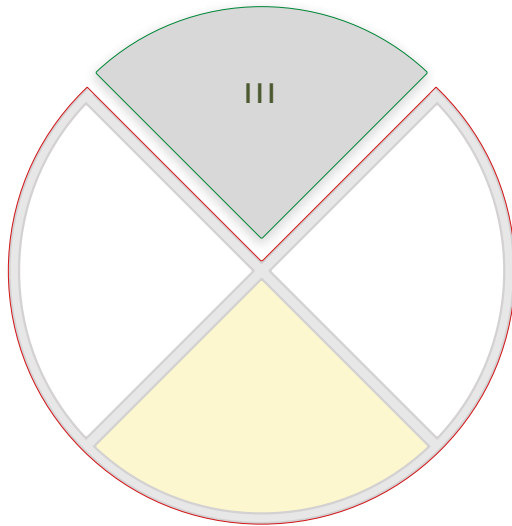



ENTEROCOCCUS SPP.

Large purple colonies

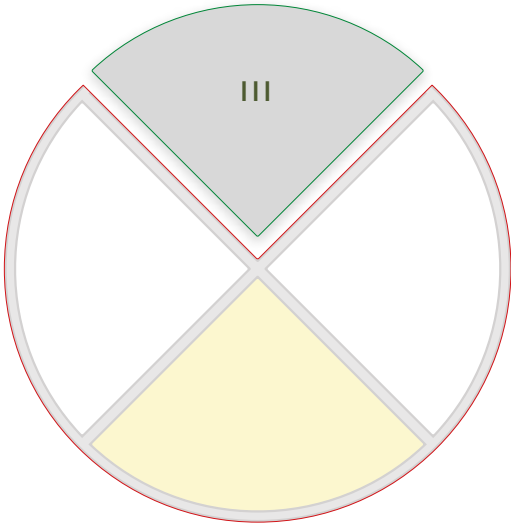


LACTOCOCCUS SPP.

Small lavender or White colonies

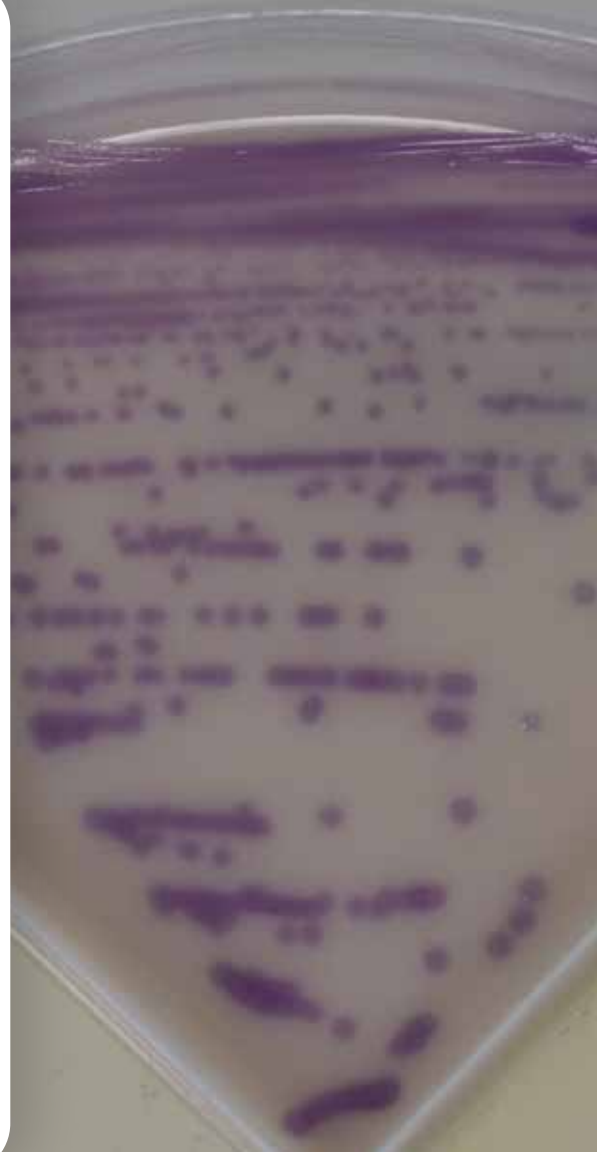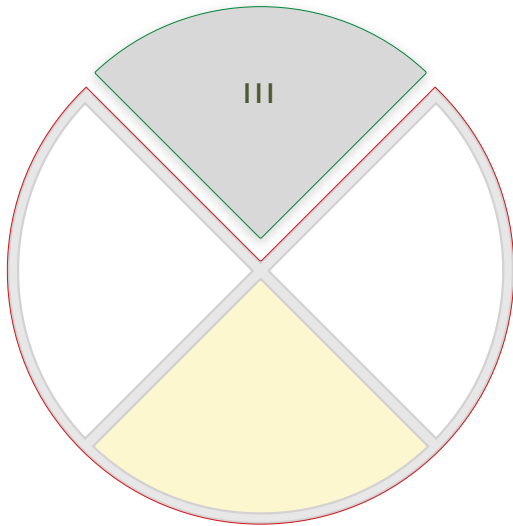


STAPHYLOCOCCUS SPP.

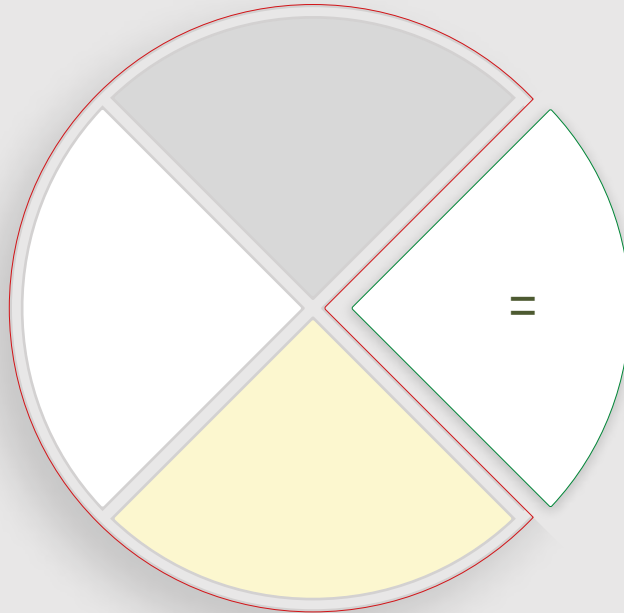


STAPHYLOCOCCUS AUREUS

Pink

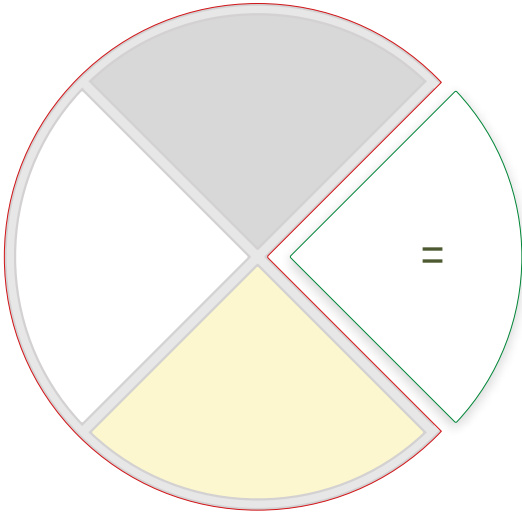


STAPHYLOCOCCUS CHROMOGENES

White

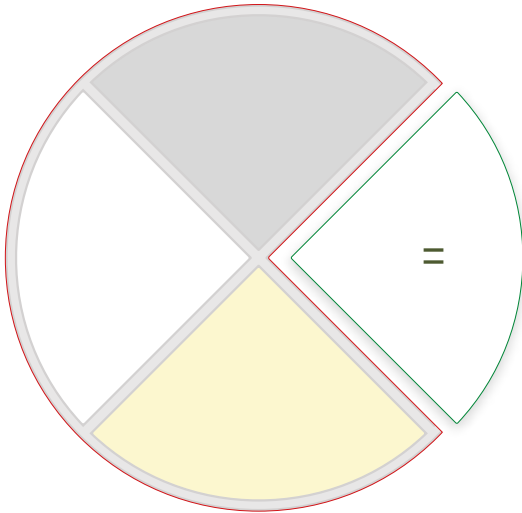


STAPHYLOCOCCUS HAEMOLYTICUS

Blue or green colonies

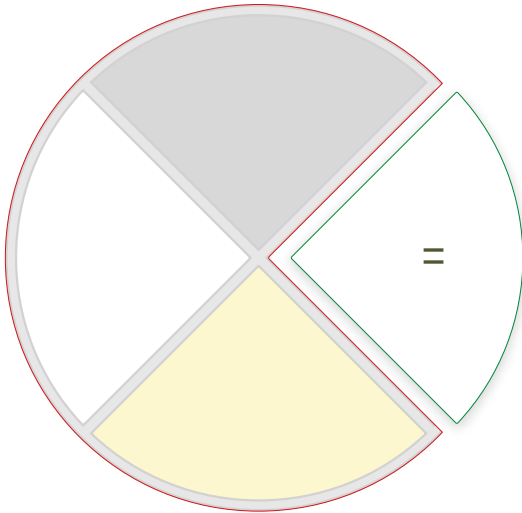


STREPTOCOCCUS PLUS

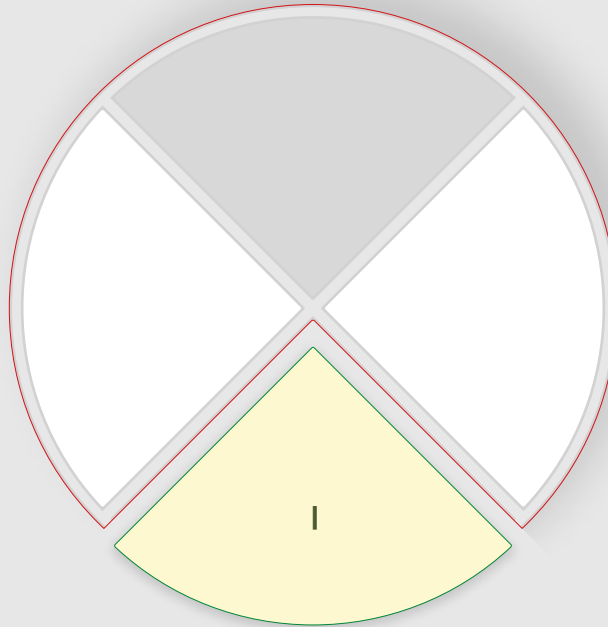


STREPTOCOCCUS AGALACTIAE

Blue

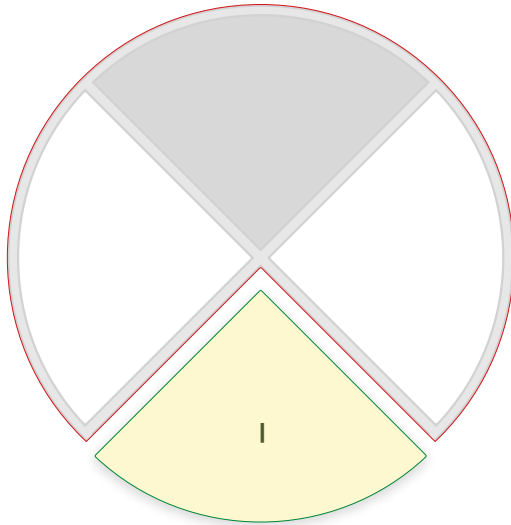


**STREPTOCOCCUS
AGALACTIAE**

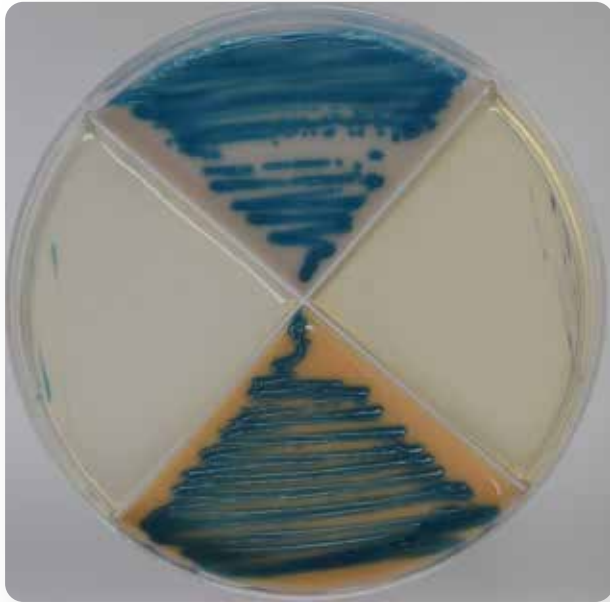


ENTEROCOCCUS & LACTOCOCCUS SPP.

Purple

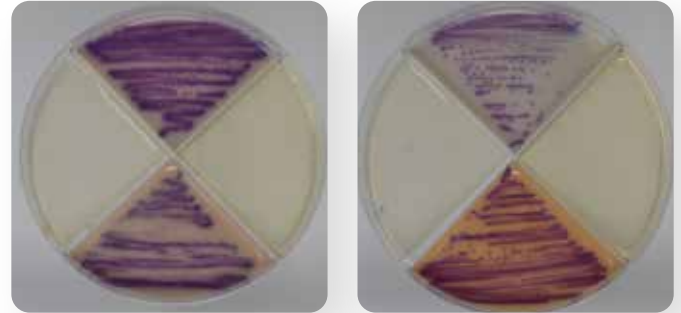

STREP. DYSGALACTIAE & UBERIS

No growth

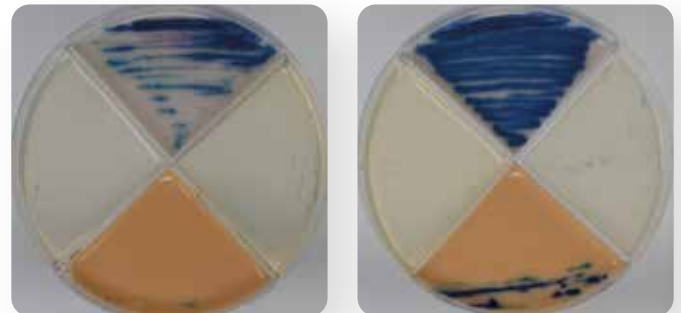



FERA

Diagnostics and Biologicals

COLOR MATCHING AND
BACTERIAL IDENTIFICATION

 **ACCUMAST^{Plus}**
The Future of Culture



A photograph of a cow in a field, overlaid with a semi-transparent yellow filter. The cow is facing left, and the background shows a field of flowers under a bright sky.


CULTURE
NEGATIVE



Culture negative is characterized by the absence of bacterial growth in the selective plates. A No Growth result may indicate that the bacterial is no longer present, most cases are due to spontaneous cure of the infection and antibiotic therapy is not likely to be beneficial.

Some plates may result in no bacterial growth, even if a cow has obvious signs of mastitis. This occurs for several reasons. We recommend to consult your veterinarian and investigate for the presence of Anaerobic bacteria, *Mycoplasma bovis*, and other unknown pathogens. Low bacteria load at the time that the sample is taken is also an obstacle to overcome, therefore we recommend to continuous sampling to catch the pathogen when it is in a shedding state.


Keep in mind that the clots, flakes, and abnormal milk secretions are a result of udder damage from the bacterial assault, not from the bacteria themselves. The artifact caused by blurring and milk clots are not bacterial colonies and are considered a negative result (i.e. plate showing growth of *Enterobacter* colonies on Gram negative media; and presence of milk clots on *Streptococcus* spp.. media)

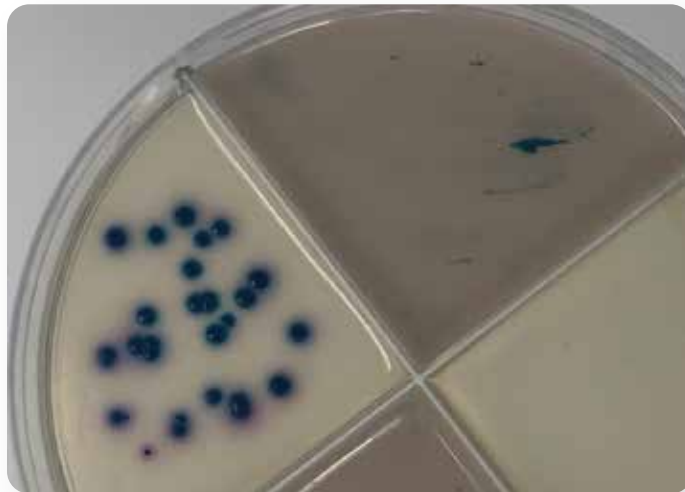

CONTAMINATED SAMPLE

In general, growth in any section of the plate indicates bacteria are present in the milk and antibiotic therapy may be indicated.

This may not be true if more than two different pathogens are present in a sample.

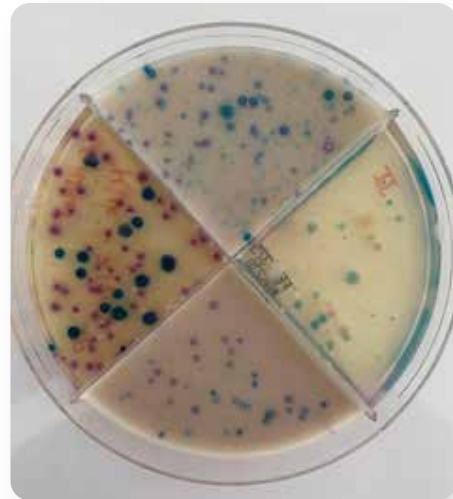

A dark horse with a yellow ear tag is shown in profile, eating from a wooden trough. The background is a blurred stable with other horses. The text is overlaid on the right side of the image.

TREATMENT

ACCORDING TO
EACH DIAGNOSTIC

HOW TO TREAT MASTITIS CAUSED BY E. COLI?

Severe cases must be treated!

Systemic antibiotics treatment is recommended to reduce the risk of bacteremia

Supportive therapy

- 2 L of hypertonic saline (7.2%)
- 500 mL of Calcium borogluconate
- 500 mL of dextrose 50%
- Anti-inflammatory

HOW TO TREAT MASTITIS CAUSED BY KLEBSIELLA?

Severe cases must be treated with systemic broad-spectrum antibiotics and supportive therapy

- 2 L of hypertonic saline (7.2%)
- 500 mL of calcium borogluconate
- 500 mL of dextrose 50%
- Anti-inflammatory

Extended IM therapy with broad-spectrum antibiotics

Probability of cure of ~60-70%

Milk production is normally impaired throughout the productive life

PSEUDOMONAS SPP.

- **Sporadic occurrence**
- **It can occur in cases of outbreaks where a single source of contamination is detected**
 - Cloth towels used in the milking routine
- **Infections become chronic and do NOT respond to antimicrobials**
- **Permanent dry off the chronic mastitic quarters or culling**
- **Severe cases must be treated with systemic broad-spectrum antibiotics and fluid therapy**
 - 2 L of hypertonic saline (7.2%)
 - 500 mL of calcium borogluconate
 - 500 mL of dextrose 50%
 - Anti-inflammatory

OTHER GRAM-NEGATIVES PATHOGENS

Pasteurella, Pseudomonas, Serratia and Enterobacter

- Supportive therapy
- Do not use intramammary antibiotics

Severe cases must be treated with systemic broad-spectrum antibiotics and fluid therapy

- 2 L of hypertonic saline (7.2%)
- 500 mL of calcium borogluconate
- 500 mL of dextrose 50%
- Anti-inflammatory

HOW TO TREAT MASTITIS CAUSED BY STREP. UBERIS?

Severe cases must be treated with systemic broad-spectrum antibiotics and fluid therapy

- 2 L of hypertonic saline (7.2%)
- 500 mL of calcium borogluconate
- 500 mL of dextrose 50%
- Anti-inflammatory

Extended intra-mammary antibiotic therapy (5-8 days)

Systemic Antibiotics are NOT necessary

HOW TO TREAT MASTITIS CAUSED BY STREP. AGALACTIAE?

USE OF INTRA-MAMMARY ANTIBIOTIC PROTOCOLS OF SHORT DURATION

5% of the cows are chronic cases → CULL

HOW TO TREAT MASTITIS CAUSED BY OTHER STREPTOCOCCUS SPP.?

**USE OF INTRA-MAMMARY ANTIBIOTIC
PROTOCOLS OF SHORT DURATION**

HOW TO TREAT MASTITIS CAUSED BY ENTEROCOCCUS?

**USE OF INTRA-MAMMARY ANTIBIOTIC
PROTOCOLS OF SHORT DURATION**

HOW TO TREAT MASTITIS CAUSED BY LACTOCOCCUS SPP.?

Extended intra-mammary antibiotic (5-8 days)

Systemic antibiotics are NOT necessary

HOW TO TREAT MASTITIS CAUSED BY STAPHYLOCOCCUS NON-AUREUS?

**USE OF INTRA-MAMMARY ANTIBIOTIC
PROTOCOLS OF SHORT DURATION**

HOW TO TREAT MASTITIS CAUSED BY STAPHYLOCOCCUS AUREUS?

Primiparous cows → USE OF INTRA-MAMMARY ANTIBIOTIC (5 – 8 days)

Multiparous cows → CULL

A dark-colored cow with a yellow ear tag is shown in profile, looking down at the ground. The background is a blurred barn interior with other cows. The text "DISCOVER OUR OTHER PRODUCTS" is overlaid in white, bold, sans-serif font, with a thin yellow horizontal line underneath the second line of text.


DISCOVER OUR
OTHER PRODUCTS

ACCUTREAT

FOR RAPID DETECTION OF MASTITIS PATHOGENS
THAT REQUIRE ANTIBIOTIC THERAPY.

**Entire herd
management for identification
of Gram-positive pathogens**


BLUE = *Streptococcus* spp.; *Lactococcus* spp. & *Enterococcus* spp.

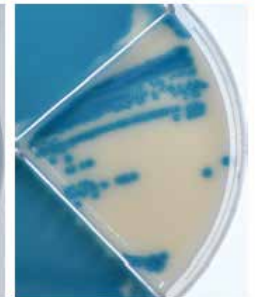

PINK: *Staph.
aureus*

ORANGE: *Staph.
chromogenes*

PURPLE: *Staph.
haemolyticus*

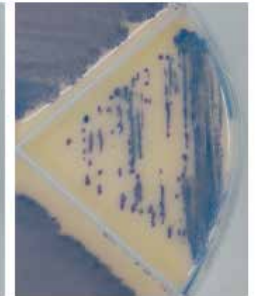

ACCUSTREP

FOR RAPID, PRECISE AND CONVENIENT IDENTIFICATION
OF STREPTOCOCCUS AGALACTIAE.

**Entire herd
control management for
Streptococcus agalactiae**



ACCUSTAPH®

DESIGNED TO KEEP STAPH AUREUS AT BAY


FOR USE WITH THE ACCUMAST® MASTITIS CULTURE SYSTEM;
SCREEN FOUR COWS FOR STAPHYLOCOCCUS AUREUS AT ONCE.

**Entire herd
control management for
Staphylococcus aureus**

PINK: Staphylococcus aureus




*BLUE: Staphylococcus
equorum*

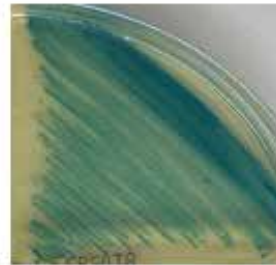

*YELLOW: Staphylococcus
warneri*



*GREEN: Staphylococcus
haemolyticus*

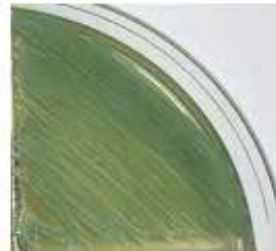

*WHITE: Staphylococcus
chromogenes*



RESOURCES

Cebra C.K., Garry F., Dinsmore R. Naturally occurring acute coliform mastitis in Holstein cattle J. Vet. Intern. Med., 10 (1996), pp. 252-257

Erskine RJ. Antibacterial therapy of clinical mastitis - part I. Drug selection. Part II Administration. North Am Vet Conf, Proc. 2003. pp. 13–16.

Ferreira JC, Gomes MS, Bonsaglia ECR, Canisso IF, Garrett EF, Stewart JL, et al. (2018) Comparative analysis of four commercial on-farm culture methods to identify bacteria associated with clinical mastitis in dairy cattle. PLoS ONE 13(3): e0194211. <https://doi.org/10.1371/journal.pone.0194211>

Fuenzalida MJ, Ruegg PL. Negatively controlled, randomized clinical trial to evaluate intramammary treatment of nonsevere, gram-negative clinical mastitis. J Dairy Sci. 2019;102(6):5438-5457. doi:10.3168/jds.2018-16156

Ganda EK, Bisinotto RS, Decter DH, Bicalho RC. Evaluation of an On-Farm Culture System (Accumast) for Fast Identification of Milk Pathogens Associated with Clinical Mastitis in Dairy Cows. PLoS One. 2016;11(5): e0155314. Published 2016 May 13. doi: 10.1371/journal.pone.0155314

Ganda, E. K. et al. Longitudinal metagenomic profiling of bovine milk to assess the impact of intramammary treatment using a third-generation cephalosporin. Sci. Rep. 6, 37565; doi: 10.1038/ srep37565 (2016).

Gröhn Y.T., Wilson D.J., González R.N., Hertl J.A., Bennett G., Schukken Y.H. Effect of pathogen-specific clinical mastitis on milk yield in dairy cows. J. Dairy Sci. 2004; 87 (15377615): 3358-3374. 10.3168/jds. S0022-0302(04)73472-4

Lago A., Godden S.M., Bey R., Ruegg P.L., Leslie K. The selective treatment of clinical mastitis based on on-farm culture results: I. Effects on antibiotic use, milk withholding time, and short-term clinical and bacteriological outcomes. *J. Dairy Sci.* 2011; 94 (21854917): 4441-4456. [10.3168/jds.2010-4046](https://doi.org/10.3168/jds.2010-4046)

Lago A., Godden S.M., Bey R., Ruegg P.L., Leslie K. The selective treatment of clinical mastitis based on on-farm culture results: II. Effects on lactation performance, including clinical mastitis recurrence, somatic cell count, milk production, and cow survival. *J. Dairy Sci.* 2011; 94 (21854918): 4457-4467. [10.3168/jds.2010-4047](https://doi.org/10.3168/jds.2010-4047)

Lago A, Godden SM. Use of Rapid Culture Systems to Guide Clinical Mastitis Treatment Decisions. *Vet Clin North Am Food Anim Pract.* 2018;34(3):389-412. doi: [10.1016/j.cvfa.2018.06.001](https://doi.org/10.1016/j.cvfa.2018.06.001)

Lima SF, Bicalho MLS, Bicalho RC (2018) Evaluation of milk sample fractions for characterization of milk microbiota from healthy and clinical mastitis cows. PLoS ONE 13(3): e0193671. <https://doi.org/10.1371/journal.pone.0193671>

Oikonomou G., Bicalho M.L., Meira E., Rossi R.E., Foditsch C., Machado V.S. Teixeira A.G., Santisteban C., Schukken Y.H., Bicalho R.C. Microbiota of cow's milk; distinguishing healthy, sub-clinically and clinically diseased quarters. PLoS One. 2014; 9 (24465777): e85904. [10.1371/journal.pone.0085904](https://doi.org/10.1371/journal.pone.0085904) Roberson J.R.

Oliver SP, Almeida RA, Gillespie BE. et al. Extended ceftiofur therapy for treatment of experimentally-induced *Streptococcus uberis* mastitis in lactating dairy cattle. J Dairy Sci. 2004; 87:3322–3329. doi: 10.3168/jds. S0022-0302(04)73468-2.

Pyörälä S, Pyörälä E. Efficacy of parenteral administration of three antimicrobial agents in treatment of clinical mastitis in lactating cows: 487 cases: (1989-1995) J Am Vet Med Assoc. 1998;212:407–412.

Rodrigues M.X., Lima S.F., Higgins C.H., Canniatti-Brazaca S.G., Bicalho R.C. The Lactococcus genus as a potential emerging mastitis pathogen group: a report on an outbreak investigation. J. Dairy Sci., 99 (2016), pp. 9864-9874, 10.3168/jds.2016-11143

Tomazi T, Freu G, Alves BG, de Souza Filho AF, Heinemann MB, Veiga dos Santos M (2019) Genotyping and antimicrobial resistance of Streptococcus uberis isolated from bovine clinical mastitis. PLoS ONE 14(10): e0223719.

Vasquez A.K., Nydam D.V., Capel M.B., Eicker S., Virkler P.D. Clinical outcome comparison of immediate blanket treatment versus a delayed pathogen-based treatment protocol for clinical mastitis in a New York dairy herd. J. Dairy Sci. 2017; 100 (28161180): 2992-3003. 10.3168/jds.2016-11614

Vasquez A.K., Ganda E.K., Capel M.B., Virkler P.D., Bicalho R.C., Nydam D.V. The microbiome of Escherichia coli and culture-negative non-severe clinical mastitis: Characterization and associations with linear score and milk production. J. Dairy Sci. 2018; 102:578-594

Warnick L.D., Moore G. Mild to moderate clinical mastitis: Efficacy of intramammary amoxicillin, frequent milk-out, a combined intramammary amoxicillin, and frequent milk-out treatment versus no treatment. *J. Dairy Sci.* 2004; 87 (15202642): 583-592

Wenz J.R., Barrington G., Garry F., McSweeney K., Dinsmore R., Goodell G., Callan R. Bacteremia associated with naturally occurring acute coliform mastitis in dairy cows. *J. Am. Vet. Med. Assoc.*, 219 (2001), pp. 976-981



FERA

Diagnostics and Biologicals

www.feraah.com