

Mepilex® Border Ag

THE EFFECTIVE ANTIMICROBIAL ABSORBENT ALL-IN-ONE FOAM DRESSING

- Minimises pain and trauma at dressing change^{1,2,3}
- Rapid and sustained antimicrobial activity⁵
- Excellent and effective exudate management
- Safe fixation but still atraumatic removal

Polyurethane backing film

- Breathable
- Water proof
- Bacteria and viral barrier*

Safetac® layer

- Reduces pain and trauma before, at and after dressing change¹
- Does not adhere to the moist wound bed but to dry skin only³
- Seals the wound margins and reduces risk of maceration⁴

Multi-layered absorbent foam pad

- The compressed foam with silver, absorbs and spreads the exudate while allowing silver to be released.
- Fibre layer for maximum retention

Rapid, sustained release of ionic silver

- Antimicrobial activity within 30 minutes up to 7 days⁵


* microbes larger than 25 nm

Safetac® technology


Safetac® technology. Less pain and less trauma.

Safetac is a proprietary adhesive technology that minimises pain to patients and trauma to wounds. Safetac technology is available exclusively on Mölnlycke Health Care dressings, including Mepilex® wound care dressings, Mepitel®, Mepiform® and Mepitac®.

For more information visit www.safetac.com

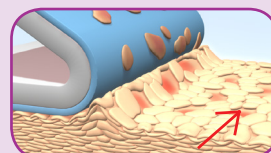


Painful skin stripping with traditional adhesive

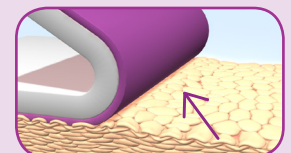


No skin stripping with Safetac technology

How to use Mepilex® Border Ag

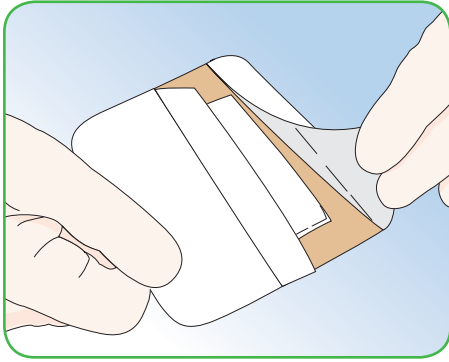

Clean the wound area.
Remove the release film.

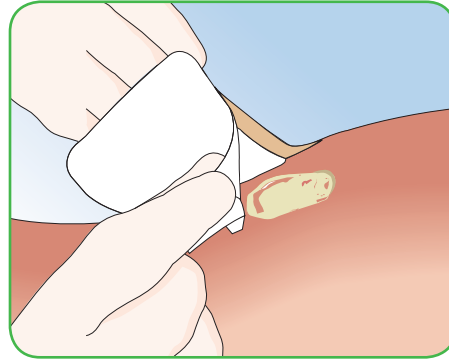

Apply the adherent side to the wound.
Do not stretch.


How Mepilex Border Ag works

Mepilex Border Ag is an antimicrobial all-in-one foam dressing that effectively absorbs and retains exudate and maintains a moist wound environment. The Safetac layer seals the wound edges, preventing the exudate to leak onto the surrounding skin, thus minimising the risk of maceration. The Safetac layer ensures that the dressing can be changed without damaging the wound or surrounding skin or exposing the patient to additional pain.

Mepilex Border Ag inactivates wound related pathogens within 30 minutes with sustained effect up to 7 days. By reducing the number of microorganisms, Mepilex Border Ag may also reduce odour.

Benefits of Mepilex Border Ag

- Minimises pain and trauma at dressing changes
- Self-adherent - no secondary fixation needed
- Shower proof
- Can be used under compression bandages
- Promotes patient comfort during wear
- Can remain in place for several days depending on the condition of the wound
- Can be lifted and adjusted without losing its adherent properties
- Low potential for skin irritation and allergy

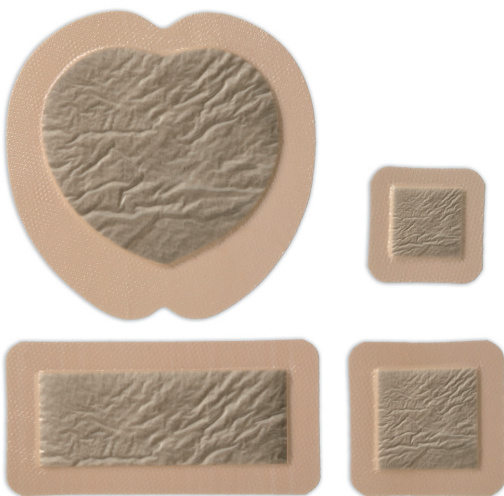

Areas of use

Mepilex Border Ag is designed for medium to high exuding wounds such as leg and foot ulcers, pressure ulcers, malignant wounds, partial thickness burns, traumatic and surgical wounds where a moist environment, exudate management, gentle fixation and an anti-microbial action is indicated.

For detailed information please see Instructions for use.



MEPILEX BORDER AG ASSORTMENT (Sterile packed)

Code	Size cm	Pcs/shelf cont.	Pcs/transp. cont.
395200	7.5 x 7.5	5	70
395300	10 x 10	5	50
395400	15 x 15	5	50
395600	15 x 20	5	45
395800	10 x 20	5	35
395700	10 x 25	5	35
395900	10 x 30	5	25
Mepilex Border Sacrum Ag			
382000	18 x 18	5	40
382400	23 x 23	5	25

Mölnlycke Health Care Pty Ltd
Suite 1.01, 10 Tilley Lane, Frenchs Forest NSW 2086
Phone 1800 005 231
www.molnlycke.com.au
New Zealand Orders & Enquiries 0800 005 231, www.molnlycke.co.nz

