

Instructions for Use



Document Ref: OPQC/IFU/014/3	Version: 3.1	Page: 1
Raised by (Name): Grant Porter Signature:	Position: Operations & Logistics Manager	Date: 01/04/2022
Authorised by (Name): Rebecca Porter Signature:	Position: Managing Director	Date: 01/04/2022

SINGLE-USE SURGICAL INSTRUMENTS: Surgical Forceps

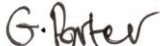

Doc No.	OPQC/IFU/014/3	Rev.	3.1	Date of Issue	20/08/22 Last doc review date: 19/12/23	Condition:	SINGLE-USE 	Condition for Use:	STERILE
Device		Single use surgical soft-tissue manipulation forceps, alligator Single use surgical soft-tissue manipulation forceps, scissors-like Single use ophthalmic soft-tissue manipulation forceps, tweezers-like Single use obstetrical forceps Single use dressing/utility forceps, tweezers-like							
Intended Use	For professional use. Intended for single use only. Surgical Forceps are hand-held, hinged instruments used for grasping and holding objects and tissues. Needle Holder Forceps are hand-held, hinged instruments used for grasping and holding needles and sutures for stitching during surgical procedures. Thumb, Dressing & Tissue Forceps are a handheld, hinged or tweezer-like instruments used for grasping and holding objects and tissues.								
Sterilisation (R)		Method:	Gamma Radiation Processed under ISO 13485: 2016 procedures		Radiation Dose Range:	25.00 Min – 35.00 Max kGy			
		Instructions for Use:	DO NOT USE if pack is damaged or open.						
Storage	<ol style="list-style-type: none"> 1. Store all instruments as packaged. 2. Do not use if packaging is damaged or open. 3. Ensure all packs are kept in a storage area free of dust, insects, chemical vapours and extreme changes in temperature and humidity. 								
Warnings/Cautions		<p>Surgical instruments may be damaged if used for unintended purposes.</p> <p>NEVER:</p> <ul style="list-style-type: none"> ▪ Misuse instruments or overstrain joints or racks. ▪ Handle instruments by their tips. These will have been cleaned by trained personnel only who will ensure the delicate working ends are adequately protected during storage or sterilisation. ▪ Attempt to clean or use abrasives on instruments as this will spoil the surface finish and compromise sterility. ▪ Store instruments with other heavier instruments on top. Always remember to store heavy instruments on the bottom and light instruments on the top. ▪ Reuse a single use device as this may result in cross infection or functional failure. 							





BlueKit Medical Limited

T 020 8224 2883 | E info@bluekitmedical.com | W www.bluekitmedical.com
 A Unit 87 Capital Business Centre | 22 Carlton Road | South Croydon | CR2 0BS
 Company Registered In England 08097913 | VAT Registration No. 141004673

Instructions for Use

BlueKit
MEDICAL

Document Ref: OPQC/IFU/014/3	Version: 3.1	Page: 2
Raised by (Name): Grant Porter Signature: 	Position: Operations & Logistics Manager	Date: 01/04/2022
Authorised by (Name): Rebecca Porter Signature: 	Position: Managing Director	Date: 01/04/2022

Instructions	<ul style="list-style-type: none"> ▪ Always store appropriately, ensuring storage is free of dust, insects, chemical vapours and extreme changes in temperature and humidity. ▪ Always check packaging before use to ensure that packaging is not damaged or open. ▪ Always check instruments for damage before use. Give special attention to microsurgical instruments: their fine tips can easily be damaged in transit or storage. ▪ Ensure all instruments are only used for the purpose for which they were designed. ▪ Handle all instruments gently. Never overstrain, drop, or misuse them. ▪ Where any serious incident occurs as a result of use, this incident must be reported to the manufacturer and competent authority.
Contraindications	<ol style="list-style-type: none"> 1. Use without Sterilisation/ Disinfection. 2. Rusty conditions visually observed. 3. Re-use of Single Use, whether used on a single patient, another patient, or another surgical procedure occurring at another time.
Shelf Life 	Refer to date of manufacture and date of expiration available on product label. Gamma Irradiation (R) provides sterility for 5 years.
Disposal	Single use surgical instruments should be disposed of after use. Used single use surgical instruments should be treated as clinical waste and should therefore be disposed of in accordance with and under the care and direction of the healthcare provider's full Policies and Procedures and must be managed strictly in accordance with local healthcare waste management policies and procedures. If decontaminated stainless steel is recyclable in ordinary steel manufacturing units.
Distributor 	BlueKit Medical Limited https://bluekitmedical.com Unit 87 Capital Business Centre, 22 Carlton Road, South Croydon, CR2 0BS, England, United Kingdom
Packing provided by	Granton Medical Limited Parkway Close, Parkway Industrial Estate, Sheffield, S9 4WJ, United Kingdom
Sterilisation provided by	Swann-Morton Limited Owlerton Green, Sheffield, S6 2BJ, United Kingdom
Regulatory Standards	This product is packed in accordance with requirements set out within ISO 14644-1 and processed under ISO 13485: 2016 procedures. All contents are CE marked.
Manufacturer  	Surgicrafts Muhammad Pura, Commissioner Road, Sialkot, Pakistan CE 1639

BlueKit Medical Limited



T 020 8224 2883 | E info@bluekitmedical.com | W www.bluekitmedical.com




A Unit 87 Capital Business Centre | 22 Carlton Road | South Croydon | CR2 0BS

Company Registered In England 08097913 | VAT Registration No. 141004673

Instructions for Use

BlueKit
MEDICAL

Document Ref: OPQC/IFU/014/3	Version: 3.1	Page: 3
Raised by (Name): Grant Porter Signature: 	Position: Operations & Logistics Manager	Date: 01/04/2022
Authorised by (Name): Rebecca Porter Signature: 	Position: Managing Director	Date: 01/04/2022

EU Representative 	Mr. Manuel Mateos CMC Medical Devices & Drugs S.L, C/ Horacio Lengo n18, C.P 29006, Málaga-Spain
Manufacturer 	S.R.R Surgical Airport Road, Kotli Behram, Sialkot 51310, Pakistan CE 1639
EU Representative 	CMC Medical Devices & Drugs S.L 19-21 Crawford Street, London, W1H 1PJ, United Kingdom

BlueKit Medical Limited

T 020 8224 2883 | E info@bluekitmedical.com | W www.bluekitmedical.com
A Unit 87 Capital Business Centre | 22 Carlton Road | South Croydon | CR2 0BS
Company Registered In England 08097913 | VAT Registration No. 141004673