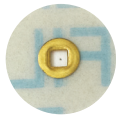
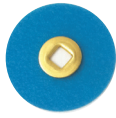


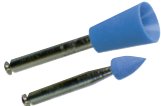

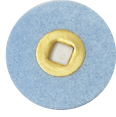



## COMPOSITE POLISHING GUIDE

Cosmedent gives you everything you need to polish composites. With our conveniently color-coded easy to use polishers, you can create a smooth and esthetic polish that mimics enamel long-term.

Cosmedent Polishing Tool	Speed Range	Use Notes
<b>Coarse FlexiDisc®</b> 	3,000 to 15,000 RPM for 3/8" disc 5,000 to 15,000 RPM for all other sizes	<ul style="list-style-type: none"> <li>• Use light pressure</li> <li>• Gross reduction</li> <li>• Initial contouring</li> <li>• Buffing out imperfections</li> </ul> <ul style="list-style-type: none"> <li>• Available Regular center sizes: 1/2" &amp; 5/8"</li> <li>• Available Mini center sizes: 3/8", 1/2", &amp; 5/8"</li> </ul>
<b>Medium FlexiDisc®</b> 	3,000 10,000 RPM around gingival margin with 3/8" disc 10,000 to 15,000 RPM for all other sizes	<ul style="list-style-type: none"> <li>• Completion of finish</li> <li>• Removal of scratches</li> <li>• Use flex and pressure</li> </ul> <ul style="list-style-type: none"> <li>• Available Regular center sizes: 1/2" &amp; 5/8"</li> <li>• Available Mini center sizes: 3/8", 1/2", &amp; 5/8"</li> </ul>
<b>Fine FlexiDisc®</b> 	10,000 to 15,000 RPM for all sizes	<ul style="list-style-type: none"> <li>• Beginning final polish</li> <li>• Use flex and pressure</li> </ul> <ul style="list-style-type: none"> <li>• Available Regular center sizes: 1/2" &amp; 5/8"</li> <li>• Available Mini center sizes: 3/8", 1/2", &amp; 5/8"</li> </ul>
<b>Superfine FlexiDisc®</b> 	10,000 to 15,000 RPM for all sizes	<ul style="list-style-type: none"> <li>• Enhancing final polish</li> <li>• Use flex and pressure</li> </ul> <ul style="list-style-type: none"> <li>• Available Regular center sizes: 1/2" &amp; 5/8"</li> <li>• Available Mini center sizes: 3/8", 1/2", &amp; 5/8"</li> </ul>
<b>Medium FlexiCups® &amp; FlexiPoints®</b> 	2,000 to 10,000 RPM MAX with Light Pressure for all cups and points	<ul style="list-style-type: none"> <li>• Use with water</li> </ul> <ul style="list-style-type: none"> <li>• Achieve labial characterization</li> </ul>
<b>Superfine FlexiCups® &amp; FlexiPoints®</b> 	2,000 to 10,000 RPM MAX with Light Pressure for all cups and points	<ul style="list-style-type: none"> <li>• Use with water</li> </ul> <ul style="list-style-type: none"> <li>• Polish gingival margin and labial surfaces</li> </ul>
<b>FlexiBuff®</b> 	15,000 RPM for all sizes	<ul style="list-style-type: none"> <li>• Add Enamelize® to tooth first and use pressure</li> </ul> <ul style="list-style-type: none"> <li>• Available Regular center sizes: 5/8"</li> <li>• Available Mini center sizes: 3/8", 1/2", &amp; 5/8"</li> </ul>
<b>NANO / Microhybrid™ Diamond Polishing Kit</b> 	10,000 to 15,000 RPM MAX with Light Pressure for all polishers	<ul style="list-style-type: none"> <li>• Use for final contouring, finishing and polishing of nanofill and microhybrid composites</li> <li>• 9 reusable diamond polishers, one point, cup and wheel of each grit</li> </ul> <ul style="list-style-type: none"> <li>• Use Coarse grit for initial contouring</li> <li>• Use Medium grit for finishing</li> <li>• Use Superfine grit to develop a final polish</li> </ul>