
Health Check Basic

Kit for self-collection of a capillary blood sample, which is then tested in the laboratory for parameters that represent, among others, an increased risk of cardiovascular diseases, diabetes, iron deficiency or gout. The test is recommended for people who, due to their lifestyle or other factors, have an increased risk of the above-mentioned diseases. The validity of the test is limited and does not replace medical consultation and specific examination. The test must not be performed by persons with blood clotting abnormalities.

Package contents

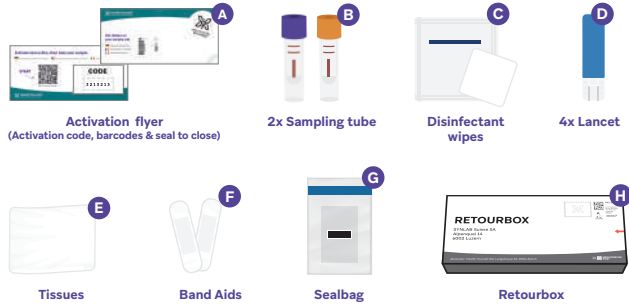Do the test in the afternoon on a weekday, so that you can drop your sample in a Swiss Post mailbox before the evening collection time.

Step 1: Online test activation

Watch the detailed video instruction in the online activation process before you start taking your sample!

Step 2: Take your sample

- 1 Wash your hands and wrists with warm water for about 4 minutes to increase the blood flow. Do NOT rub or clap your hands! This can falsify your results.
2 Take the sampling materials out of the test package and place them on a table in front of you. Unscrew both sample tubes B, attach their lids to the bottom of the tubes and place them on the table.
3 Use a disinfectant wipe C to disinfect the fingertip of the finger you intend to prick with the lancet. We recommend the ring finger of your non-dominant hand, but you can also use any other finger aside from your thumbs and pinky fingers.
4 Remove the white protective pin from the lancet D by turning the pin to the left and then pulling it off (see picture).
5 Stand up and press the nail of the disinfected finger firmly onto the edge of a table (fingertip pointing upwards). With the thumb of the same hand apply pressure below the tip of the ring finger so that the skin of the fingertip tightens. With your free hand, press the lancet against the tightened fingertip (see image for correct placement of the lancet). Press firmly on the lancet until it "clicks".
A small cut is made, which may hurt a little. Stop applying pressure to the fingertip with your thumb and wipe off the first drop of blood with the tissue E.

- 6 Hold the finger with the cut pointing downwards over the test tube. The palm of your hand should continue to face you, so you avoid blood running into your nail. Wait until a drop of blood has formed and then wipe it off on the inner edge of the tube.
7 Ensure continuous blood flow using the "fingertip pump": Continue to hold the bleeding finger over the sample tube. Pinch the fingertip of your bleeding finger between the thumb & index finger of your free hand (see picture). You should be able to apply pressure to the fingertip (skin should tighten). When you apply pressure on the fingertip, a drop of blood should form. Release the pressure, wipe the blood drop off at the inner edge of the tube.
Important: Wait 5-10 seconds to allow more blood to flow into your finger tip. Then apply pressure again.
Not clear? Check out the instruction video in the activation process.

- 8 Repeat the "Finger Tip Pumping Process" until initially the purple tube is filled between the two marked lines. Seal the tube with the cap. You will hear a "click" when lid is properly closed. Immediately turn it upside down 5 times. Continue with the yellow tube: fill it also between the two marked lines, seal it, and turn it upside down 5 times. If the finger stops bleeding before both tubes are filled, repeat the same process on another finger with another lancet. Cover the cut with a band aid F as soon as you are finished. Dispose of lancets in a suitable place, preferably for sharp objects.


- 9 Wipe off any blood on the outside of your closed vials. Take the two barcode stickers from the back of your activation flyer A and stick them onto the sampling vials as shown in the image.
10 Put each sample vial in one of the small pouches made of white absorbing material in the Sealbag G. Close the bag by pulling off the adhesive strip. Write your last name on the barcode label on the Sealbag. Put the Sealbag into the Retourbox H and close the Retourbox with the seal sticker from the back of your activation flyer. Put the Retourbox in a yellow mailbox (no stamp necessary) before the last collection of the day.

In extreme heat or extreme cold, use a protected mailbox or post office for shipping to the laboratory.

Step 3: Laboratory results

Your results will normally be available 2-3 business days after submitting your sample. You can log in to access your results here: https://yourself.health/pages/mytests

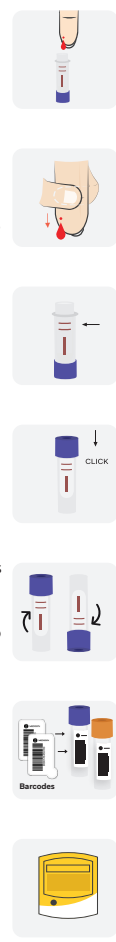This product is intended for single use only. Do not use the test if any content of the package, in particular the lancets or the sampling tubes, is damaged or if the expiry date of the test (see outer package) has passed. Ask for a new test by contacting service@yourself.health. Underfilling or overfilling of the sampling tubes leads to an incorrect ratio of blood to preparation/additive and can lead to incorrect analysis results. In case of direct contact of a third person with the blood sample or a puncture injury of a third person with a used lancet, consult a physician, as this may transmit HIV, HCV, HBV or other infectious diseases. Store and transport at room temperature. Any serious incident related to the device must be reported to the manufacturer and the competent authority of the Member State where the user and/or the patient resides.

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