



"I use DynaBlast® for all my bone grafting procedures. Both the DynaBlast® putty and paste consistently produce predictable clinical results"

– Courtesy: Timothy M. Blieden, DDS, MS, PhD - Rochester, New York

Pre-Clinical Performance - Demonstrated bone formation in a large load-bearing animal model

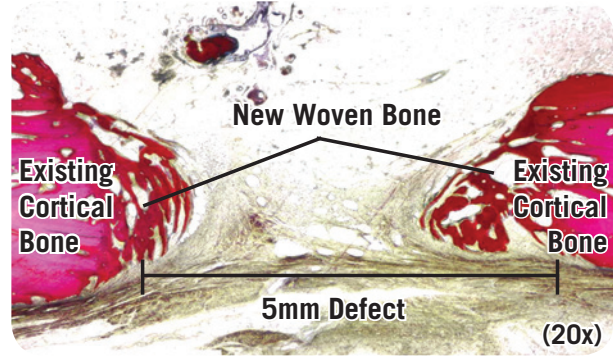

Minimal bone regeneration was observed within the defect. Healing was limited to the area adjacent to the existing cortical bone.

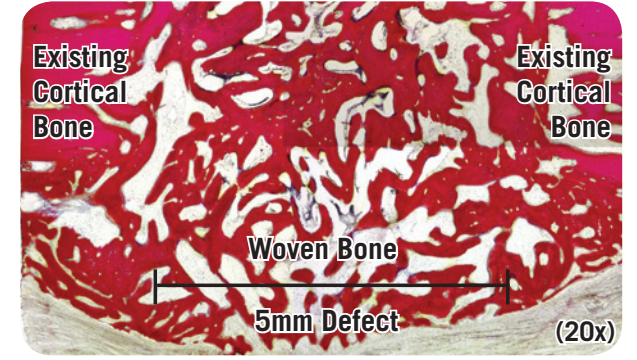

Prolific woven bone was seen bridging the defect at 8 weeks. Active remodeling was evident with no adverse inflammatory response.



16 weeks - Healing of defect demonstrated by the transformation of woven bone to new cortical bone

DynaBlast® Safety

Sterilization procedure for bone graft products
 Keystone Dental bone graft products are produced in HEPA filtered clean rooms with environmentally/microbiologically monitored air. After final packaging, each product is sterilized to an assurance level of 10⁻⁴ using a validated procedure involving the use of electron beam (e-beam) irradiation. By incorporating e-beam sterilization as the last step of the manufacturing process, Keystone Dental provides bone graft products that feature the highest standards of safety—a primary concern of surgeons and patients alike.

Safety of Demineralized Bone Matrix
 Keystone Dental takes many steps to ensure that the demineralized bone matrix (DBM) provided is of the highest quality and safety. This includes obtaining tissue from tissue banks that comply with all tissue banking standards as provided by the Food and Drug Administration (FDA) and the Clinical Laboratory Improvement Amendments (CLIA).

Donor Screening & Testing
 Each lot of DBM is obtained from a single human donor and is not pooled with other donors. Potential tissue donors are rigorously screened by the tissue bank in compliance with FDA requirements, including serological testing performed by a CLIA certified laboratory.

REFERENCES

1. Wientroub S, Reddi AH. Influence of irradiation on the osteoinductive potential of demineralized bone matrix. *Calif Tissue Int* 1988; 42:255-60.
2. Scarborough NL, While EM, Hughes JV, et al. Allograft safety: viral inactivation with bone demineralization. *Contemp Orthop* 1995;31(4):257-61.
3. Mellonig JT, Prewett AB, Moyer MP, et al. HIV inactivation in a bone allograft. *J Periodontal* 1992;63(12):979-83.

HISTOLOGY

"Histological Screening of Graft Materials Using An Ovine Model Phase 1 – 6 treatment 4 and 4A vs. Empty and Autograft Filled Defects"
 PI - Donna L. Wheeler, PhD Department of Mechanical Engineering/Biomedical Engineering Colorado State University
 I (surgical aspect) – A Simon Turner, B.V.Sc., MS. Orthopedic Bioengineering Research Lab, Colorado State University. November 22, 2004



154 Middlesex Turnpike • Burlington, MA 01803
 1.866.902.9272 • www.keystonedental.com
 MK30133 Rev B 11/2017



DynaBlast®

A powerful 2-in-1 solution

DynaBlast® is an allograft material that contains both demineralized and mineralized bone in a proprietary, Reverse Phase Medium carrier. DynaBlast® stimulates new bone growth while providing a natural structural scaffold to support new bone formation.

Dependable Bone Growth

DynaBlast® is a powerful 2-in-1 combination of osteoconductive and osteoinductive elements¹⁻³ to promote new bone growth. Demineralized bone provides active signals to stimulate bone regeneration and mineralized cancellous chips provide a natural structural scaffold to encourage the attachment of osteogenic precursor cells.

Indications for Use


- Filling of defects after tooth resection, apicoectomy and cystectomy
- Filling of extraction sockets to enhance preservation of the alveolar ridge
- Elevation of maxillary sinus floor

Designed to Stay in Place



The proprietary Reverse Phase Medium carrier thickens at body temperature, resists both irrigation and suction and allows for guided clot formation and easy vascularization. DynaBlast® stays in place, with no bone chips migrating from the site—alleviating possible irritation and patient concern.



Implant Dehiscence

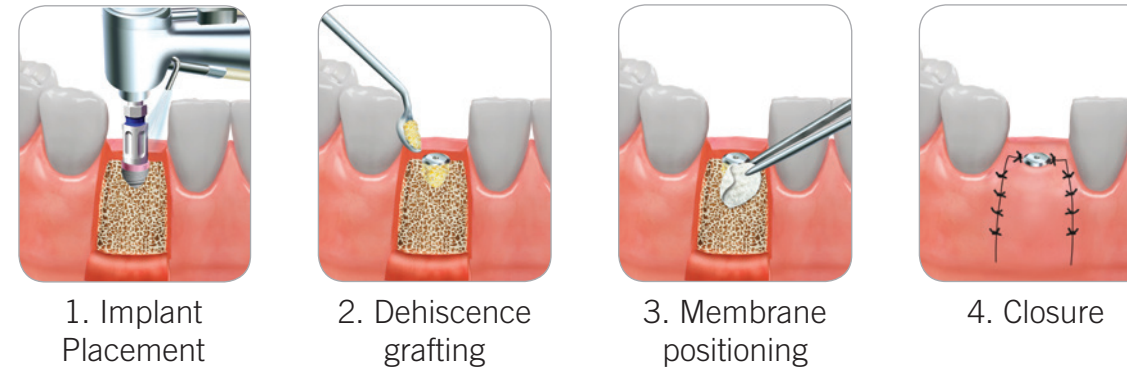

Extraction Socket

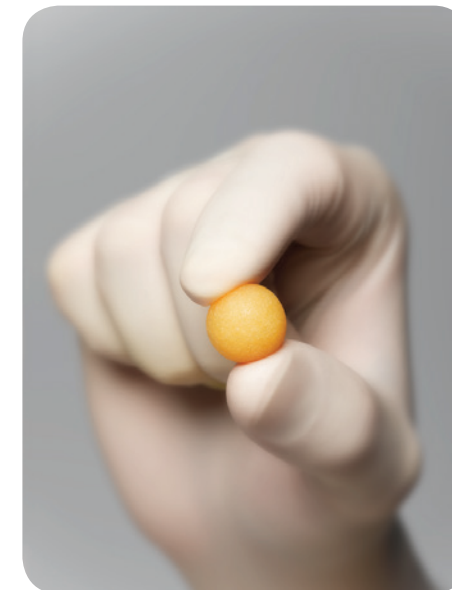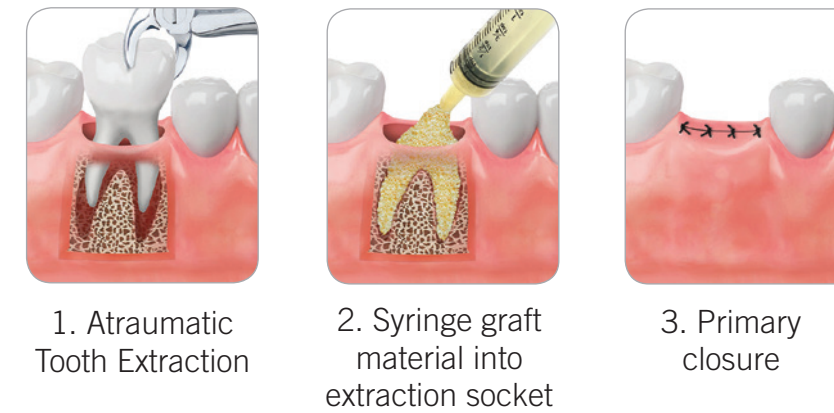

Easy Handling

DynaBlast® is available as a moldable putty or a flowable paste in a prefilled syringe. Both are premixed for easy application without refrigeration—no additional time required for thawing, mixing or setting the material.

Cells

Patient's own osteoblasts, osteoclasts and mesenchymal stem cells integrate with DynaBlast graft material

Signals

Proven osteoinductive properties of demineralized bone stimulates the growth of new bone. ¹⁻³

Scaffolds

Mineralized cancellous chips in this unique formulation provide the osteoconductive elements to encourage the attachment of osteogenic precursor cells. Scaffold maintains volume of the defect site while anchoring surrounding cells.

DynaBlast® Provides a Unique Bone Growth Dynamic

DynaBlast®

Dematerialized Bone Matrix with Cancellous Bone

DynaBlast® paste (in syringe) 0.5cc

DynaBlast® paste (in syringe) 1.0cc

DynaBlast® paste (in syringe) 3.0cc

DynaBlast® putty (in vial) 1.0cc

DynaBlast® putty (in vial) 2.5cc

DynaBlast® putty (in vial) 5.0cc

Catalog Number

10.210.1050

10.210.1060

10.210.1070

10.220.1030

10.220.1040

10.220.1050

1. Clokie C, Urist MR. Bone Morphogenetic Protein Excipients: Comparative Observations on Poloxamer. *Plast Reconstr Surg* 2000;105: 628.

2. Coulson R, Clokie C, Peel S. Collagen and a Thermally Reversible Poloxamer Deliver Demineralized Bone Matrix (DBM) and Biologically Active Proteins to Sites of Bone Regeneration. *Proc Portland Bone Symp* 1999. p. 619-637.

3. Kay, JF. Validated Assay for Measuring Osteoinductivity of Human Demineralized Bone Matrix, *Integra OrthoBiologics* Publication.