



RNA and DNA from soil

User manual

NucleoBond[®] RNA Soil Mini




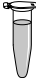








DNA Set for NucleoBond[®] RNA Soil
Mini

December 2017 / Rev. 01

RNA and DNA from soil

Protocol at a glance (Rev. 01)

NucleoBond® RNA Soil

1 Prepare sample		Fill tubes up to the 1 mL graduation mark
2 Add lysis buffers		800 µL Buffer E1 <i>Optional: 100 µL OPT</i> 100 µL Phenol:Chloroform:Isoamylalcohol Vortex 2 s
3 Bead beating		Vortex horizontally 5 min at RT
4 Lysate clearing	 	20,000 x g, 2 min Transfer the supernatants to a 1.5 mL tube
5 Adjust binding conditions	 	125 µL Buffer E2 per mL supernatant Vortex 5 s 2 min RT 20,000 x g, 2 min
6 Column setup and equilibration		1.5 mL Buffer EQU
7 Load column		Load supernatant on NucleoBond® RNA Column
8 1 st washing step		250 µL Buffer E3
9 2 nd washing step		2 mL Buffer E4
10 RNA elution		1 mL Buffer ERNA

RNA and DNA from soil

Protocol at a glance (Rev. 01)

NucleoBond® RNA Soil



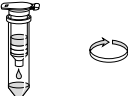
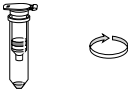

11 Precipitate nucleic acids		700 µL Isopropanol Vortex 5 s
12 Bind nucleic acids		Load sample 11,000 x <i>g</i> 15 s
13 Washing		500 µL Buffer E5 11,000 x <i>g</i> , 15 s Discard flow through Repeat this step once.
14 Drying		11,000 x <i>g</i> , 1 min
15 Elute nucleic acids		50–100 µL RNase-free H ₂ O 11,000 x <i>g</i> , 1 min

Table of contents

1 Components	4
1.1 Kit contents	4
1.2 Reagents, consumables and equipment to be supplied by user	5
2 Product description	6
2.1 Basic principle	6
2.2 Kit specifications	7
2.3 Bead tubes and lysis	7
3 Storage conditions and preparation of working solutions	8
4 Safety instructions	9
5 Protocol	10
6 Appendix	14
6.1 Troubleshooting	14
6.2 Ordering information	16
6.3 Product use restriction / warranty	17

1 Components

1.1 Kit contents

REF	NucleoBond® RNA Soil Mini		DNA Set for NucleoBond® RNA Soil Mini	
	10 preps 740142.10	50 preps 740142.50	10 preps 740143.10	50 preps 740143.50
Lysis Buffer E1	13 mL	70 mL	-	-
Buffer OPT	3 mL	10 mL	-	-
Equilibration Buffer EQU	30 mL	100 mL	-	-
Binding Buffer E2	10 mL	10 mL	-	-
Wash Buffer E3	13 mL	30 mL	-	-
Wash Buffer E4	30 mL	125 mL	-	-
RNA Elution Buffer ERNA	13 mL	60 mL	-	-
DNA Elution Buffer EDNA	-	-	13 mL	60mL
Wash Buffer E5 (Concentrate)*	6 mL	2 x 6 mL	6 mL	2 x 6 mL
RNase-free H ₂ O	13 mL	13 mL	13 mL	13 mL
NucleoSpin® Bead Tubes Type A	10	50	-	-
NucleoBond® RNA Mini Columns	10	50	-	-
NucleoSpin® RNA Columns	10	50	-	-
NucleoSpin® DNA Columns	-	-	10	50
Collection Tubes	-	-	10	50
Plastic Washers	10	10		
User manual	1	1		

* For preparation of working solutions and storage conditions see section 3.

1.2 Reagents, consumables and equipment to be supplied by user

Reagents

- Phenol:Chloroform:Isoamylalcohol (25:24:1 v/v), e.g., VWR Cat. No. 0883-100ML
- 96–100% ethanol
- Isopropanol

Consumables

- 1.5 mL centrifuge tubes
- 2.0 mL centrifuge tubes
- RNase-free pipette tips
- NucleoBond® Rack Small, alternatively, 50 mL centrifuge tubes or standard laboratory flask

Equipment

- Manual pipettors
- Centrifuge for microcentrifuge tubes
- Equipment for sample disruption and homogenization (see section 2.3)
- Personal protection equipment (e.g., labcoat, gloves, goggles)

2 Product description

2.1 Basic principle

The **NucleoBond® RNA Soil Mini Kit** is a downscaled version of the NucleoBond® RNA Soil Kit (740140.20) and the corresponding DNA Elution Set (740141.20). A smaller anion exchange column reduces both, capacity and dead volume, resulting in a decreased elution volume. The final purification, concentration and desalting step can be performed in centrifuge for microcentrifuge tubes.

For downstream applications like PCR or RT-PCR, high yields of nucleic acid might not be necessary. In contrast to the larger NucleoBond® RNA Soil kit a single NucleoSpin® Bead Tube Type A, filled to the 1 mL mark with 250 to 500 mg of soil sample is sufficient and will result in ¼ of the larger kit's yield (0.25–2.5 µg RNA). Nucleic acid yield is mainly influenced by the type of soil and varies strongly from soil to soil.

Organisms in a soil sample are lysed by bead beating in the presence of Lysis Buffer E1 (check for precipitated SDS!) and Phenol:Chloroform:Isoamylalcohol. Buffer OPT reduces the adsorption of nucleic acids to clay and mineralic soil components, but will also increase the contamination with humic acids if present. Performance cannot be predicted and depends on the soil composition, but as a rule of thumb a sample with a high content of organic and **humic components**, e.g., forest soil, should be lysed **without Buffer OPT** while predominantly **mineralic soils**, e.g., river sediments and clay, should be lysed **with Buffer OPT**.

Unlysed components are sedimented by centrifugation. The supernatant of the bead tube is transferred into one fresh 1.5 mL tube and mixed with Binding Buffer E2. The incubation and centrifugation steps that follow can be used to prepare the NucleoBond® RNA Mini Columns. The columns are fixed on a NucleoBond® Rack Small (not supplied, see ordering information) or by the use of supplied Plastic Washers on top of a standard laboratory flask or a 50 mL centrifuge tube. Columns are equilibrated by Buffer EQU to prepare the anion exchange chromatography columns.

Once the columns are equilibrated, the clear supernatant of the centrifuged samples is loaded onto the purification columns.

Wash Buffer E3 is used as a first washing buffer for the column matrix. Nucleic acids and soluble polyanionic molecules including some fraction of humic substances (if present) will bind to the anion exchange surface, resulting in a brown silica matrix. Buffer E4 removes contaminants in a second washing step.

RNA is eluted by Buffer ERNA and afterwards DNA can be eluted by Buffer EDNA, supplied in the additional **DNA Set for NucleoBond® RNA Soil Mini**. If DNA isolation is not necessary, the DNA elution step may be skipped.

Isopropanol is added to the eluates and the mixture is bound to NucleoSpin® RNA Columns or NucleoSpin® DNA Columns by centrifugation. The columns are washed with Buffer E5, dried and eluted in RNase-free H₂O. Elution with 50 µL of RNase-free H₂O will result in highly concentrated nucleic acids but total yield might be decreased. Elution with 100 µL of RNase-free H₂O will result in a high total yield but decreased concentration.

Eluted RNA and DNA fractions might contain a fraction of the other nucleic acid. If necessary, digest DNA or RNA enzymatically. Contact MACHEREY-NAGEL (tech-bio@mn-net.com) for RNase A and DNase Buffer Sets and protocols if required.

2.2 Kit specifications

The **NucleoBond® RNA Soil Mini** kit is designed for the isolation of pure RNA from soil, the additional **DNA Set for NucleoBond® RNA Soil Mini** for the isolation of pure DNA from soil.

Kit specifications at a glance		
Parameter	NucleoBond® RNA Soil	DNA Set for NucleoBond® RNA Soil
Format	Anion exchange chromatography and silica spin purification	Elution Set (NucleoBond® RNA Soil Mini kit required)
Processing	Manual	Manual
Sample material	0.25–0.5 g soil	NucleoBond® RNA Soil Mini kit required
Elution volume	50–100 µL	50–100 µL
Preparation time	60 min/12 preps	15 min/12 preps
Typical yield	0.25–2.5 µg	1.25–12.5 µg

2.3 Bead tubes and lysis

Bead beating can be performed with any dedicated machine (e.g., FastPrep® instrument - MP Biomedicals, Precellys – Bertin Corp., Retsch mill or the like) which can accommodate 2 mL tubes. Alternatively, the MN Bead Tube Holder can be used in combination with a Vortex Genie II – Scientific Industries for 2 mL tubes.

Bead beating for 5 minutes at full speed is sufficient to lyse the cells when using a vortexer. Parameters for other beat beating machines need to be established individually. The following parameters are a good starting point for further optimization:

- 5 m/s for 30 s on a FastPrep® instrument (twice to three times)
- 30 Hz for 5 min on a Retsch mill
- 2700 rpm for 5 min on a vortexer (horizontally)

Do not overfill the NucleoSpin® Bead Tubes! Lysis is most effective when having a ratio of one part beads to one part sample to two parts of lysis buffer in a 2 mL tube. Further increase of sample material will result in a reduced lysis activity and therefore in a reduced yield. For convenient processing the 1 mL mark of the skirted bead tubes can be used as a fill marking. NucleoSpin® Bead Tubes Type A filled to the 1 mL mark will contain the optimal range of 250–500 mg of sample, depending on soil density and water content.

3 Storage conditions and preparation of working solutions

Attention: Buffer E4 contains guanidine thiocyanate! Wear gloves and goggles!

CAUTION: Buffer E4 contains guanidine thiocyanate which can form highly reactive compounds when combined with bleach (sodium hypochloride). DO NOT add bleach or acidic solutions directly to the sample preparation waste.

Attention: Phenol:Chloroform:Isoamylalcohol is poisonous! Read the safety data sheet, supplied with the Phenol:Chloroform:Isoamylalcohol, and follow the instructions for use and disposal carefully! Take appropriate safety measures when working with Phenol:Chloroform:Isoamylalcohol and discard waste according to legal guidelines!



- All kit components can be stored at room temperature (18–25 °C) and are stable for at least one year.
- **Sodium dodecyl sulfate (SDS) in Buffer E1 may precipitate if stored at temperatures below 20 °C.** If a precipitate is observed in Buffer E1, incubate bottle at 30–40 °C for several minutes and mix well.

Before starting the NucleoBond® RNA Soil protocol prepare the following:

- Add 24 mL of 96–100 % ethanol to each bottle of **Wash Buffer E5**.

4 Safety instructions

Only harmful features do not need to be labeled with H and P phrases up to 125 mL or 125 g. *Mindergefährliche Eigenschaften müssen bis 125 mL oder 125 g nicht mit H- und P-Sätzen gekennzeichnet werden.*

Component	Hazard contents	GHS symbol	Hazard phrases	Precaution phrases
<i>Inhalt</i>	<i>Gefahrstoff</i>	<i>GHS-Symbol</i>	<i>H-Sätze</i>	<i>P-Sätze</i>
Wash Buffer E4	ethanol 5–20 % <i>Ethanol 5–20 %</i> CAS 64-17-5d	 WARNING <i>ACHTUNG</i>	226	210
RNA Elution Buffer ERNA	ethanol 5–20 % <i>Ethanol 5–20 %</i> CAS 64-17-5d	 WARNING <i>ACHTUNG</i>	226	210



The symbol shown on labels refers to further safety information in this section. *Das auf Etiketten dargestellte Symbol weist auf weitere Sicherheitsinformationen dieses Kapitels hin.*

Hazard phrases

H 226 Flammable liquid and vapour.
Flüssigkeit und Dampf entzündbar.

Precaution phrases

P 210 Keep away from heat, hot surfaces, sparks, open flames and other ignition sources. No smoking.
Von Hitze, heißen Oberflächen, Funken, offenen Flammen sowie anderen Zündquellenarten fernhalten. Nicht rauchen.

5 Protocol

Before starting the preparation:

- Check Lysis Buffer E1 for precipitated SDS according to section 3.
- Check if Wash Buffer E5 was prepared according to section 3.
- Check if Phenol:Chloroform:Isoamylalcohol (25:24:1 v/v) is available.
- Equilibrate Phenol:Chloroform:Isoamylalcohol to room temperature and mix well!



1 Prepare Sample

Fill **1 NucleoSpin® Bead Tube Type A** with sample material up to the **1 mL graduation mark**.

Usually: 250–500 mg sample per tube.



fill tubes up to the 1 mL graduation mark

2 Add lysis buffers

Add **800 µL Lysis Buffer E1**.

Optional: Add 100 µL Buffer OPT.

Buffer OPT is recommended for soil with a predominantly mineralic matrix (e. g., clay, sediment) but negative for predominantly organic soil. It cannot be identified in advance, if addition of Buffer OPT increases yield or decreases purity.

Add **100 µL Phenol:Chloroform:Isoamylalcohol** (25:24:1 v/v).

Mix Phenol:Chloroform:Isoamylalcohol well before use! Take appropriate safety measures for working with Phenol:Chloroform:Isoamylalcohol!

Vortex bead tube for 2 s to distribute beads, sample and buffers homogeneously.



+ 800 µL Buffer E1

optional:
+ 100 µL Buffer OPT

+ 100 µL Phenol:Chloroform:Isoamylalcohol

Vortex 2 s

3 Bead Beating

Lyse cells by **bead beating**

- **5 m/s for 30 s** on a FastPrep® Instrument (twice to three times)
- **30 Hz for 5 min** on a Retsch mill
- **2700 rpm for 5 min** on a vortexer (horizontally)



Vortex RT, 5 min

See section 2.3 for details.

Recommended: use MN Bead Tube Holder in combination with a Vortex-Genie II (Scientific Industries Inc.).

4 Lysate clearing

Centrifuge samples for **2 min** at **full speed** (best: 20,000 x g).



20,000 x g
2 min

Transfer the supernatant into a fresh 1.5 mL tube (not supplied) and determine the approximate volume.



Transfer the supernatants to a 1.5 mL tube

Avoid transferring debris or the organic phase which is a thin layer between the aqueous phase and the solid debris and usually much darker than the aqueous supernatant!

5 Adjust binding conditions

For each mL of supernatant, add **125 µL Buffer E2**.



+ 125 µL E2 per mL supernatant

Example:

to 0.6 mL supernatant add 75 µL Buffer E2

to 0.7 mL supernatant add 87.5 µL Buffer E2

to 0.8 mL supernatant add 100 µL Buffer E2

Vortex for 5 s.

Incubate for **2 min** at **room temperature** and **centrifuge for 2 min** at **full speed** (best: 20,000 x g) afterwards.



Vortex 5 s
2 min RT

20,000 x g
2 min

Proceed to step 6 (column setup and equilibration) during the incubation and centrifugation time.

6 Column setup and equilibration

For each soil sample, combine a **NucleoBond® RNA Mini Column** with a Plastic Washer (reusable) and place the combination on top of a laboratory flask or 50 mL centrifuge tube.

Equilibrate columns with **1.5 mL Buffer EQU**.

All steps involving the NucleoBond® RNA Column are performed with gravity flow. Do not use vacuum!

See ordering information for a convenient NucleoBond® Rack Small instead of Plastic Washers and centrifuge tubes or flasks.



+ 1.5 mL Buffer EQU

7 Load column

Load supernatant from centrifugation (step 5) into the column.



Load supernatant

Avoid transferring debris.

Attention: Flow through contains remaining PCI! Discard separately according to local guidelines!

8 1st washing step

Wash column with **250 µL Buffer E3**.



**+ 250 µL
Buffer E3**

9 2nd washing step

Wash the NucleoBond® RNA Mini Column with **2 mL Buffer E4**.



**+ 2 mL
Buffer E4**

10 RNA elution

Place a fresh 2.0 mL tube below the NucleoBond® RNA Mini Column.

Elute RNA with **1 mL Buffer ERNA**.



**+ 1 mL
Buffer ERNA**

Optional: DNA elution (DNA Set for NucleoBond® RNA Soil Mini required)

Place a fresh 2.0 mL tube below the NucleoBond® RNA Mini Column.

Elute DNA with **1 mL Buffer EDNA**.

11 Precipitate nucleic acids

Add **700 µL isopropanol** (not supplied) to each eluate.

Vortex 5 s.



**+ 700 µL
Isopropanol**

12 Bind nucleic acids

Load **600 µL** of the resulting mixture into a **NucleoSpin® RNA Column** (blue, RNA) or a **NucleoSpin® DNA Column** (green, DNA) in a Collection Tube.

Centrifuge NucleoSpin® Columns for **15 s at 11,000 x g**

Discard flow through and place the NucleoSpin® Columns back into the Collection tubes.

Repeat this step until all of the mixture has passed the silica of the NucleoSpin® Columns.



**11,000 x g
15 s**

13 Washing

Add 500 µL Buffer E5 to the center of the silica.

Centrifuge NucleoSpin® Columns for **15 s** at **11,000 x g**

Discard flow through, place the NucleoSpin® Columns back into the Collection tubes.

Repeat this step once.



**+ 500 µL
Buffer E5**



**11,000 x g
15 s**

14 Drying

Centrifuge empty NucleoSpin® Columns for **1 min** at **11,000 x g**.

Transfer NucleoSpin® Columns to a fresh 1.5 mL centrifuge tube (not supplied).

Avoid contact between the outlet of the columns and the flow through.



**11,000 x g
1 min**

15 Elute nucleic acids

Add **50 µL** (high concentration) to **100 µL** (high yield) **RNase-free H₂O** to the center of the silica.

Centrifuge NucleoSpin® Columns for **1 min** at **11,000 x g**



**+ 100 µL
RNase-free H₂O**



**11,000 x g
1 min**

6 Appendix

6.1 Troubleshooting

Problem	Possible cause and suggestions
No or low nucleic acid yield	<p data-bbox="333 309 978 331"><i>Nucleic acid adsorbed to soil matrix</i></p> <ul data-bbox="333 347 978 496" style="list-style-type: none"> <li data-bbox="333 347 978 400">• Use Buffer OPT or increase volume of Buffer OPT. Reduce volume of Buffer E1 when increasing volume of Buffer OPT <li data-bbox="333 416 978 496">• Increase volume of Phenol:Chloroform:Isoamylalcohol. Some soils do not require Phenol:Chloroform:Isoamylalcohol, but in most cases addition is mandatory to isolate long, intact RNA.
	<p data-bbox="333 512 978 534"><i>Storage of soil sample</i></p> <ul data-bbox="333 550 978 608" style="list-style-type: none"> <li data-bbox="333 550 978 608">• Storage conditions might alter RNA integrity and total amount of nucleic acids.
	<p data-bbox="333 624 978 646"><i>Sample amount</i></p> <ul data-bbox="333 662 978 815" style="list-style-type: none"> <li data-bbox="333 662 978 742">• Do not fill the NucleoSpin[®] Bead Tubes Type A higher than the 1 mL mark as this will result in decreased bead beating efficiency. “Less is more”. <li data-bbox="333 758 978 815">• Check if the beads are moving vigorously in the tubes during the bead beating.
	<p data-bbox="333 831 978 853"><i>Room temperature</i></p> <ul data-bbox="333 869 978 927" style="list-style-type: none"> <li data-bbox="333 869 978 927">• Best performance is achieved at 18–25 °C. Do not store buffers in the fridge or on ice!
	<p data-bbox="333 943 978 965"><i>Sample material</i></p> <ul data-bbox="333 981 978 1114" style="list-style-type: none"> <li data-bbox="333 981 978 1114">• Usually the RNA content in a soil sample is much lower than the DNA content. DNA might originate from dead organic material as well. If the isolated RNA concentration is too low, try the larger version of this kit (see ordering information: REF 740140 and 740141).

Problem

Possible cause and suggestions

Poor downstream performance

Humic acid contamination

- To verify PCR inhibition by humic substances as root cause for poor downstream application results, check if a dilution series (e.g., 1:10, 1:100, 1:1000) of the eluates will restore functionality in downstream application. If not, other reasons must be taken into account.
 - Buffer OPT will release PCR inhibitors from soils. Skip addition or decrease volume.
 - Reduce bead beating time.
-

DNA contamination

- Separation of DNA and RNA is not complete. Degraded DNA from dead cells will elute with intact RNA of living cells. If residual DNA in the RNA fraction or RNA in the DNA fraction is influencing downstream application results, digest RNA and DNA enzymatically. Contact MACHEREY-NAGEL (techbio@mn-net.com) for enzymes and protocols if necessary.
-

6.2 Ordering information

Product	REF	Pack of
NucleoBond® RNA Soil	740140.20	20 preps
DNA Set for NucleoBond® RNA Soil	740141.20	20 preps
NucleoBond® RNA Soil Mini	740142.10 / .50	10 / 50 preps
DNA Set for NucleoBond® RNA Soil Mini	740143.10 / .50	10 / 50 preps
NucleoSpin® Soil	740780.10 / .50 / .250	10 / 50 / 250 preps
NucleoSpin® RNA Stool	740130.10 / .50	10 / 50 preps
NucleoSpin® DNA Stool	740472.10 / .50 / .250	10 / 50 / 250 preps
rDNase Set	740963	500 reactions
Liquid RNase A	740397	500 reactions
MN Bead Tube Holder	740469	1 piece
NucleoSpin® Bead Tubes Type A (0.6–0.8 mm ceramic beads)	740786.50	50 pieces
NucleoBond® Rack Small	740562	1 piece

6.3 Product use restriction/warranty

NucleoBond® RNA Soil Mini and **DNA Set NucleoBond® RNA Soil Mini** components are intended, developed, designed, and sold FOR RESEARCH PURPOSES ONLY, except, however, any other function of the product being expressly described in original MACHEREY-NAGEL product leaflets.

MACHEREY-NAGEL products are intended for GENERAL LABORATORY USE ONLY! MACHEREY-NAGEL products are suited for QUALIFIED PERSONNEL ONLY! MACHEREY-NAGEL products shall in any event only be used wearing adequate PROTECTIVE CLOTHING. For detailed information please refer to the respective Material Safety Data Sheet of the product! MACHEREY-NAGEL products shall exclusively be used in an ADEQUATE TEST ENVIRONMENT. MACHEREY-NAGEL does not assume any responsibility for damages due to improper application of our products in other fields of application. Application on the human body is STRICTLY FORBIDDEN. The respective user is liable for any and all damages resulting from such application.

DNA/RNA/PROTEIN purification products of MACHEREY-NAGEL are suitable for IN VITRO-USES ONLY!

ONLY MACHEREY-NAGEL products specially labeled as IVD are also suitable for IN VITRO-diagnostic use. Please pay attention to the package of the product. IN VITRO-diagnostic products are expressly marked as IVD on the packaging.

IF THERE IS NO IVD SIGN, THE PRODUCT SHALL NOT BE SUITABLE FOR IN VITRO-DIAGNOSTIC USE!

ALL OTHER PRODUCTS NOT LABELED AS IVD ARE NOT SUITED FOR ANY CLINICAL USE (INCLUDING, BUT NOT LIMITED TO DIAGNOSTIC, THERAPEUTIC AND/OR PROGNOSTIC USE).

No claim or representations is intended for its use to identify any specific organism or for clinical use (included, but not limited to diagnostic, prognostic, therapeutic, or blood banking). It is rather in the responsibility of the user or - in any case of resale of the products - in the responsibility of the reseller to inspect and assure the use of the DNA/RNA/protein purification products of MACHEREY-NAGEL for a well-defined and specific application.

MACHEREY-NAGEL shall only be responsible for the product specifications and the performance range of MN products according to the specifications of in-house quality control, product documentation and marketing material.

This MACHEREY-NAGEL product is shipped with documentation stating specifications and other technical information. MACHEREY-NAGEL warrants to meet the stated specifications. MACHEREY-NAGEL's sole obligation and the customer's sole remedy is limited to replacement of products free of charge in the event products fail to perform as warranted. Supplementary reference is made to the general business terms and conditions of MACHEREY-NAGEL, which are printed on the price list. Please contact us if you wish to get an extra copy.

There is no warranty for and MACHEREY-NAGEL is not liable for damages or defects arising in shipping and handling (transport insurance for customers excluded), or out of accident or improper or abnormal use of this product; defects in products or components not manufactured by MACHEREY-NAGEL, or damages resulting from such non-MACHEREY-NAGEL components or products.

MACHEREY-NAGEL makes no other warranty of any kind whatsoever, and SPECIFICALLY DISCLAIMS AND EXCLUDES ALL OTHER WARRANTIES OF ANY KIND OR NATURE WHATSOEVER, DIRECTLY OR INDIRECTLY, EXPRESS OR IMPLIED, INCLUDING, WITHOUT LIMITATION, AS TO THE SUITABILITY, REPRODUCTIVITY, DURABILITY, FITNESS FOR A PARTICULAR PURPOSE OR USE, MERCHANTABILITY, CONDITION, OR ANY OTHER MATTER WITH RESPECT TO MACHEREY-NAGEL PRODUCTS.

In no event shall MACHEREY-NAGEL be liable for claims for any other damages, whether direct, indirect, incidental, compensatory, foreseeable, consequential, or special (including but not limited to loss of use, revenue or profit), whether based upon warranty, contract, tort (including negligence) or strict liability arising in connection with the sale or the failure of MACHEREY-NAGEL products to perform in accordance with the stated specifications. This warranty is exclusive and MACHEREY-NAGEL makes no other warranty expressed or implied.

The warranty provided herein and the data, specifications and descriptions of this MACHEREY-NAGEL product appearing in MACHEREY-NAGEL published catalogues and product literature are MACHEREY-NAGEL's sole representations concerning the product and warranty. No other statements or representations, written or oral, by MACHEREY-NAGEL's employees, agent or representatives, except written statements signed by a duly authorized officer of MACHEREY-NAGEL are authorized; they should not be relied upon by the customer and are not a part of the contract of sale or of this warranty.

Product claims are subject to change. Therefore please contact our Technical Service Team for the most up-to-date information on MACHEREY-NAGEL products. You may also contact your local distributor for general scientific information. Applications mentioned in MACHEREY-NAGEL literature are provided for informational purposes only. MACHEREY-NAGEL does not warrant that all applications have been tested in MACHEREY-NAGEL laboratories using MACHEREY-NAGEL products. MACHEREY-NAGEL does not warrant the correctness of any of those applications.

Last updated: 07/2010, Rev. 03

Please contact:
MACHEREY-NAGEL GmbH & Co. KG
Tel.: +49 24 21 969-270
tech-bio@mn-net.com

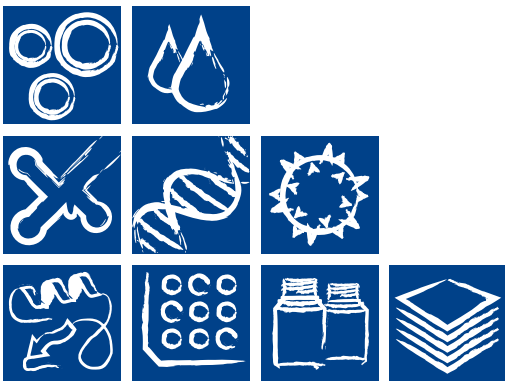
Trademarks:

NucleoBond and NucleoSpin are registered trademarks of MACHEREY-NAGEL GmbH & Co KG

Fastprep is a registered trademark of MP Biomedicals, LLP

Vortex-Genie is a registered trademark of Scientific Industries

All used names and denotations can be brands, trademarks, or registered labels of their respective owner – also if they are not special denotation. To mention products and brands is only a kind of information (i.e., it does not offend against trademarks and brands and can not be seen as a kind of recommendation or assessment). Regarding these products or services we can not grant any guarantees regarding selection, efficiency, or operation.



www.mn-net.com

MACHEREY-NAGEL



MACHEREY-NAGEL GmbH & Co. KG · Neumann-Neander-Str. 6-8 · 52355 Düren · Germany

DE / International:

Tel.: +49 24 21 969-0
 Fax: +49 24 21 969-199
 E-mail: info@mn-net.com

CH:

Tel.: +41 62 388 55 00
 Fax: +41 62 388 55 05
 E-mail: sales-ch@mn-net.com

FR:

Tel.: +33 388 68 22 68
 Fax: +33 388 51 76 88
 E-mail: sales-fr@mn-net.com

US:

Tel.: +1 484 821 0984
 Fax: +1 484 821 1272
 E-mail: sales-us@mn-net.com