



Chair *Caddie*

RETRO-FIT



Turns...Rolls...and Brakes for Safety!



Enables a standing person get seated... without the fear of falling.



Enables a caregiver to move a seated person up-to the table...with ease, grace, and safety!

- Extends life of chairs by 4-6 years
- Eliminates damage to floors
- Attaches to existing chairs



- Reduces caregiver fatigue
- Reduces workplace injuries
- Provides dignity to persons with limited mobility

 www.comfortek.com/cc

Chair *Caddie*

RETRO-FIT

1. **Attaches to existing chairs** (wood or metal)
2. Low profile design – chair is raised $\frac{3}{4}$ " above the floor
3. Fits inside the footprint of an existing chair
4. Medical grade casters:
 - Chair Caddie A (3" front / 2" rear)
 - Chair Caddie B (4" front / 3" rear)
5. Anti-tip glides
6. Front wheel tension control
7. Foot activated locking rear caster
8. Rated for 450 lbs.
9. Attaches to arm or armless chairs
10. Adjustable - Custom sizing available
11. Eliminates **pushing, pulling, shoving,** and **twisting**



PATENT
PENDING



The Chair CADDIE's mobility features (**turn...roll...and brake** for safety) help to;

- **reduce** a standing person's risk of falling
- **reduce** caregiver's risk of incurring a work / personal injury!





FEBRUARY 2023

The Chair CADDIE attaches to the base of a kitchen/dining chair and eliminates the physical demands expected of the caregiver when assisting a standing person get seated and moved up-to the table!

The Chair CADDIE's mobility features (**turn...roll...and brake** for safety) help to;

- **reduce** a standing person's risk of falling
- **reduce** caregiver's risk of incurring a work / personal injury!

With the Chair CADDIE, caregivers are able to guide a standing person to a chair, help them get seated and then ROLL them to the table with **ease, grace and safety**. This eliminates the awkward task of pushing-pulling-shoving-twisting when moving the seated person up-to the table!




1. Getting the person seated in the chair

For the safety of the person being seated, ask the person to **REMAIN STANDING**, supported either by the table, crutches, walker/rollator or another caregiver.

With the person standing, a caregiver ROLLS the Chair/CADDIE directly up behind the standing person. When the chair contacts the back of the person's legs, the brakes of the Chair/CADDIE should be engaged.

It is now safe for the standing person to sit down!

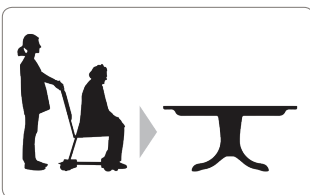

2. Moving the seated person up-to the table

For the safety of the caregiver, with the person comfortably seated in the chair, the caregiver ROLLS the seated person up-to the table. Once positioned at the table, the brakes should be engaged.



3. Moving the seated person away-from the table

For the safety of the caregiver, the caregiver releases the brakes and gently ROLLS and TURNS the seated person away from the table.

For the safety of the seated person, when the walker/rollator has been re-positioned directly in front of the chair, the seated person rises from the chair and walks directly away from the table.



DESIGN ELEMENTS OF THE CHAIR CADDIE



1. The Chair CADDIE is a fully integrated “undercarriage unit” which strengthens and supports the frame of the chair to which it is attached.
2. Side rails hold the legs of the chair in place preventing them from splaying outwards. The connecting arms are adjustable enabling the Chair CADDIE to attach to most existing armed and armless chairs. The Chair CADDIE is rated to carry 450 lbs.
3. Casters (oversized) attached directly to the “undercarriage unit” (not to the legs of the chair) enable the chair with the person seated in the chair to move easily on tile, laminate or carpeted surfaces, while raising the seat height of the chair only 3/4”.
4. For stability, the front casters are positioned in-line with the front legs of the chair. The wide wheel-base stance provides unparalleled stability preventing the chair from tipping sideways. Anti-tip glides installed under the front legs prevent the chair from tipping forwards.
5. For mobility, the “undercarriage unit” is fitted with rotating rear casters.
6. For safety, both rear casters are fitted with foot activated brakes. It is recommended the brakes are engaged whenever the chair is not attended by a caregiver.
7. Additionally, the entire “undercarriage unit” fits inside the “foot print” of the chair to which it is attached and does not pose a tripping hazard to the caregiver or other guests.



Online Reference: www.comfortek.com/cc

OPERATING GUIDELINES



info@comfortek.com

1.888.678.2060