

**vive<sup>®</sup>**

# **STAND ALONE TOILET SAFETY RAIL**

User Manual

LVA1O23



**vivehealth.com**

## PRODUCT DESCRIPTION


The Vive stand alone toilet safety rail is a bathroom safety product designed to protect the disabled and elderly. It easily fits around most toilets to provide the user with support that makes sitting down and standing up safe and easy. Its armrests are padded for extra comfort which is built to withstand up to 300 lbs. Its freestanding design makes it easy to remove and store as needed. Protected by a lifetime warranty.

## WHAT'S INCLUDED

- 1x Stand Alone Toilet Safety Rail
- Lifetime Warranty

## PARTS LIST

- A) 2x Frames with Armrests
- B) 2x U-Shape Support Bars
- C) 1x Magazine Basket
- D) 8x Hand Knob Screws

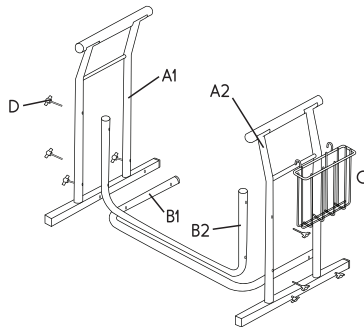

## SPECIFICATIONS

- Length: 26.5"
- Height: 25.5"
- Depth: 19"
- Width Between Armrests: 21"
- Product Weight: 9 LBS
- Max User Weight: 300 LBS


## WARNING ⚠

- Check periodically to make sure all fixings are tight and secure before every use.

## ASSEMBLY INSTRUCTIONS

No Tools Required! 

1. Remove product from carton and read instructions carefully before attempting to install the product. No tools required.
2. Attach frame with armrest (A1) to U-shape support bar (B1) using 2 hand knob screws (D). Insure the holes on the vertical post are facing forward.
3. Attach the other frame with armrest (A2) to the other side of the U-shape support bar (B1). Fasten with 2 hand knob screws.
4. Attach the last U-shape support bar (B2) to the frame with armrest (A2) using 2 hand knob screws for each side like steps 2 and 3. Insure the support bar is in a U-shape.
5. Attach optional magazine basket to either side of the frame with armrest.

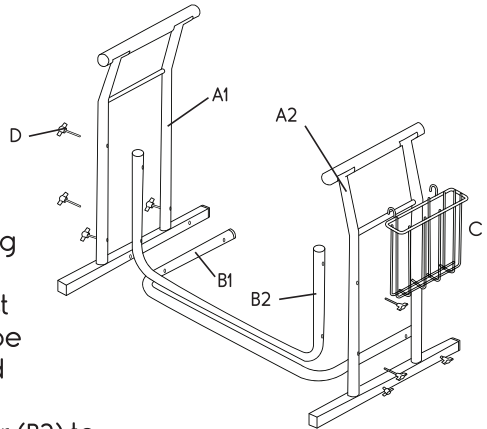

## INSTALLATION INSTRUCTIONS

Once fully assembled, position the product up against the toilet so the center supports are wrapping around the base of the toilet. The armrests are located on either side for support when sitting and standing up.

## WARRANTY REGISTRATION

This product is guaranteed to be free from manufacturing defects. We will replace or exchange any product that is defective to your complete satisfaction. This warranty extends for the lifetime of the purchaser, as long as the product has not been misused or abused.

**To ensure your warranty,  
please register your product at  
[www.vivehealth.com/register](http://www.vivehealth.com/register).**

**Product Code: LVA1O23**



We sincerely appreciate your business.  
We strive to provide you with the best  
quality products at great value.

If you have any questions please contact us.



[service@vivehealth.com](mailto:service@vivehealth.com)



1-800-487-3808

Monday - Friday 9am - 5pm EST