

REORDER NO. ZZAN602



## **ANALOG TENS UNIT** Instruction Manual



## TABLE OF CONTENTS

Index .....	2
1. Introduction .....	3
2. Cautions .....	4
3. Warnings .....	5
4. General Description.....	6
5. Construction.....	6
6. Technical Specifications .....	7
7. Replaceable Parts .....	8
8. Accessories.....	8
9. Graphic Symbols .....	9
10. Parameter Controls .....	9
11. Attachment of Electrodes Lead Wires.....	10
12. Lead Wire Maintenance .....	10
13. Electrode Options .....	11
14. Electrode Placement .....	11
15. Tips For Skin Care .....	13
16. Application of Re-usable Self Adhesive Electrodes .....	14
17. Adjusting the Controls .....	15
18. Battery Information .....	18
19. Maintenance, Transportation and Storage of TENS Device .....	19
20. Safety Control .....	19
21. Malfunctions .....	20
22. Conformity to Safety Standards.....	20
23. Warranty .....	20

### **EXPLANATION OF PAIN**

Pain is a warning system and the body's method of telling us that something is wrong. Pain is important; without it abnormal conditions may go undetected, causing damage or injury to vital parts of our bodies.

Even though pain is a necessary warning signal of trauma or malfunction in the body, nature may have gone too far in its design. Aside from its value in diagnosis, long-lasting persistent pain serves no useful purpose. Pain does not begin until a coded message travels to the brain where it is decoded, analyzed, and then reacted to. The pain message travels from the injured area along the small nerves leading to the spinal cord. Here the message is switched to different nerves that travel up the spinal cord to the brain. The pain message is then interpreted, referred back and the pain is felt.

### **EXPLANATION OF TENS**

Transcutaneous Electrical Nerve Stimulation is a non-invasive, drug free method of controlling pain. TENS uses tiny electrical impulses sent through the skin to nerves to modify your pain perception. TENS does not cure any physiological problem; it only helps control the pain. TENS does not work for everyone; however, in most patients it is effective in reducing or eliminating the pain, allowing for a return to normal activity.

### **HOW TENS WORKS**

There is nothing "magic" about Transcutaneous Electrical Nerve Stimulation (TENS). TENS is intended to be used to relieve pain. The TENS unit sends comfortable impulses through the skin that stimulate the nerve (or nerves) in the treatment area. In many cases, this stimulation will greatly reduce or eliminate the pain sensation the patient feels. Pain relief varies by individual patient, mode selected for therapy, and the type of pain. In many patients, the reduction or elimination of pain lasts longer than the actual period of stimulation (sometimes as much as three to four times longer). In others, pain is only modified while stimulation actually occurs. You may discuss this with your physician or therapist.

### 1. Precautions:

Isolated cases of skin irritation may occur at the site of electrode placement following long-term application. Effectiveness is highly dependent upon patient selection by a person qualified in the management of pain patients.

### 2. Contradictions:

TENS devices can affect the operation of demand type cardiac pacemakers. TENS is not recommended for patients with known heart disease without physical evaluation of risk. Do not use TENS on the carotid sinus(neck) region. Do not apply TENS for undiagnosed pain syndromes until etiology is established. Do not stimulate on the site that may cause current to flow transcerebrally (through the head).

### 3. Adverse Reactions

Possible allergic to gel, skin irritation and electrode burn are potential adverse reactions.

4. Read operation manual before use of TENS.
5. We emphasize that patient with an implanted electronic device (for example, a pacemaker) should not undergo TENS treatment without first consulting a doctor. The same applies to patients with any metallic implants.
6. If TENS therapy becomes ineffective or unpleasant, stimulation should be discontinued until its use is reevaluated by the physician or therapist.
7. This device should not be used while driving, operating machinery, close to water, or during any activity in which involuntary muscle contractions may put the user at undue risk of injury.
8. Turn the T.E.N.S. off before applying or removing electrodes.
9. T.E.N.S. devices have no AP/APG protection. Do not use it in the presence of explosive atmosphere and flammable mixture.


1. Caution should be used in applying TENS to patients suspected of having heart disease. Further clinical data is needed to show there are no adverse results.
2. The safety of TENS devices for use during pregnancy or birth has not been established. **Do not use TENS during pregnancy.**
3. TENS is not effective for pain of central origin. (This includes headache.)
4. TENS devices should be used only under the continued supervision of a physician.
5. TENS devices have no curative value.
6. TENS is a symptomatic treatment and as such suppresses the sensation of pain which would otherwise serve as a protective mechanism.
7. Electronic monitoring equipment (such as ECG monitors and ECG alarms) may not operate properly when TENS stimulation is in use.
8. There should be a prominently placed statement warning that stimulus delivered by this device may be sufficient to cause electrocution. Electrical current of this magnitude must not flow through the thorax because it may cause a cardiac arrhythmia.
9. Do not place electrodes on the front of the throat as spasm of the Laryngeal and Pharyngeal muscle may occur.
10. Electrodes should not be placed over the eyes, in the mouth, or internally.
11. Keep this device out of the reach of children.
12. Caution: Federal law restricts this device to sale by or on the order of a physician

## CHAPTER 4 : GENERAL DESCRIPTION

The ZZAN602 is a battery operated pulse generator that sends electrical impulses electrodes to the body and reach the nerves causing pain. The device is provided with two controllable output channels, each independent of each other. An electrode pair can be connected to each output channel.

The electronics of the ZZAN602 create electrical impulses whose Intensity, duration, number per second and modulation may be altered with the controls or switches. Dial controls are very easy to use and the slide cover prevents accidental changes in the setting.

## CHAPTER 5 : CONSTRUCTION

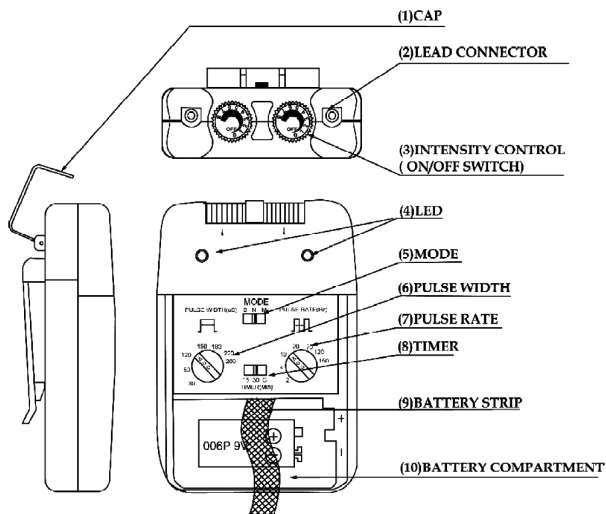

The technical specification details of ZZAN602 are as follows.

	MECHANISM	TECHNICAL DESCRIPTION
1	Channel	Dual, isolated between channels
2	Pulse Amplitude	Adjustable, 0-80 mA peak into 500 ohm load each channel
3	Pulse Rate	Adjustable, from 2 to 150 Hz
4	Pulse Width	Adjustable, from 30 to 260 microseconds
5	Modulation Mode	Pulse rate is automatically varied in a cyclic pattern over an interval of nominally 10 Seconds.(in max 150Hz) Pulse rate decreases linearly over a period of 4 seconds from the control setting value to a value which is 40%less. The lower pulse rate will continue for 1 second. Then increase linearly over a 4 seconds period to its original value. The original pulse rate will continue for 1 second. The cycle is then repeated.
6	Burst Mode	Bursts occur twice very second. Pulse width(adjustable), frequency = 100 Hz
7	Wave Form	Asymmetrical Bi-Phasic Square Pulse
8	Timer	15, 30 minutes or Continue
9	Voltage	0 to 40 V (Load : 500 ohm)
10	Max. Charge per pulse	20 micro-coulombs
11	Power Supply	Approximately 50 hours at nominal settings.
12	Battery Life	Approximately 50 hours at nominal settings.
13	Size	95(H) x 65(W) x 23.5(T) mm
14	Weight	115 grams(battery included)

## CHAPTER 7 : REPLACEABLE PARTS

The replaceable parts and accessories of T.E.N.S. devices are as given below-Except leads, electrodes and battery, battery case cover, please do not try to replace the other parts of a device.

NO.	PARTS
01	ELECTRODES LEADS
02	ELECTRODES
03	9V BATTERY , TYPE 6F22
04	BELT CLIP
05	BATTERY CASE COVER
06	LEAD CONNECTOR
07	MAIN PCB
08	INTENSITY KNOB
09	B-N-M SWITCH
10	PULSE WIDTH KNOB
11	PULSE RATE KNOB

## CHAPTER 8 : ACCESSORIES

Each set ZZAN602 are completed with standard accessories and standard label as given below

Ref. NO.	DESCRIPTION	Q'TY
01	40 X 4 0 mm Adhesive Electrodes	4 pcs.
02	Electrodes Leads	2 pcs.
03	9V Battery, type 6F22	1 pcs.
04	Instruction Manual	1 pcs.
05	Carrying Case	1 pcs.



## CHAPTER 9 : GRAPHIC SYMBOLS



Note Operating Instructions



Degree of Electrical Protection BF



Do not insert the plug into AC power supply socket.



Direct Current (DC power source)

## CHAPTER 10 : PARAMETER CONTROLS

### PULSE DURATION

Wider pulse duration settings will deliver stronger stimulation for any given intensity setting. As mentioned in the Controls section, by using a combination of intensity and pulse duration, it is felt that various pulse widths are capable of stimulating different groups of nerve fibers.

The choice of which pulse duration to use is partially dependent upon the Treatment Mode and Protocol selected (refer to the appropriate section).

### PULSE RATE

The Pulse Rate (hertz or pulses per second) chosen depends greatly upon the type of electrode placement given to the patient.

When using contiguous and dermatome electrode placements (i.e.

stimulating directly through the area of pain or localized enervation), a quick pulse rate (setting greater than 80Hz on the Pulse Rate Control) is desired. The patient should not perceive individual pulses but rather have the sensation of steady continuous stimulation.

Despite above recommendations, these individual patients may require slight variations of the above settings, according to the nature of their condition.

### **TREATMENT MODE**

Normal or Conventional TENS offers the practitioners complete control over all the various treatment parameters of the instrument.

Burst Mode is analogous to the Low Rate TENS technique except the low frequency individual pulses are replaced by individual "bursts" of 7-10 individual pulses. It is thus a combination of Conventional TENS and Low Rate TENS. In Burst Mode, the treatment frequency is fixed by the instrument and is not adjustable with the Frequency Rate control.

Modulated Mode attempts to prevent nerve accommodation by continuously cycling the treatment intensity. When using Modulated Mode increase the intensity only when the unit is at the maximum intensity of the modulation cycle. If the intensity is increased during a low intensity period of the cycle, the patient may turn up the control very slowly, so that they may feel the intensity any higher.


### **TIME DURATION**

The onset of pain relief should occur shortly after the intensity setting has been determined. However, in some cases, pain relief may take as long as 30 minutes to achieve. TENS units are typically operated for long periods of time, with a minimum of 20 – 30 minutes and in some post-operation protocols, as long as 36 hours.

In general, pain relief will diminish within 30 minutes of the cessation of stimulation.

## CHAPTER 11 : ATTACHMENT OF ELECTRODE LEAD WIRES

The wires provided with the system insert into the jack sockets located on top of the device. Holding the insulated portion of the connector, push the plug end of the wire into one of the jacks (see drawing); one or two sets of wires may be used.

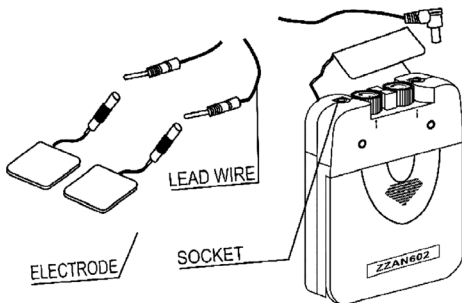

After connecting the wires to the stimulator, attach each wire to an electrode. Use care when you plug and unplug the wires. Jerking the wire instead of holding the insulated connector body may cause wire breakage.

### CAUTION

Do not insert the plug of the patient lead wire into the AC power supply socket.

## CHAPTER 12: LEAD WIRE MAINTENANCE

Clean the wires by wiping with a damp cloth. Coating them lightly with talcum powder will reduce tangling and prolong life.

## CHAPTER 13 : ELECTRODE OPTIONS

The electrodes are disposable and should be routinely replaced when they start to lose their adhesive nature. If you are unsure of your electrode adhesive properties, order replacement electrodes. Replacement electrodes should be re-ordered through or on the advice of your physician to ensure proper quality. Follow application procedures outlined in electrode packing, to maintain optimal stimulation and to prevent skin irritation.

## CHAPTER 14 : ELECTRODE PLACEMENT

The placement of electrodes can be one of the most important parameters in achieving success with TENS therapy. Of utmost importance is the willingness of the clinician to try the various styles of electrode placement to find which method best fits the needs of the individual patient.

Every patient responds to electrical stimulation differently and their needs may vary from the conventional settings suggested here. If the initial results are not positive, feel free to experiment. Once an acceptable placement has been achieved, mark down the electrodes sites and the settings on the patient's reference sheet of this manual, so the patient can easily continue treatment at home.

### **CONTIGUOUS PLACEMENT**

This is the most common placement technique. It involves placing the electrodes alongside the area of localized pain site, in such a way as to direct the flow of current through or around the area of pain.

In a single channel application, this would involve placing each pad on either side of the pain site if the pain is localized on a limb and deep within the tissue. Pad placement on the posterior and anterior aspects of the affected limb will allow the current to flow completely through the limb and thus through the endogenous pain site.

With a two channels application, the clinician may either direct the current flow to cross through the pain site or, in what is called the "bracket" method allowing the current flow on either side of the painful area, generally through the nerve branches that feed into the pain site.

## CHAPTER 15 : TIPS FOR SKIN CARE

**To avoid skin irritation, especially if you have sensitive skin, follow these suggestions:**


1. Wash the area of skin where you will be placing the electrodes, using mild soap and water before applying electrodes, and after taking them off. Be sure to rinse soap off thoroughly and dry skin well.
2. Excess hair may be clipped with scissors; do not shave stimulation area.
3. Wipe the area with the skin preparation your clinician has recommended. Let this dry. Apply electrodes as directed.
4. Many skin problems arise from the "pulling stress" from adhesive patches that are excessively stretched across the skin during application. To prevent this, apply electrodes from center outward; avoid stretching over the skin.
5. To minimize "pulling stress", tape extra lengths of lead wires to the skin in a loop to prevent tugging on electrodes.
6. When removing electrodes, always remove by pulling in the direction of hair growth.
7. It may be helpful to rub skin lotion on electrode placement area when not wearing electrodes.
8. Never apply electrodes over irritated or broken skin.

### Application

1. Clean and dry the skin at the prescribed area thoroughly with soap and water prior to application of electrodes.
2. Insert the lead wire into the pin connector on the pre-wired electrodes.
3. Remove the electrodes from the protective liner and apply the electrodes firmly to the treatment site.

### Removal

1. Lift at the edge of electrodes and peel; do not pull on the lead wires because it may damage the electrodes.
2. Place the electrodes on the liner and remove the lead wire by twisting and pulling at the same time.

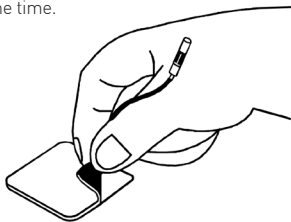

### Care and Storage


1. Between uses, store the electrodes in the resealed bag in a cool dry place.
2. It may be helpful to improve repeated application by spreading a few drops of cold water over the adhesive and turn the surface up to air dry. Over Saturation with water will reduce the adhesive properties.

### Important

1. Do not apply to broken skin.
2. The electrodes should be discarded when they are no longer adhering.
3. The electrodes are intended for single patient use only.
4. If irritation occurs, discontinue use and consult your clinician.
5. Read the instruction for use of self-adhesive electrodes before application.

## 1. Slide Cover:

A slide-on panel cover covers the controls for Pulse Width, Pulse Rate, Mode Selector and Modulation Selector. Your medical professional may wish to set these controls for you and request that you leave the cover in place.

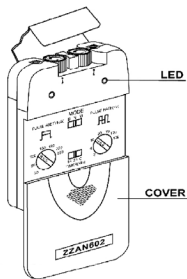

## 2. Display Led


Each of the leds illuminates whenever the electronics of the device create a current impulse. Due to the capacity of the human eye, the illumination of the lamp can only be recognized up to a frequency of approximately 30 Hz. At higher frequencies, the lamp will appear to be constantly illuminated.

## 3. On/Off Switch and Intensity Control :

If both controls are in the off-position (white markings on the housing), the device is switched off. By turning the controls clockwise, the appropriate channel is switched on and the impulse display led will illuminate and begin to pulse according to the frequency set.



The current strength of the impulses transmitted to the electrodes increases further when the control is turned clockwise.

To reduce the current strength or switch the device off, turn the controls counter clockwise to the required setting or off- position.

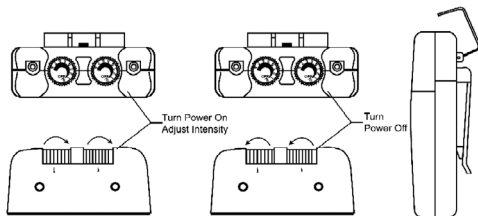

## 5. Mode Control

Expose the controls by sliding front cover down from top of unit. This switch has 3 positions: B for Burst stimulation, N for Constant stimulation, and M for modulation stimulation. Push the Mode Selector until engaged in position desired.

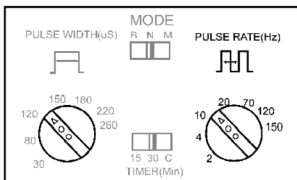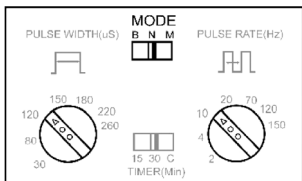


## 6. Pulse Rate Control:


This dial determines how many electrical impulses are applied through the skin each second. By turning these controls, the number of current impulses per second(Hz) for both channels can be continually adjusted. Unless otherwise instructed, turn the pulse rate control to the 70-120 Hz range.

## 7. Pulse Width Control:

This dial adjusts the length of time each electrical signal is applied through the skin, which controls the strength and sensation of the stimulation. If no instructions regarding the pulse width are given in therapy, set the control to the suggested 70-120  $\mu$ s setting.

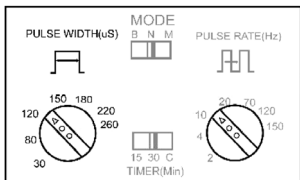



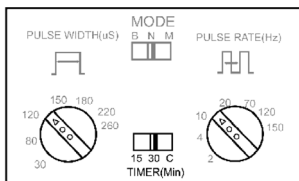


## 8. Timer Control

Treatment time of TENS can be preset with Timer Control. This switch has 3 positions, 15, 30 and C(Continue). . Push the Timer Control until engaged in position desired.

## 9. Check/Replace the Battery:

Over time, in order to ensure the functional safety of TENS, changing the battery is necessary.

1. Make sure that both intensity controls are switched to off position.
2. Slide the battery compartment cover and remove.
3. Remove the battery from the compartment.
4. Insert the battery into the compartment. Note the polarity indicated on the battery and in the compartment.
5. Replace the battery compartment cover and slide to close

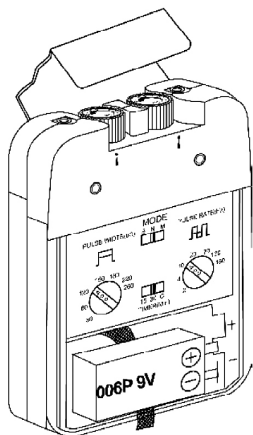

ZZAN602 can be used with a rechargeable battery when necessary. If you use rechargeable batteries, please follow the instructions.

### **RECHARGEABLE BATTERIES (NOT INCLUDED)**

Prior to the use of a new unit, the rechargeable battery should be charged according to the battery manufacturer's instructions. Before using the battery charger, read all instructions and cautionary markings on the battery and in this instruction manual. After being stored for 60 days or more, the batteries may lose their charge. After long periods of storage, batteries should be charged prior to use.

### **BATTERY CHARGING**

- (1) Plug the charger into any working 110 or 220/240v mains electrical outlet. The use of any attachment not supplied with the charger may result in the risk of fire, electric shock, or injury to persons.
- (2) Follow the battery manufacturer's instructions for charging time.
- (3) After the battery manufacturer's recommended charging time has been completed, unplug the charger and remove the battery.
- (4) Batteries should always be stored in a fully charged state. To ensure optimum battery performance, follow these guidelines:
  - (a) Although overcharging the batteries for up to 24 hours will not damage them, repeated overcharging may decrease useful battery life.
  - (b) Always store batteries in their charged condition. After a battery has been discharged, recharge it as soon as possible. If the battery is stored more than 60 days, it may need to be recharged.
  - (c) Do not short the terminals of the battery. This will cause the battery to get hot and can cause permanent damage. Avoid storing the batteries in your pocket or purse where the terminals may accidentally come into contact with coins, keys or any metal objects.

### **(d) WARNINGS:**

1. Do not attempt to charge any other types of batteries in your charger, other than the nickel-cadmium rechargeable batteries. Other types of batteries may leak or burst.
2. Do not incinerate the rechargeable battery as it may explode!

## CHAPTER 19: MAINTENANCE, TRANSPORTATION AND STORAGE OF TENS DEVICE

1. Non-flammable cleaning solution is suitable for cleaning the device. Note: Do not smoke or work with open lights (for example, candles, etc.) when working with flammable liquids.
2. Stains and spots can be removed with a cleaning agent.
3. Do not submerge the device in liquids or expose it to large amounts of water.
4. Return the device to the carrying box with sponge foam to ensure that the unit is well-protected before transportation.
5. If the device is not to be used for a long period of time, remove the batteries from the battery compartment (acid may leak from used batteries and damage the device). Put the device and accessories in carrying box and keep it in cool dry place.
6. The packed TENS device should be stored and transported under the temperature range of  $-20^{\circ}$  ~  $+60^{\circ}$ , relative humidity 20% ~ 95%, atmosphere pressure 500 hPa ~ 1060 hPa.

## CHAPTER 20: SAFETY-TECHNICAL CONTROLS

**For safety reasons, check your ZZAN602 each week based on the following checklist.**

1. Check the device for external damage.
  - deformation of the housing.
  - damaged or defective output sockets.
2. Check the device for defective operating elements.
  - legibility of inscriptions and labels.
  - make sure the inscriptions and labels are not distorted.
3. Check Led
  - led must be illuminated when switched on.
4. Check the usability of accessories.
  - patient cable undamaged.
  - electrodes undamaged.

Please consult your distributor if there are any problems with device and accessories.

## CHAPTER 21 : MALFUNCTIONS

### **Should any malfunctions occur while using the TENS, check**

- whether the switch/control is set to the appropriate form of therapy. Adjust the control correctly.
  - whether the cable is correctly connected to the device. The cables should be inserted completely into the sockets.
  - whether the impulse display led is illuminated. If necessary, insert a new battery.
  - for possible damage to the cable. Change the cable if any damage is detected.
- \* If there is any other problem, please return the device to your distributor. Do not try to repair a defective device.

## CHAPTER 22: CONFORMITY TO SAFETY STANDARDS

The ZZAN602 devices are in compliance with the EN 60601-1-2:2001 and EN 60601-1:1990+A1:1993+A2:1995+A13: 1996 safety standards.

## CHAPTER 23 : WARRANTY

All ZZAN602 models carry a warranty of one year from the date of delivery. The warranty applies to the stimulator only and covers both parts and labor relating thereto. The warranty does not apply to damage resulting from failure to follow the operating instructions, accidents, abuse, alteration or disassembly by unauthorized personnel.



Manufactured for BodyMed®

Hudson, Ohio 44236 1-866-528-2152

Made in Taiwan. Rev. 082410