

## Sterile SatPax® 1000

Sterile Pre-wetted 55% Cellulose / 45% Polyester Nonwoven Cleanroom

Sterile SatPax® 1000 combines Durx® 670 nonwoven wipers composed of a hydroentangled nonwoven blend of 55% cellulose and 45% polyester with a high saturation level of 70% IPA and 30% DI water.

This pre-wetted format provides a cost effective and easy to use solution versus traditional bulk handling of solvents, maintenance of squirt bottles and inconsistent wetting and cleaning associated with wetting a dry wiper.

Each package is then gamma irradiated sterile to a 10<sup>-6</sup> Sterility Assurance Level making the product ideal for cleaning critical sterile processing environments.



### Key Attributes

- 55% cellulose / 45% polyester hydroentangled nonwoven blend
- No chemical binders in base material
- Pre-wetted with consistent 70% IPA/ 30% DI Water to a high saturation level
- Re-sealable solvent resistant packaging
- Wipers in c-folded configuration for single withdrawal
- Gamma irradiated and sterile validated to a 10<sup>-6</sup> Sterility Assurance Level per AAMI guidelines
- Certificate of Sterility available to download on the web by lot number 24/7. Details: Date of Manufacture, Radiation Dosage, Lot Information, Meets cGMP requirements for traceability

### Benefits

- Low extractables and fiber and particle counts
- Smooth and durable with good wet strength
- Reduces alcohol usage and preparation / handling costs
- Reduces VOC emissions
- Increases cleaning efficiency
- Increases cleaning protocol consistency

### Other Sterile Pre-wetted wipers

- Sterile SatPax® 670
- Sterile SatPax® 1200
- Sterile SatPax® -HA

### Environmental

- Reduces VOC (Volatile Organic Compound) emissions
- Complies with RoHS and pre-registered under REACH

### Applications

- Designed for use in ISO Class 5 and higher sterile cleanroom environments
- Designed for use in wet cleaning of critical surfaces where control of flammable solvents and flammable solvent concentrations is required
- Final cleaning of surfaces or products prior to manufacturing or packaging
- High saturation level is ideal for removing cleaning and disinfecting residues in sterile environments

### Saturation Levels

The amount of solution contained in each wiper can be varied according to customer requirements. Higher saturation levels apply more solution to the surface during cleaning.

### Validated Sterile

ANSI/AAMI/ISO procedures have been followed for verification of irradiation dose with audits being performed to monitor bioburden levels and dose. The Sterility Assurance Level (SAL) or the probability of a viable organism being present on a product unit after sterilization is a SAL of 10<sup>-6</sup> (which is the chance of 1 non-sterile unit in 1,000,000).

[www.berkshire.eu.com](http://www.berkshire.eu.com)

Contact: Tel 44 (0) 870 757 2877

[enquiries@berkshire.uk.com](mailto:enquiries@berkshire.uk.com)

Americas	Tel 1 413 528 2602	Fax 1 413 528 2614	<a href="mailto:info@berkshire.com">info@berkshire.com</a>
Europe	Tel 44 (0) 870 757 2877	Fax 44 (0) 870 757 2878	<a href="mailto:enquiries@berkshire.uk.com">enquiries@berkshire.uk.com</a>
SE Asia	Tel 65 6252 4313	Fax 65 6252 4312	<a href="mailto:enquiries@berkshire.com.sg">enquiries@berkshire.com.sg</a>
Japan	Tel 81 3 5827 2380	Fax 81 3 5827 2382	<a href="mailto:master@berkshire.co.jp">master@berkshire.co.jp</a>

## Technical Data (In Dry State)

Attribute		Units	Value	Test Method
<b>Basis Weight</b>		g/m <sup>2</sup>	<b>68.0</b>	TAPPI T-410
<b>Caliper</b>		µm	<b>264</b>	TAPPI T-411
<b>Fibers</b>	≥100µm	fibers/cm <sup>2</sup>	<b>160</b>	IEST-RP.CC004.3, Sec 6.1.3 / Sec 6.2.2
<b>Particles</b>	≥0.5µm	x10 <sup>9</sup> /cm <sup>2</sup>	<b>10</b>	IEST-RP.CC004.3, Sec 6.1.3 / Sec 6.2.1
<b>Sorbency</b>	Capacity	mL/m <sup>2</sup>	<b>320</b>	IEST-RP.CC004.3, Sec 8.1 modified / Sec 8.2 modified
	Efficiency	mL/g	<b>4.7</b>	
	Rate	seconds	<b>2</b>	
<b>Non-Volatile Residue</b>	DI Water	g/m <sup>2</sup>	<b>0.028</b>	IEST-RP.CC004.3, Sec 7.1.2
	IPA	g/m <sup>2</sup>	<b>0.0038</b>	
<b>Ions</b>	Na <sup>+</sup>	ppm	<b>62</b>	IEST-RP.CC004.3, Sec 7.2.2
	K <sup>+</sup>	ppm	<b>5.9</b>	
	Ca <sup>++</sup>	ppm	<b>22</b>	
	Mg <sup>++</sup>	ppm	<b>5.0</b>	
	Cl <sup>-</sup>	ppm	<b>31</b>	

### Notes:

- Technical data represented in this table are typical values at the time of publication. These should not be used as product specifications.
- Due to differences in test methods applied and equipment utilized by different wiper manufacturers, valid product comparisons may only be obtained through side-by-side testing in the same test facility, under similar conditions
- Third party testing can be performed upon request

### Order Information:

Product	Number	Size	Shts/pk	Pks/cs	IPA/DI	Saturation	VOC % by	Style
Sterile SatPax® 1000	SSP1000.001.24	9x9" (23x23cm)	30	24	70/30	70%	46.2%	C-fold

### Other Berkshire products




**Wipers**



**Glove Liners**



**Mop Systems**

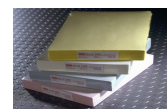

**Documentation Systems**



**Face Masks**



**Swabs**