

TIVA2 Washers

Thermal High Disinfector Line

Dental and medical practices need strict cleaning, disinfection, and sterilization of instruments before they can be used on patients. Thermal High Disinfection employs moist heat above 194° for a minimum of 5 minutes to kill bacteria and viruses in the washing phase. The Tiva line provides:

- Thermal High Disinfection to achieve a complete decontamination cycle.
- Minimizes consumption and operation costs while maximizing sustainability
- Improves user experience with an advanced control system and an easy to use interface



TIVA2-H TD
Under Counter Washer



TIVA2 TD
Counter Top Washer



High end washers for the decontamination of complex instruments

TIVA2 Thermal High Disinfector washers provide a flexible turnkey solution for washing, disinfecting and drying diverse medical instruments for a complete decontamination cycle.



Pre-wash
Cold water removes soiled material



Detergent wash
Cleaning with warm demineralized water & suitable alkaline detergents



1st Rinse
Rinse with warm water to remove all cleaning agent residues, second rinse when required



Thermal disinfection
Water at a temperature above 194oF for a minimum of 5 minutes



Hot air drying
The dried instruments are ready for the next sterilization phase

TIVA2/TIVA2-H TD are fully-automatic thermal high disinfectors for all practices.

Features:

- Touch panel, color graphic display and real-time monitoring
- A total of 40 programs: Preset and programmable cycles
- Air drying prefilter + HEPA H14 filter
- Integrated water softener
- 2 dosing pumps with flowmeter for alkaline detergent and neutralizer (additional pump optional for TIVA2-H TD)
- USB/Printer for data documentation
- Provides complete & total evacuation of previous water (ETS)


Utility Requirements:

Water Connections

- Cold water connection for TIVA2 (default)
- Connection for demineralized water
- Integrated Water Softener with Salt Box

Electrical Connections

- 208+240V 1 phase 60Hz

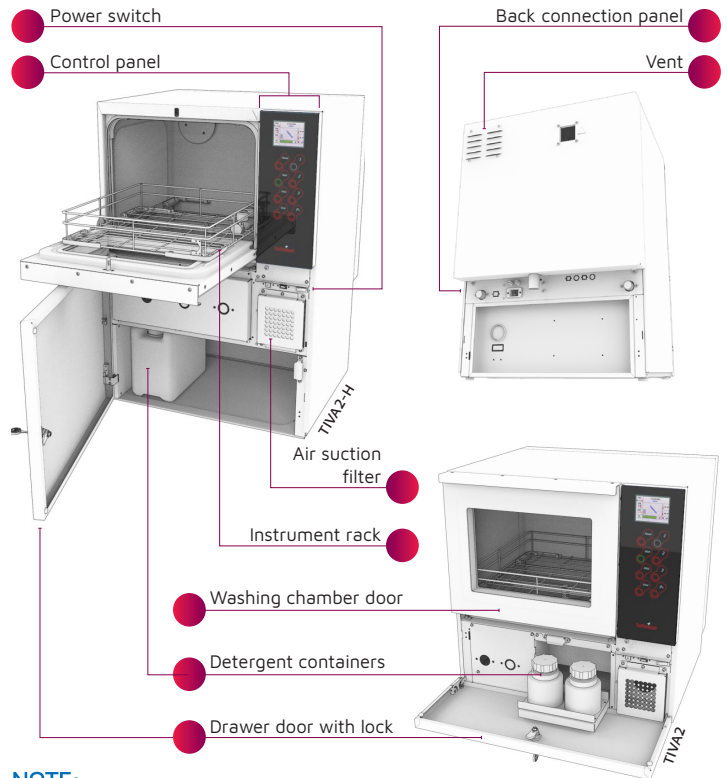


NOTE:


The front and back panel of the TIVA2 and TIVA2-H are identical

Comes with the Following Configurations:

Dental:

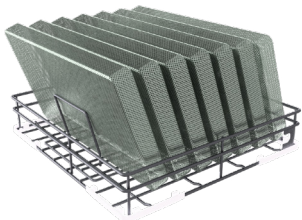

Medical:

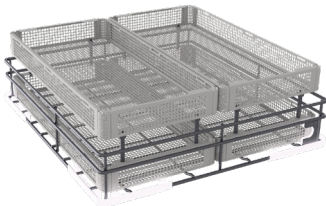

*Unit does not come with any cassettes

Dimensions:

TIVA2-H TD Under Counter Washer

- (WxDxH) 23.4"x21.3"x32.2"
- Usable Chamber Volume: 65L

TIVA2 TD Counter Top Washer

- (WxDxH) 23.4"X21.3"X23.6"
- Usable Chamber Volume: 65L

T-Clean Tuttnauer Line of Detergents & Neutralizers

-User friendly for use in your Washer

T-Clean TIVA Detergent



T-Clean TIVA Neutralizer:



T-Clean Prep Foam
PF0006
500 ml with sprayer



Indicators for Cleaning Efficiency Testing

Cleaning Indicator (High Level)



WTL198-0062

Cleaning Indicator Holder for Washer



WTL198-0065



Learn from our experts.
Visit our blog: tuttnauer.com/blog



Tuttnauer USA Co. Ltd.
25 Power Drive, Hauppauge, NY 11788
Tel: +800 624 5836, +631 737 4850, Fax: +631 737 0720
E-mail: info@tuttnauerUSA.com, www.tuttnauerUSA.com

