

# 1<sup>st</sup> Brushless ENDO MOTOR

by Woodpecker

## MotoPex

Super mini head & Slim neck  
Stronger power  
Precise rotation  
Lower noise  
Balances Weight




GUILIN WOODPECKER MEDICAL INSTRUMENT CO., LTD.




## High performance brushless motor

- Adjustable reciprocating angle  
(10° interval; adjustment range: 20° - 340° );  
Accurate rotating angle reduces risk of instrument fracture and enables higher preparing efficiency.

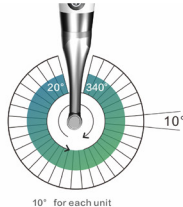

- Smooth rotation and lower noise make treatment more endurable for patients.

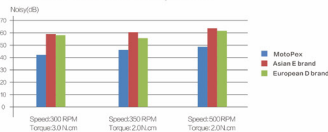

- Strong power, together with accurate and stable torque output, helps to realize efficient root canal preparation.




Built-in 30 mainstream file systems.



2000mAh Lithium battery


## Super Mini Contra-angle

Super mini head with a diameter of 8mm and a height of 9.7mm brings clear vision during operation.

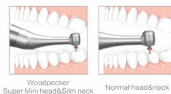

## Built-in Apex locator

Integrated length determination makes endodontic treatment safer and more efficient. Realize real-time display of file location in canal. Automatically stops or reversely rotates as soon as the file reaches the apical stop, so as to prevent perforation. Auto Start and Auto Stop functions respectively enable motor to automatically start while the file gets in canal and stop when it is out.

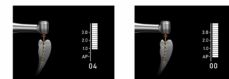

Front center-of-gravity



360° rotatable contra-angle



6:1 Contra-angle



Three times longer service life



The choice of millions of dentists worldwide

Guilin Woodpecker Medical Instrument Co., Ltd.  
Information Industrial Park, Guilin National High-Tech Zone, Guilin, Guangxi, 541004 P. R. China

Tel: Europe Sales Dept.: +86-773-5873196, +86-773-2125222  
North America, South America & Oceania Sales Dept.: +86-773-5873198, +86-773-2125123  
Asia & Africa Sales Dept.: +86-773-5855350, +86-773-2125896  
Fax: +86-773-6822450  
E-mail: woodpecker@glwoodpecker.com, sales@glwoodpecker.com  
Website: http://www.glwoodpecker.com

ECREP MedNet GmbH  
Borkstrasse 10 · 48163 Münster · Germany