



Product Code

Description	Size	Code
N5B Nasal Mask	Extra large	10020336
N5B Nasal Mask	Large	10020319
N5B Nasal Mask	Medium	10020320
N5B Nasal Mask	Small	10020321
N5B Nasal Mask	All Size Kit (L, M & S)	10020322

Description	Size	Code
Cushion Assembly	Extra large	60380390
Cushion Assembly	Large	60380242
Cushion Assembly	Medium	60380243
Cushion Assembly	Small	60380244

Description	Size	Code
Headgear Assembly	/	61030060

Technical Specifications

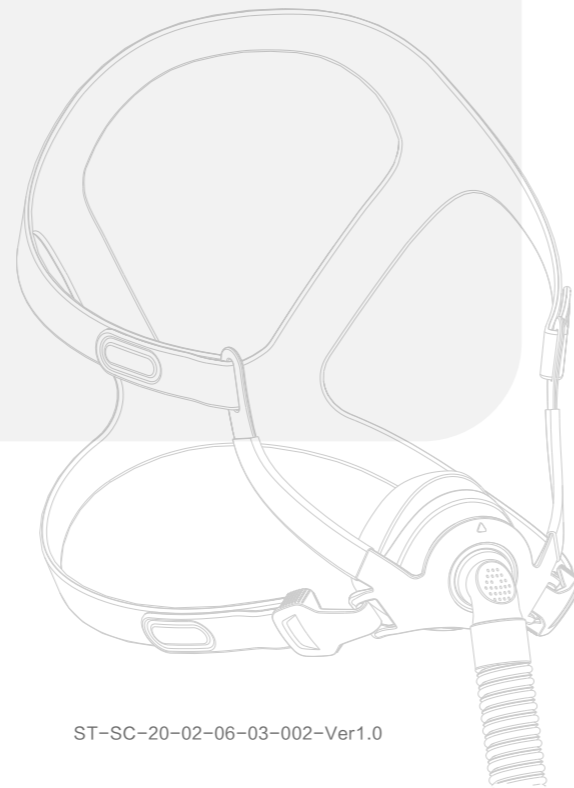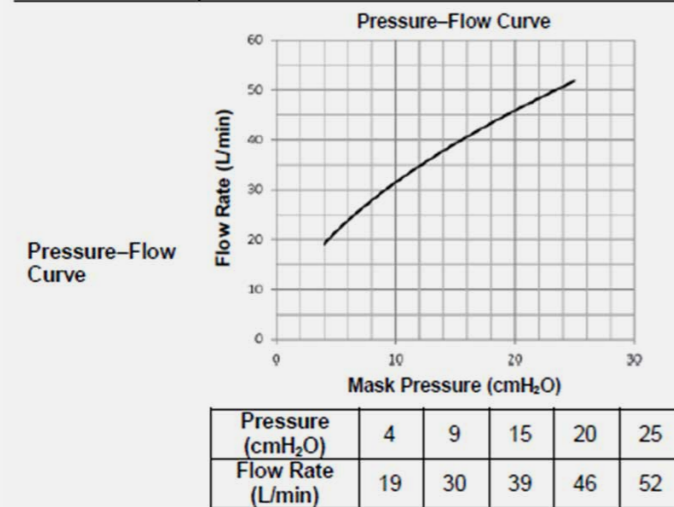

Widened vision, improved seal

The design without forehead frame is suitable for the patient who suffers from the nasal septal deviation. Much comfort plus dynamic seal makes the therapy easier.



family cares

N5B Nasal Mask



BMC Medical Co., Ltd.

Room 110, Tower A Fengyu Building, No.115 Fucheng Road, Haidian, 100036 Beijing, PEOPLE'S REPUBLIC OF CHINA

Tel: +86-10-5166 3880 / Fax: +86-10-5166 3880 ext. 810

en.bmc-medical.com


N5B Nasal Mask

Minimal contact

Thanks to the minimal contact design, less pressure is imposed on the face.

Real comfort

The groove design on the cushion plus the ball and socket elbow cares about your feelings.

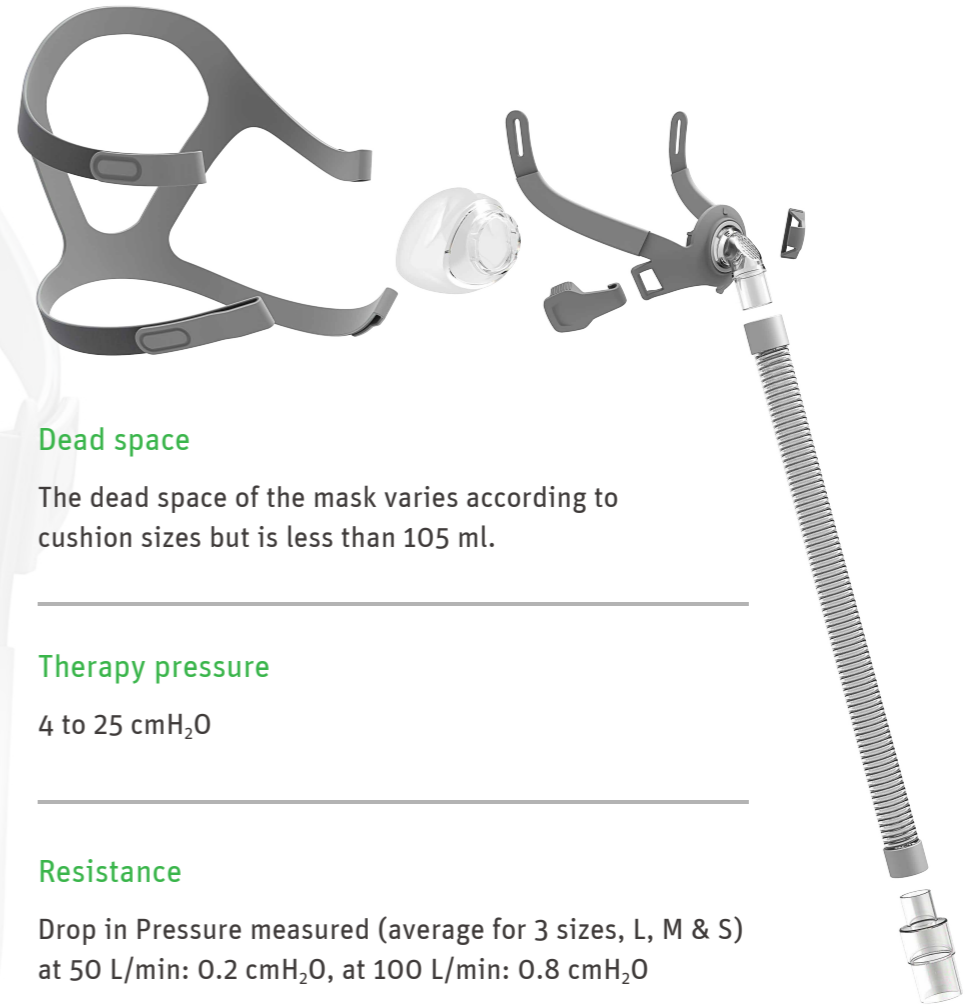

Dead space

The dead space of the mask varies according to cushion sizes but is less than 105 ml.

Therapy pressure

4 to 25 cmH₂O

Resistance

Drop in Pressure measured (average for 3 sizes, L, M & S) at 50 L/min: 0.2 cmH₂O, at 100 L/min: 0.8 cmH₂O

Sound

DECLARED DUAL-NUMBER NOISE EMISSION VALUES in accordance with ISO 4871. The A-weighted sound power level of the mask is 28 dBA, with uncertainty 3 dBA. The A-weighted sound pressure level of the mask at a distance of 1 m is 20 dBA, with uncertainty 3 dBA.



Platform-based frame

The mask frame is compatible with different cushion sizes by one-click lock.

Crown-style headgear

The crown-style headgear provides the user with more ease of use.