

# RHYTHM MULTI

Provides you, the caregiver, with peace of mind, and your patient with world class care and comfort




**HAPPY FEEL™** Heel Relief  
Patent Pending

**PRIUS HEALTHCARE USA**  
*Creating innovation and difference*

[priushcusa.com](http://priushcusa.com)

Perfect for patients and caregivers

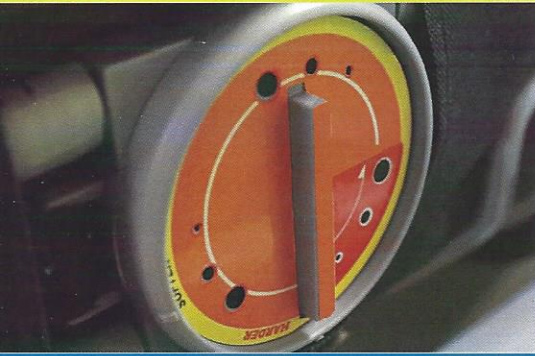


**HAPPYHEEL™** provides extra comfort settings for the heels



Low slip nylon under mattress base cover to reduce shearing and friction between mattress and bed frame



Nylon woven fabric with PU coating finish

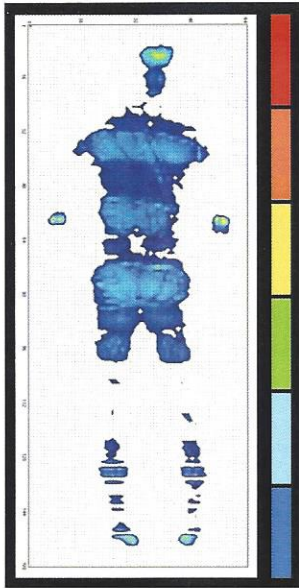

## PRESSURE MAPPING

Rhythm Multi provides constant low pressure by redistributing the patient's pressure over a greater area than a standard mattress.




Foam Safety Mat to prevent bottoming out in case of power loss



10 adjustable comfort settings provides different levels of comfort depending on the patient's therapy required

3 different modes of clinical therapy: static/pulsation/alternation

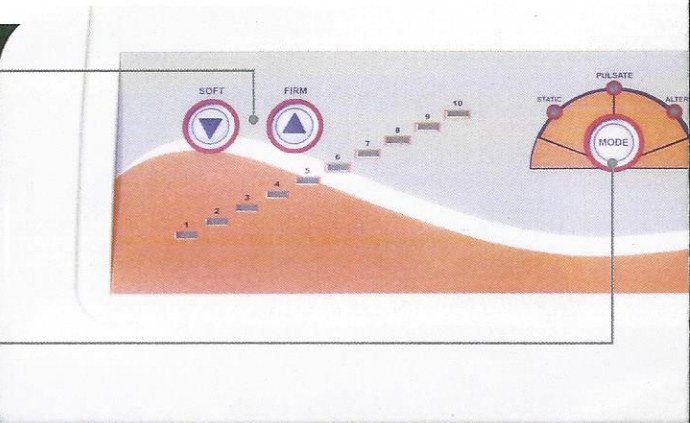


# RHYTHM MULTI

- The Rhythm Multi replacement mattress can accommodate patient weights up to 1000 lbs
- 3-in-1 alternating pressure covers 66% of patient's body to provide a more comfortable yet still effective therapy
- 60 seconds of pulsation simulates a massage to assist in maximizing a patient's comfort
- Modular structure design ensures ease of use and cleaning
  - True low air loss provides extra comfort for moisture and temperature control



- provides 3 selectable softness settings in the last 4 air cells to provide protection for patient's fragile heels

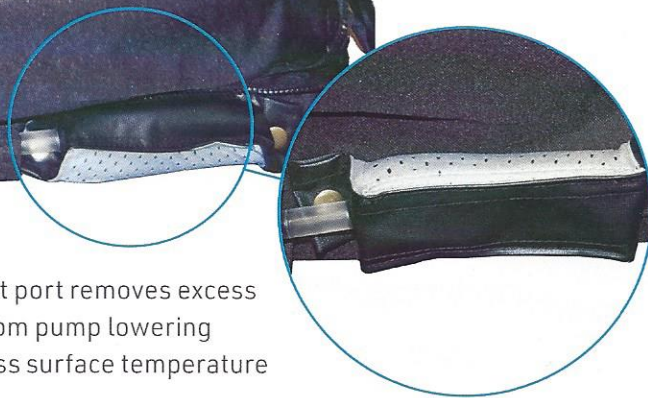

Exhaust port removes excess heat from pump lowering mattress surface temperature


Rhythm Multi is your partner in protection and treatment of pressure ulcers




Cable management keeps electrical cables off the floor



CPR quick release system deflates the system within 20 seconds

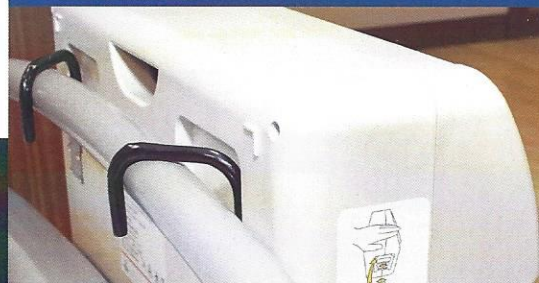

Hanging bracket fits on all type of beds

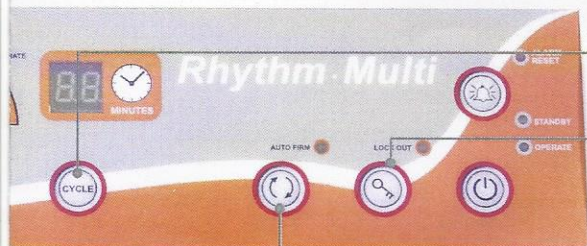

Alternation cycle times from 3-95 minutes


Control Panel Lock Out

Auto Firm/Patient Care control



# RHYTHM MULTI

**Rhythm Multi** Mattress System is a medium to high risk system for patients who are at risk of developing pressure ulcers according to your sound clinical judgment. Rhythm Multi can also be used for patients who have an existing stage 1, 2, 3 and 4 pressure ulcers, in conjunction with your hospital's policy on pressure area management.

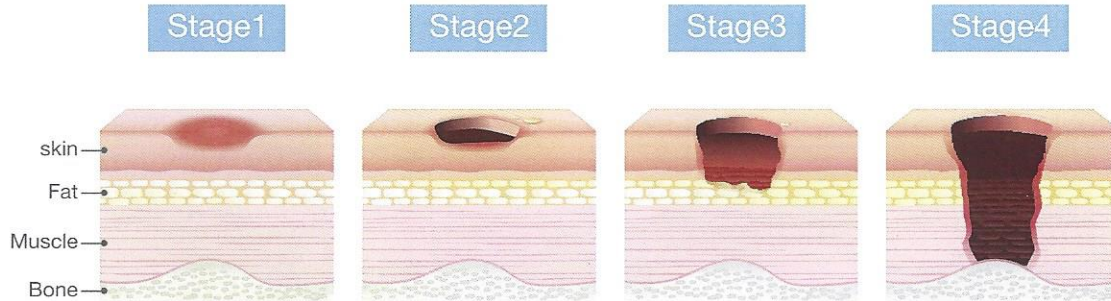

Optional Raised Bolster Cover available for added patient safety. Bolster cover has center access providing easy ingress and egress.

## MASTER CONTROL UNIT

<b>Dimension:</b>	(L)17.7"x (W)6.8"x (H)10.8"
<b>Weight:</b>	13.2 lbs
<b>Rating Voltage:</b>	110-120 volts
<b>Max Current:</b>	5A
<b>Rating Frequency:</b>	60Hz
<b>Cycle Time:</b>	3~95 mins 3 in 1 alternation
<b>Min/Max Pressure:</b>	10~61 mmHg +/- 3~4 mmHg
<b>Protection Type:</b>	Class I, Type BF Power outage visual and audio alarm Static / Pulsation / Alternating Pressure Auto-firm for fast mattress inflation Keypad lock-out function
<b>Air flow rate:</b>	> 1200 liters per minute

## MATTRESS REPLACEMENT

<b>Model No:</b>	Rhythm Multi
<b>Dimensions:</b>	(L)80"x (W)36"x (H)10"
<b>Bariatric Size:</b>	42, 48 inch width
<b>Weight:</b>	28.6 lbs
<b>Weight Capacity:</b>	600, 800, 1000 lbs
<b>Cell Material:</b>	Black Nylon with PU coating
<b>Cell Number:</b>	18 cells
<b>Cover Material:</b>	Blue Nylon with PU coating finish
<b>Base Material:</b>	Black Polyester with PVC coating CPR Quick Release System Cable Management System on bottom cover
<b>HCPSC Code</b>	E0277