

## Product Description

Pure Hold Sanitising Hand Gel – disposable bottle for use in the Pure Hold Hygiene Handle

## Direction for Use

Only to be used with the Pure Hold Hygiene Handle

Suitable for hand disinfectants and surgical hand disinfectants

Apply to hands, rub in, ensuring contact with fingernails and in between fingers. Do not use on damaged, sensitive or irritated skin

## Replacement Gel Bottles

Pure Hold Sanitising Hand Gel comes in a 850ml disposable bottle

PH-GL-001 — 850ml Pure Hold Sanitising Hand Gel

PH-GL-009 — Pack of 9 x 850ml Pure Hold Sanitising Hand Gel

## Specifications

- **Active Ingredients:** Ethanol, Propanol 1, Isopropanol
- **Alcohol Content:** 70%
- **Appearance:** Clear gel
- **Odour:** Alcoholic
- **pH Value:** 5.5

## Benefits

- Independently proven effective against bacteria and viruses
- Dermatologically tested to be kind on skin
- pH of 5.5 to balance and be complementary to the natural pH of skin
- Suitable for use in medical, food and industrial environments
- Leaves hands feeling clean and refreshed, with no sticky residue

## Microbiological Testing:

Pure Hold Sanitising Hand Gel has been tested to the highest European standards, ensuring it was subject to the most stringent and comprehensive testing regimes. The microbiological tests selected covered those organisms which typically represent the highest risk elements in the medical, food and industrial environments

Bactericidal according to:

- BS EN1500
- BS EN12054
- BS EN13727

Virucidal according to:

- BS EN14476

Dermatologically tested - proven to be kind of skin and not cause irritation

## Recommended Shelf Life

36 months

## Order Information


To order Pure Hold Hygiene Handles and/or Pure Hold Sanitising Hand Gel, or to arrange for a trial, simply call us on: +44 (0)2392 299100, e-mail us:

[info@purehold.co.uk](mailto:info@purehold.co.uk) or visit our website [www.purehold.co.uk](http://www.purehold.co.uk)




## Gel Bottle Replacement Instructions


The disposable gel bottles are extremely fast and simple to replace in the Pure Hold Hygiene Handle by following the simple instructions below:

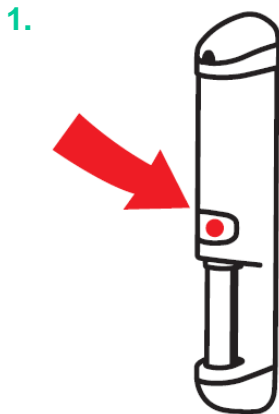


When the red flotation ball becomes fully visible in the reservoir window it indicates that the gel bottle is empty and that it needs to be replaced

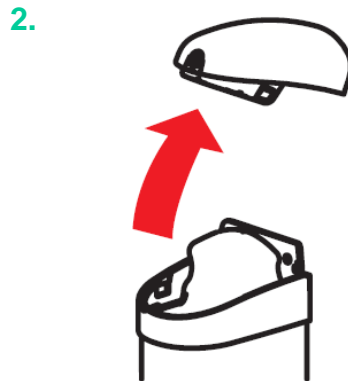


Insert the Pure Hold allen key into the handle lid and removed the lid to reveal the empty bottle

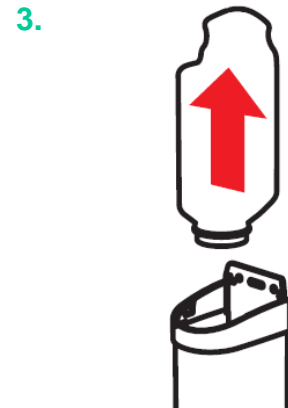

Remove the empty bottle

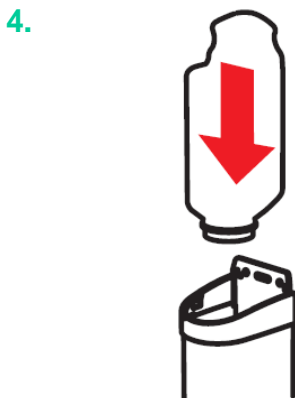

Replace with a new bottle, ensuring bottle clicks into place

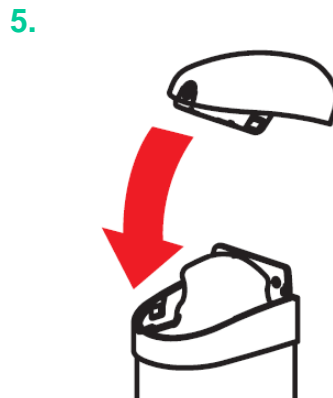

Replace the handle lid by 'snapping' back into position, checking that it is securely locked



Recycle empty bottle

**IMPORTANT NOTE** : The disposable gel bottle unit is not designed to be refilled. Attempts to refill the disposable gel bottle may cause the Hygiene Handle to malfunction and will void the warranty of the Hygiene Handle

## Additional Information

For full safety information see Pure Hold Sanitising Hand Gel Safety Data Sheet (SDS) – this gives guidance on how to safely use the product and also any first aid measures or environmental disposal measures in the event of an accident

For full microbiological testing information see the Pure Hold Sanitising Hand Gel independent tests and performance report pack

Pure Hold Sanitising Hand Gel is also available in a portable, personal-use tottle (PH-GL-050 — 50ml Tottle)