





STAYERの由来





STAYER literally translates to "long-distance running horse". You can use it for a long time with attachment (running), powerful (durable) I named it with this image. A horse is an auspicious animal that is depicted on votive tablets. The logo is colored with a seat color with a riding ground as a motif.

—Wheelchair that can coexist for a long time, STAYER—



STAYERとは



Improved usability





Improvement of Durability and maintainability

1) Ease of use






②Durability /Maintenance



Assistance brake lever and interlocking brake are new mechanisms Reduced wire trouble (improved durability)

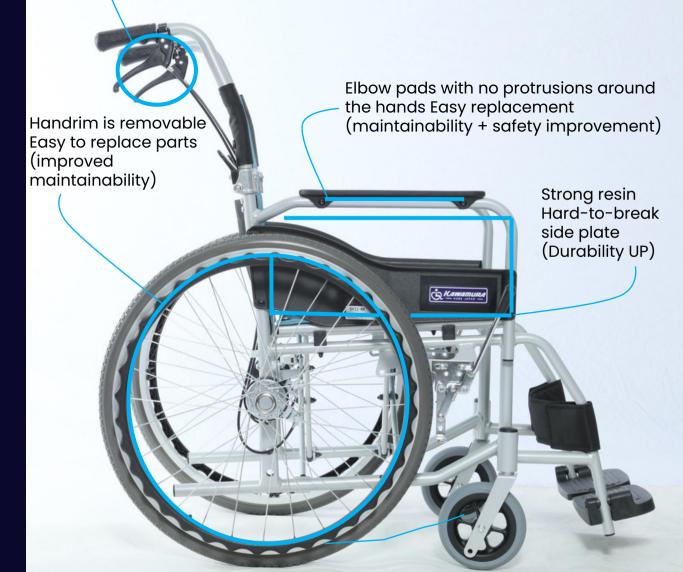






- I want a product that can be disinfected with 70% alcohol (ethanol) and 50ppm hypochlorous water. Or, please let us know which existing products can be used.
- Corona countermeasure products We would like to propose wheelchairs that take measures to prevent infection. "There was a voice." Specifically, a wheelchair that can handle hypochlorous acid.
- Wheelchair or preventive goods for preventing corona infection "Are there any products that specialize in prevention?" Or sell infection control goods "Do you?" I have an inquiry like this.



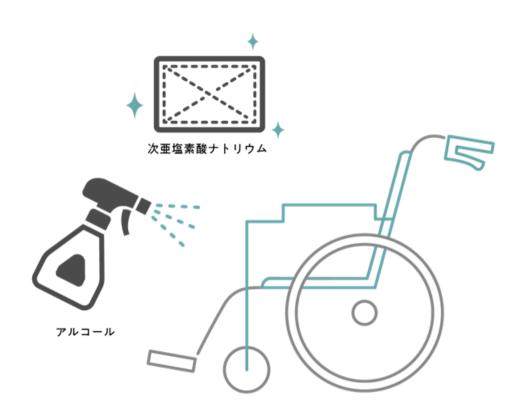

STAYERの特徴


The two new sheet colours are easy to care for. It is generally said that the alcohol concentration at which the sterilization effect is exhibited is 60 to 70%, but mint has been confirmed to be resistant to 80% alcohol. It can also be cleaned with diluted sodium hypochlorite for hygienic cleaning.

*Regarding the handling of chemicals, please refer to the product Please follow the instructions on the package and website.



• 寸法

SY22-45

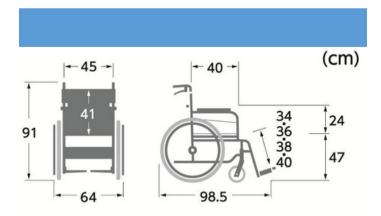


*各寸法には、シート・シートクッション等の厚みは含まれません。

Folded dimensions (cm) ...

Overall Length: 95.5 x Width 32 x Height 71 Wheels: 22in high polymer (no puncture)


Weight: 13.3kg

Sheet colour selection

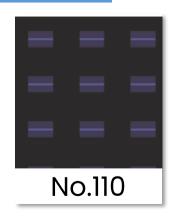



SY16-45

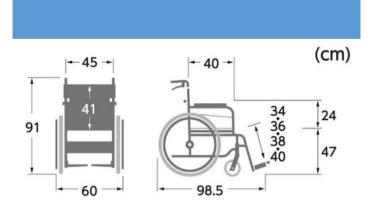

*各寸法には、シート・シートクッション等の厚みは含まれません。

Folding dimensions (cm)...

Length: 86.5 x Width 28 x Height: 71

Wheels: 16in high polymer (no puncture)

Weight: 12.2kg



Flbow Pads

Hand Rim

Foot Step

Interlocking Brake

End Cap

Side Plate



	BM-S	SY	
image Front seat Height	43cm	47cm	
Pusher Height	86.5cml	91.5cm	
Tipping Bar	11cm	14cm	
Step Height Adjustment	3 level	4 level	
Aspect		Adv	vant



The seat height is high, so even tall people can ride.

The height of the pusher makes it easier for a careaiver to push.



The tipping bar becomes longer, allowing you to go over bumps smoothly.

The adjustment range is expanded, and the comfort of sitting is improved.

1-S
۷

When rowing the wheelchair, the pipe around the shoulder blade is angled, making it easy to row. **Back Pipe**

> Since it is fixed from the side, the bolt does not come into contact with the hand and is safe. It is also easy to replace parts.

Removable and easy-to-replace parts. (For rental)

The gap between the steps must be within 2 cm for wide steps. No need to worry about your feet falling.

The assistance brake lever and interlocking brake mechanism are new, reducing wire trouble.


It has a drainage hole and is easy to dry when washed. (For rental)

A side plate that is hard to break using strong resin.

Back seat design for easy parts replacement. (For rental) **Back Seat**

予告NEW MODEL



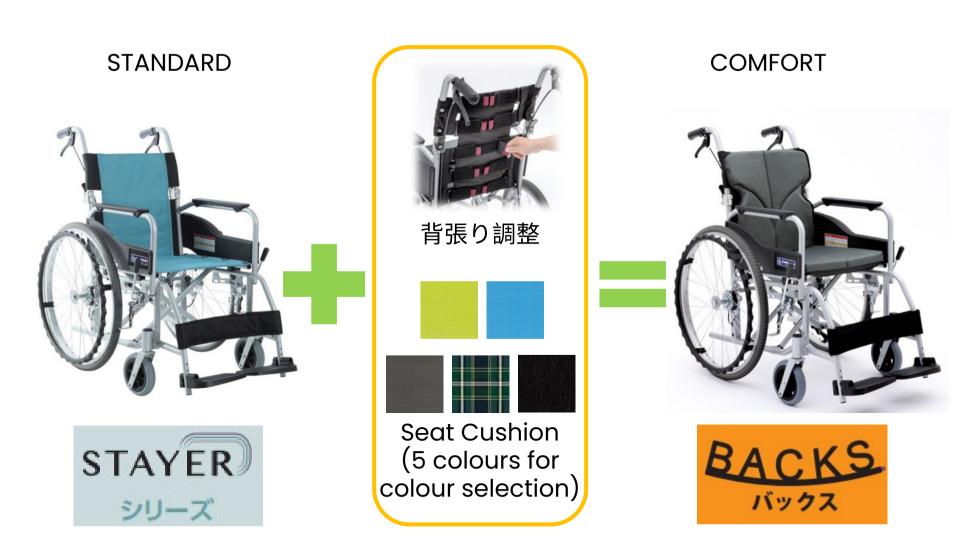

STAYER - Based on this, we will also develop BACKS with backrest adjustment and seat cushion added.