


FreeWheel Rigid Frame Wheelchair

FreeWheel Compatibility Chart

The following rigid frame wheelchairs will work with the FreeWheel Wheelchair Attachment with footrest types similar to those noted in Picture A to the right. —————>

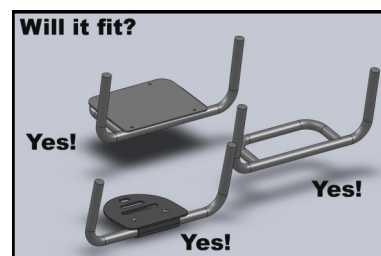

Picture A

Manufacturer	Wheelchair Model	Notes
TiLite	TR-S Series 3, TRA, Aero T, ZR, Aero R	Angle adjustable footrests will require a modified clamp and pinning screws.
	Twist	With non-moving footrest. Some adjustment will be required to footrest mounting.
Sunrise Medical	Quickie 7 series, 5R, GT, GP series, Quickie Helium, Argon 2, Quickiet Life R	With the GT any spacers between the clamp and the footrest will need to be removed. The Q7 will require a custom perch clamp to fit the oval back bar for the FreeWheel storage perch. Argon 2 must have a footrest that is flush with the front tubing (no overhang)
Ki Mobility	Rogue, tsunami ALX	The Rogue requires a custom perch clamp to fit the larger back bar for the FreeWheel storage perch.
	Little Wave Klik, Little Way	Klik footrest must be flush with front tubing. Little Wave will need need a non-moving foorest.
InvaCare	Top End Reveal, terminator, Crossfire T7A, A-4	
Colours	Razor Blade, Eclipse, Shock Blad, Spazz, Saber, Spazz G, Zepher	On Zepher, or other modes with a footrest depth over 6-3/8" may need a custom rear end or some modification to footrest.
	Chump, Razor Blade Jr., Little Dipper, Chump G	Footrests over 4-3/4" will require a custom frame FreeWheel. Footrests on risers will require some modification to footrest (modified D install under front leading edge), and pinning of bottom clamp.

FreeWheel Folding Wheelchair FreeWheel Compatibility Chart*

The following folding wheelchairs will work with the FreeWheel Wheelchair Attachment WITH the ADAPTOR for Folding Wheelchairs. The ADAPTOR requires 2.5” of unobstructed side, vertical tubing to attach the ADAPTOR with single footrests. Folding chairs rarely have back tubing for FreeWheel storage on back of wheelchair using the included Perch Post.

Manufacturer	Wheelchair Model	Notes
TiLite	2GX2, 2Gx, TX	
Sunrise Medical	Zippie/Youngter 3, Quickie Life, Life I, Easy 200, 300 and 160i.	Side vertical tubing under 13”, center-to-center, may require modification of Adaptor.
	Quickie 2, Quickie QXi	
Invacare	Tracer SX5, Tracer, IV, Tracer, EX2, 9000 SL, 9000 Jymni Pediatric	Side vertical tubing on 9000 Jymni Pediatric wheelchair under 13”, center-to-center, may require modification of Adaptor.
Colours	ActiFold	

This is not a complete list of all wheelchair models that the FreeWheel and FreeWheel Folding Chair Adaptor will work with. Please contact us if you have any questions regarding FreeWheel’s ability to fit a certain type of chair. We can be reached at info@gofreewheel.com— Thank you!