



Molift Duo

A free standing and flexible gantry system





Molift Duo

A flexible, free standing gantry system with no need of installation in wall, ceiling or floor. The perfect solution for most temporary hoist situations.

The easy solution for all bed and room sizes





Molift Duo offers a system that is easy to install without being fixed to the ceiling or walls. It is ideally suited to temporary hoist applications in hospitals, nursing homes or private homes. In combination with Molift Air or Molift Nomad hoist and a sling, Molift Duo solves most overhead lifting situations, whatever the environment. With a lifting capacity of 300 kg, it is ideal for lifting and transferring a person from or to; bed, floor, chair, wheelchair, toilet etc.

Lightweight construction

Molift Duo has a lightweight aluminium construction and a functional design which makes the gantry system easy to install, dismantle and move elsewhere.

This 2-post gantry system has a rail length between 2000 mm and 3500 mm and has 4 step adjustable heights up to 2550 mm, making this free standing hoist unit suitable for all bed and room sizes.



Product name	Item no.	Weight	
Molift Duo 2000	1600013	Maximum user weight (SWL)	300 kg
Molift Duo 2500	1600015	Total weight	38 - 44,3 kg
Molift Duo 3000	1600017	Material	
Molift Duo 3500	1600019	Anodized aluminium, Plastic	
Weight of unit		 	
Post (complete)	14,8 kg		
Post (without foot)	11 kg		
Boom	4,2 kg/m		
Foot	3,8 kg		
			
		Clean / disinfect the product using a standard non-abrasive cleaning agent with a pH between 5 and 9, or with a 70% disinfectant solution. The product can be safely cleaned in a cabinet washer-disinfector at 85°C for 3 minutes.	

Features



Easy to install

Due to the lightweight system, it is easy to install the gantry and move it to a new location.



Supreme stability

The gantry system has supreme stability and a functional design.




Adjustable height

Height positions can be set at four fixed positions from 225 cm up to 255 cm. The height is adjustable with the quick release pin.



Charging of lift motor

Both Molift Air and Nomad motors can be charged with the hand control which is placed in the charger mounted on the post.

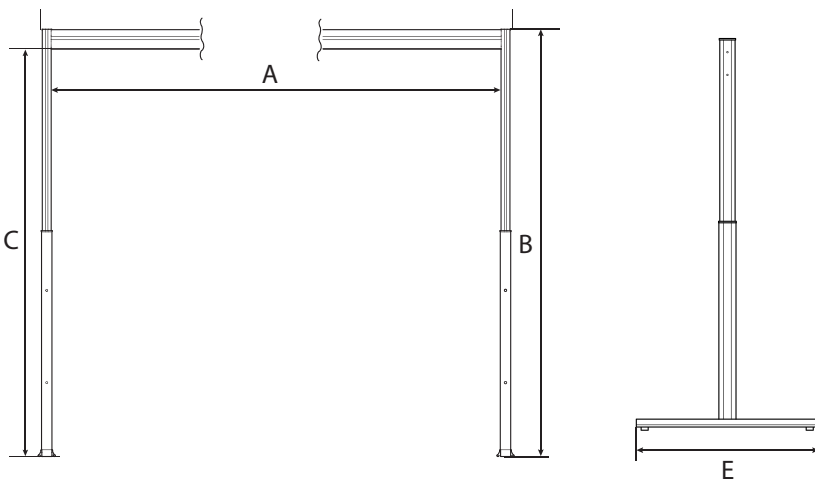

Measurement

A: Rail lengths	2000, 2500, 3000, 3500 mm
B: Height	2250 - 2550 mm
C: Internal height	2400 mm
E: Depth	1200 mm
Max internal space	3500 x 2400 mm





“Everyone should be able to live a free and independent life
and pursue their dreams regardless of any physical circumstances”

For the latest news and continuously updated
product information – visit: www.etac.com

Etac AB
Kista Science Tower
SE-164 51 Kista, Sweden
Tel +46 371 - 58 73 30
Fax +46 371 - 58 73 90
info@etac.se www.etac.com


Creating Possibilities