

Molift Quattro

The temporary hoist solution with excellent stability

A freestanding traverse gantry

Molift Quattro is a freestanding hoist unit that offers a (four post) Traverse System which is the best solution for temporary hoisting needs or when the facility is not built for fixed installations. The gantry is easy to install without the need to mount additional fixtures on walls or ceilings. It is ideally suited for temporary hoist applications in nursing or private homes.




Solves most overhead lifting situations

The traverse has a lightweight aluminum frame and four steel posts that provides excellent stability and covers the surface of the hoisting requirement. Molift Quattro solves all overhead lifting situations together with the Molift Air or Molift Nomad hoist motor and a sling. With a lifting capacity of max 300 kg, it is ideal for lifting and transferring a person from or to bed, floor, chair, wheelchair, toilet etc.

Molift Quattro can be delivered in two sizes: 3x3 meters and 4x4 meters. The length of the posts are fixed. The lowest allowed distance from floor to ceiling is 2.40 meters.

Molift Air or Molift Nomad with trolleys are ordered separately.

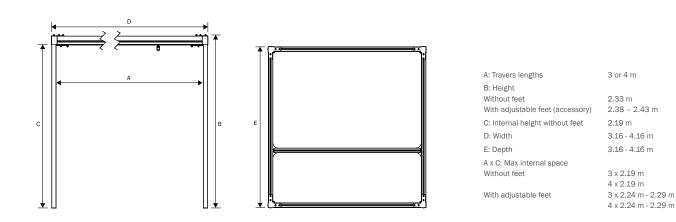

Product name	Internal space	Item no.	Facts		
Molift Quattro 3	3x3 m	1600030	Safe Working Load: 300 kg	max R max	
Molift Quattro 4	4x4 m	1600035	Material: Steel, Anodized Aluminum, Plastic Superior Expected Lifetime: 10 years $300 \atop lbs$		
Facts Weight of unit:			Cleaning		
Skinner Kg/m: 4.2 kg/m/8.8 lbs/m Leg with adjustable feet: 8.5 kg/18.7 lbs Corner: 2.7kg/6 lbs			Clean/disinfect the product using a standard non-abrasive cleaning agent with a pH between 5 and 9, or with a 70% disinfectant solution.	Alk. 70% 5-9	ϵ
Total weight (Complete u Total weight (Complete u	, , , ,				

Molift Quattro

Molift | a part of Etac


Features



Supreme stability and manoeuvrability

persons.

Molift Quattro is a gantry system which has supreme stability and provides access in all lifting situations.

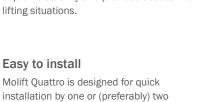



Charging of hoist motor

Molift hoist motors can easily be charged through the hand control or through the in-rail charging option.



Accessory

Adjustable Foot

The Molift Quattro can be optimized with a height adjustable foot on each leg.. Item no 1610200



For the latest news and continuously updated product information – visit www.etac.com

Etac AB
Kista Science Tower
SE-164 51 Kista, Sweden
Tel +46 371 - 58 73 30
Fax +46 371 - 58 73 90
info@etac.se www.etac.com

