

DEVICE PRODUCT LIST

HISTOLAB making your beauty more beautiful

OMEGA LIGHT

Multipurpose near-infrared LED equipment
resolving complicated skin problems!

What is LLLT?

(Low Level Laser Therapy)

UV light harmful for the skins or unnecessary IR are not released, even though it uses a powerful high brightness special LED. It has some advantages such as fewer side effects and fewer damages of the tissues or eyes due to low energy.

Principle

OMEGA LIGHT facilitates regeneration, soothing, easing inflammation and relief of pains by irradiating lights with various wavelengths ranging from 423nm to 640nm. Since it is available Blue Light and Red Light with cross mode every half second, it is optimized LLLT near-infrared LED equipment for Acne PDT treatment (Photo Dynamic Therapy).

OMEGA LIGHT Treatment Effect

* The latest special LED application (287 LEDs)



TYPE OF LIGHT SOURCE



[Red]

Light Dosage (wave) 640nm

Effect of Light

- Atopic, Tinea, Psoriasis
- Activation of the skin

Wave Length 1-6mm




[Green]

Light Dosage (wave) 532nm

Effect of Light

- Improvements of sensitive skins
- Stabilization of body and mind
- Soothing effects

Wave Length 0.5-2mm

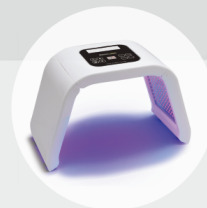

[Blue]

Light Dosage (wave) 423nm

Effect of Light

- Soothing skin troubles
- Skin balance control
- Suppression of sebaceous glands

Wave Length 1mm



[Yellow]

Light Dosage (wave) 583nm

Effect of Light

- Flush / Erythema
- Relief of redness

Wave Length 1-2mm

CRYO LAB

Skin Cooling System, Nutrient Infiltration, and Facelift

Multifunctional Post Laser Aesthetic Device



Skin Cooling and Healing System

Rapidly soothe and heal compromised skin with wide temperature control (-10 °C to 45 °C)



Microcurrent Anti-aging Treatment

This safe and painless facial erases fine lines and wrinkles, lifts muscles and reduces puffiness. The overall results are a younger and healthier skin.



Nutrient Infiltration

Iontophoresis treatment delivers nutrients (such as vitamins, anti-oxidant, peptides) to deeper skin layers



Cryo and Ionto Hand Piece



Micro Current Hand-Piece



Iontophoresis Ground Patch

CRYO LAB BENEFITS

CRYO + NUTRIENT INFILTRATION



Soothe and calm irritated skin



Vitalize and brighten aged dull skin



Hydration and rejuvenation



Pain relief



De-puffiness and facelift

MICROCURRENT

VITAL LAB

Medical Skin Care for various
regeneration programs!



Peeling without irritation
Skin care for even sensitive skin with fine ultrasonic vibration



Various skin care programs
Provides a range of customized skin care with this 1 device



Excellent regeneration effect
Excellent whitening and regeneration effect with stable electrical current!

Principle

It is designed a multi-medical skin care device which provides scrub care for dead skin cell and ultrasonic care for soothing effect and iontophoresis care for infiltration of ionized nutrition into the skin with electrical current

VITALION III EFFECT



Regeneration care post laser treatment



Dead skin cell and waste removal



Moisturizing for dry skin



FUNCTION



VITALION III, SONOZYME common use

SKIN SCRUBBER

Helps lifting and scaling function and infiltration of nutrition with 28kHz powerful vibration.



VITALION III, SONOZYME common use

ULTRASONIC WAVE

Fine ultrasonic vibration helps skin elasticity, skin regeneration and body shape management with infiltration of nutrition

SONO LAB

Premium Medical Skin Care with a convenient and effective lonzyme system!



Non-irritation
Safe treatment with stable electrical current



Ease of use
Uniform infiltration of nutrition for whole face at once



Excellent effect for treatment time
Excellent skin care effect with upgraded Sonozyme function

Principle

It is designed a medical skin care device with the best optimized system including infiltration of nutrition function and chemical peeling using lonzyme principle and soothing care for damaged skin and dead skin cell removal with fine ultrasonic vibration

SONOZYME EFFECT



Soothing, regeneration care, after squeezing acne



Elasticizing care for aged skin



Whitening care for dark and somber skin



VITALION III, SONOZYME common use

IONTOPHORESIS

Promotes skin regeneration following laser treatment and vitalizing on the skin for infiltration of ionized water-soluble nutrition such as IONTO C, V, P products



only SONOZYME use

CROCODILE CLIP

Helps infiltration of nutrition for whole face that is not a part of face with stable electrical current.

AQUA LAB

Smart pore cleaner sucking up every waste in pores without irritation!



Effective sebum removal in pores
Various waste removal effect such as sebum in pores and Acarus folliculorum that are difficult to remove



Simple treatment without side effect
Ideal peeling treatment for busy contemporary people that doesn't affect their daily lives




Peeling and Moisturizing at once
Helps peeling and moisturizing at the same time with one treatment


Principle

AQUA LAB is a smart brightening peeling device which combines all together for chemical peeling, mechanical peeling and physical peeling. Chemical peeling helps effectively dissolves dead skin cells and sebum in pores with AHA and BHA peeling solution with AQUA LAB. Mechanical peeling helps suctioning operation for the wastes in the skin with powerful tornado vacuum principle. Physical peeling helps remove excessive dead skin cell with patented particular tip of AQUA LAB

In Vacuum attachment (AQUA LAB TIP)?

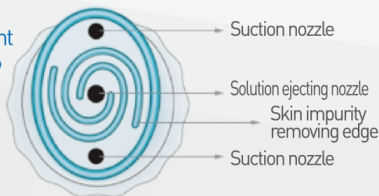

AQUA LAB EFFECT

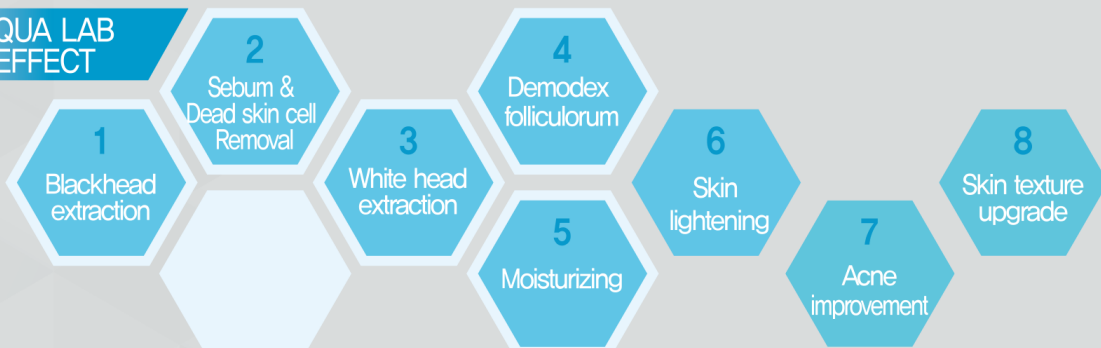

RF-LAB

Radio Frequency + Vacuum Suction + LED Color Therapy

Triple-action mechanism
for maximized FACIAL & BODY contouring



Body slimming effect
Excellent lipolysis and slimming effect with multi function of high-frequency current and vacuum suction



Bi-Polar Radiofrequency
Safe and convenient treatment without Plate!



Penetration of nutrition
Maximizes penetration of nutrition and skin regeneration with cell activation by color therapy from LED

Principle

Radio frequency is non-surgical, non-invasive treatment which helps the skin look younger by safely heating up the deep collagen rich layers of skin, which tighten existing collagen and stimulates the formation of new collagen. It also liquefies and shrinks fat cells to achieve body sculpting and smoother contours. At the same time, vacuum suction improves lymphatic circulation and remove waste products, and LED therapy helps skin regeneration and nutrient penetrations. These functions work together to achieve the best outcome of a more lifted, firmer and younger skin.



BODY



FACE & BODY



FACE



RF-LAB EFFECT

EXCELLENT;
LIFTING CARE!
OBESITY CARE!
RELIEF OF PAINS!



- Face : V line Care, Lifting, the synthesis of Collagen
- Toxin elimination, Neck wrinkle care, Swelling care post plastic surgery
- Neck • Shoulder • Waist muscle pain care, Pain reduction
- Extra flab removal such as love handles, Diet effect
- Self-healing power increase with body temperature increase
- Discharge waste from body, lymphokinesis swelling, and whole body slimming care

Lo
m ips



HISTOLAB
HISTOPLANT CELL LABORATORY

CATLAS INTERNATIONAL TRADING INC.

TEL | 1-647-949-8015 EMAIL | info@histolab.ca

WEBSITE | www.histolab.ca