

Extensions - Toe Length

Information Provided By Kevin Orthopedic Institute



Toe Extension

Toe Extension Myolite 4.5 3 1.5 (mm)

Myolite material frame-extension to toes

FUNCTIONS:

- Forefoot cushion
- Metatarsal punch out material

CLINICAL INDICATIONS:

- Leg length discrepancy
- Sensitive feet
- Accommodations requiring material

A toe extension is a Myolite forefoot padding that extends from the superior dorsal edge of the orthotic frame to approximately 12mm past the toe of the patient's foot. The extension is oversized along with corresponding top and bottom covers, so the orthotic may be trimmed down for appropriate size of patient shoe gear.

Note: If a met head punch is requested, a forefoot extension is required.

CLINICAL PEARL


Extensions add bulk and thickness to an orthotic. Consider topcover thickness and amount of space available in shoe gear when selecting extension thicknesses of 3mm or greater.



Superior view




4.5mm

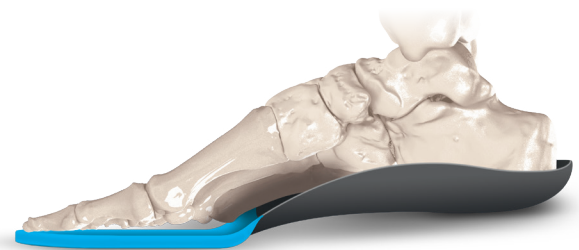

4.5mm




3mm

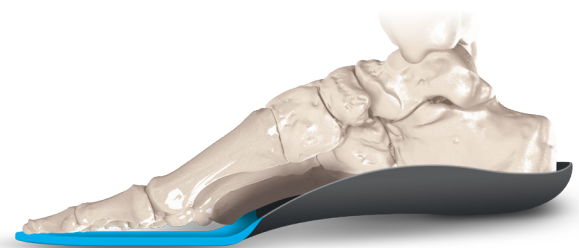

3mm



1.5mm

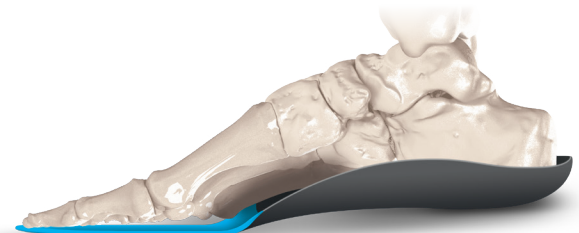

1.5mm

Anterior view

Medial view

Notes: All illustrations and diagrams are of right foot
Colors on illustrations are for visual purposes and will vary on final product