

Extensions - Sulcus Length

Information Provided By Kevin Orthopedic Institute



Sulcus Extension

Sulcus Extension Myolite 4.5 3 1.5 (mm)

Myolite material frame-extension to sulcus

FUNCTIONS:

- Metatarsal head cushion
- Metatarsal punch out material

CLINICAL INDICATIONS:

- Low volume shoes
- Accommodations requiring material
- Forefoot fat pad atrophy
- Hammertoes
- Mallet toes
- Claw toes
- Leg length discrepancy

A sulcus extension is a Myolite forefoot padding that extends from the superior distal edge of the orthotic frame and terminates just distal to the metatarsal heads.


Note: If a met head punch is requested, a forefoot extension is required.



Superior view




4.5mm

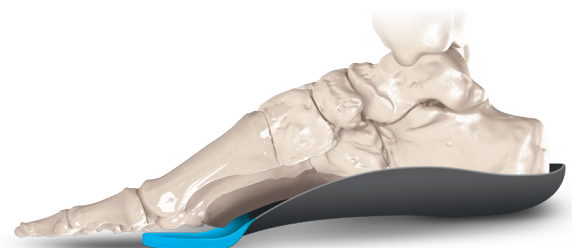

4.5mm




3mm

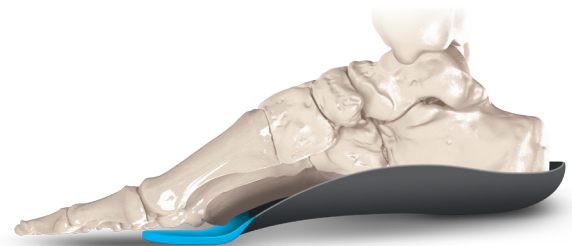

3mm



1.5mm

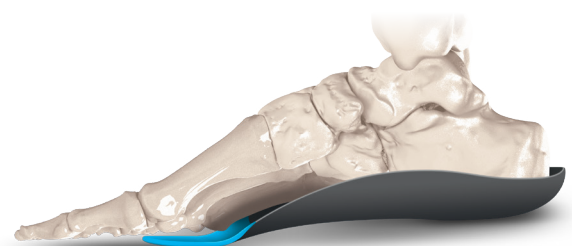

1.5mm

Anterior view
(phalanges excluded)

Medial view

Notes: All illustrations and diagrams are of right foot
Colors on illustrations are for visual purposes and will vary on final product