



**ORMESA®**  
• MADE IN ITALY SINCE 1980 •

**BUG**

BLU

# BUG.

## Safe, modular & versatile.

Bug by Ormesa is a stroller with a tilt-in-space function incorporated in the seat.

This aid is highly engineered and versatile; the wide range of adjustments, makes Bug a **long-lived pushchair**, able to accompany the child during the entire growth process.

Bug is available in two sizes: **Small and Medium.**



BUG with base 869 passed CRASH TEST according to the **ISO 7176/19** and **ANSI RESNA WC/19**.

PRODUCT CERTIFIED BY  
**CERAH**  
FRANCE

PRODUCT TESTED BY  
**BERLIN CERT GmbH**  
AND LISTED IN THE GERMAN HMV (Technical Aids Listing)

PRODUCT REGISTERED **FDA** USA

PRODOTTO MARKED **CE**



# Modular Seating System

## Quick and simple to adjust for optimum posture

The seating adjustment hinges are close to the user's joints (coxo-femoral, knee & ankle joints). This structural feature is a major contributing factor to the child's well-being, as **the posture is not changed during tilting and reclining operations.**

## Tilt in space and reclining

Tilting and reclining of the backrest make it possible to get the **ideal support position**. They are regulated by two independent gas springs that makes adjustment operations particularly gentle & smooth. **At full recline, the chair assumes a nearly supine position.**



# Reversible Seat

## face to face with the child

The seat of the Bug measure Small with base 869 is **reversible**. Can be used in **two directions: forward facing or facing you position** (mother direction).

The forward direction will favor the visual **exploration of the surrounding environment**, the mother direction instead the **visual contact between the child and the adult**.

forward facing



# Adjustable on 4 bases

## outdoor use



### 869 four-wheel outdoor base

Light solid wheels, sliding on all surfaces, with direction locks, independent drum brakes.

Front wheels: 17 cm.

Rear wheels: 25 cm (Small) and 30cm (Medium).



BUG with base 869 passed CRASH TEST according to the **ISO 7176/19** and **ANSI RESNA WC/19**.



# and on any flat surfaces

## indoor use

### 856 indoor high-low feeding base

With 10 cm solid pivoting wheels and brakes. Simply and safely height adjustable, with pedal gas spring. Height goes from 34 to 61 cm.



### 851 indoor high-low compact base

With cm.7.5 solid pivoting wheels and brakes. Simply and safely Height adjustable, with gas spring. Height goes from 23 to 39 cm.



### 850 indoor compact floor base

It allows the development of activities on the ground, maintaining maximum stability and safety.










# Technical features and adjustments



**ATTENZIONE!**  
ACCERTARSI CHE LA LEVA DI BLOCCAGGIO SIA INSERITA, PRIMA DI USARE IL PASSEGGINO

**ATTENZIONE**  
LEGGERE E COMPRENDERE LE ISTRUZIONI PER L'USO E LA MANUTENZIONE  
ATTENTION: READ AND UNDERSTAND THE USE AND MAINTENANCE INSTRUCTIONS  
BEFORE USING THE EXERCISE EQUIPMENT  
ATTENTION: READ AND UNDERSTAND THE INSTRUCTIONS FOR USE AND MAINTENANCE BEFORE USING THE EXERCISE EQUIPMENT

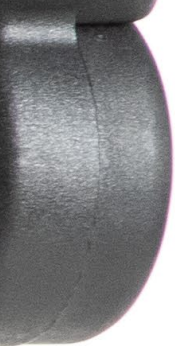

## Headrest

incorporated in the backrest. Your child's head is always supported during the tilt adjustments of the seat or of the backrest, also without accessories.



## Locking/unlocking system

on a flat surface that ensures a stable positioning well visible and extremely safe. The release lever has a safety button to avoid the accidental release of the seat from the frame.



## Push handle

adjustable in tilt through red buttons. It helps the care-giver in handling it even when Bug is tilting and/or reclining. (Only for 869 four-wheel outdoor base Bug).



## Backrest

concave, rigid, reclinable and adjustable in height.



## Tilting seat

with gas spring and adjustable in depth, width and rotation. The seat of the Bug measure Small with base 869 is reversible.



## Leg rest

adjustable in tilt, with hinges close to the knee joints.

## Tip-up foot rest

adjustable in height and tilt for the foot flex-extension.



# Padded, ergonomic, re fixed with buttons and made by “made in It

dark grey  
pink



dark grey  
pink & light blue



dark grey  
green



dark grey  
green & light blue



light grey  
pink



# movable upholstery

## aly" high quality fabric

light grey  
pink & light blue



light grey  
green



light grey  
green & light blue



- › **Fire-proof:** its fire-proof properties are inherent to the yarn and they do not come from any additional treatments.
- › **Recyclable:** production cycle with low environmental impact.
- › **Martindale method:** maximum resistance to rubbing: 100.000 cycles. UNI ENI ISO 12947.
- › **Oeko-Tex Certificate:** for all fabrics that are in direct contact with the skin, to guarantee quality and safety for health.
- › **Long-lasting and highly resistant colours:** they are resistant to sweat, light, dry cleaning and washing machine (see its Use and Instruction manual).



CERTIFICATE  
MED annex  
A.1/3.19-20-21

# Components



**863 Headrest**  
with parietal supports. Adjustable in height, width, tilt and backward forward.



**852 Headrest**  
with parietal supports. Adjustable in height and width.



**942 Ergonomic headrest**  
Adjustable in height, tilt and forward-backward with adjustable sides. Size L (large) available for special case.



**835 Shaped headrest**



**934 Suspension headrest attachment**  
Available only for 863 headrest.



**838 Lateral supports**  
adjustable in height and width.



# Components



**868 Wrappable and flexible trunk supports**  
Adjustable in height and width.



**834 Padded abduction block**  
Adjustable in depth.



**834R Padded abduction block**  
Adjustable in depth and abduction.



**853 Vest harness**



**903 Five point vest harness**



**906 Five point harness**

# Components



**894 45° pelvic belt**



**920 Four Point belt**



**827 Foot straps**



**839 Front handle**  
Adjustable in tilt.



**824 Transparent tray**  
wrappable, with recess. Adjustable in tilt.



**824B Wide tray for indoor base**  
It is made of resistant and hygiene material, no edge design, it does not retain any residue of food.

# Components



**905 Hand brake**  
for 869 base.



**935 Oxygen bottle holder for 869 base**  
Maximum dimensions of the oxygen bottle: height 600 mm (valve included) and diameter 140 mm.  
**The frame can be folded without removing the accessory.**



**911 Vent tray**  
for carrying medical equipment. For 869 base. **The frame can be folded without removing the accessory.**



**819 Canopy**  
in combination with the upholstery. Resistant and easily adjustable.



**825 Rain cover**  
(to fix at the canopy 819) for 869 base.



**818 Thermal cover**

# Components



**892 Covering for footrest**  
padded and removable.



**854 Parasol**  
For 869 base, in combination with the upholstery.



**858 Net basket for 869 base**



**921 Sibling rear footboard**  
Placed behind the stroller. For children of at least 15 months. Max. load 20 kg (44 Lbs.) For 869 base.



**922 Bottle holder**  
For 869 base.




**891 Set of 4 tie down hooks**  
for the forward-facing transport of the pushchair in a motor vehicle (private/cars, buses...). Only for 869 base.  
*See User and Maintenance Manual provided with the aid.*



# Bug | 869 base sizes and weight

Modular postural system on four-wheel outdoor base.

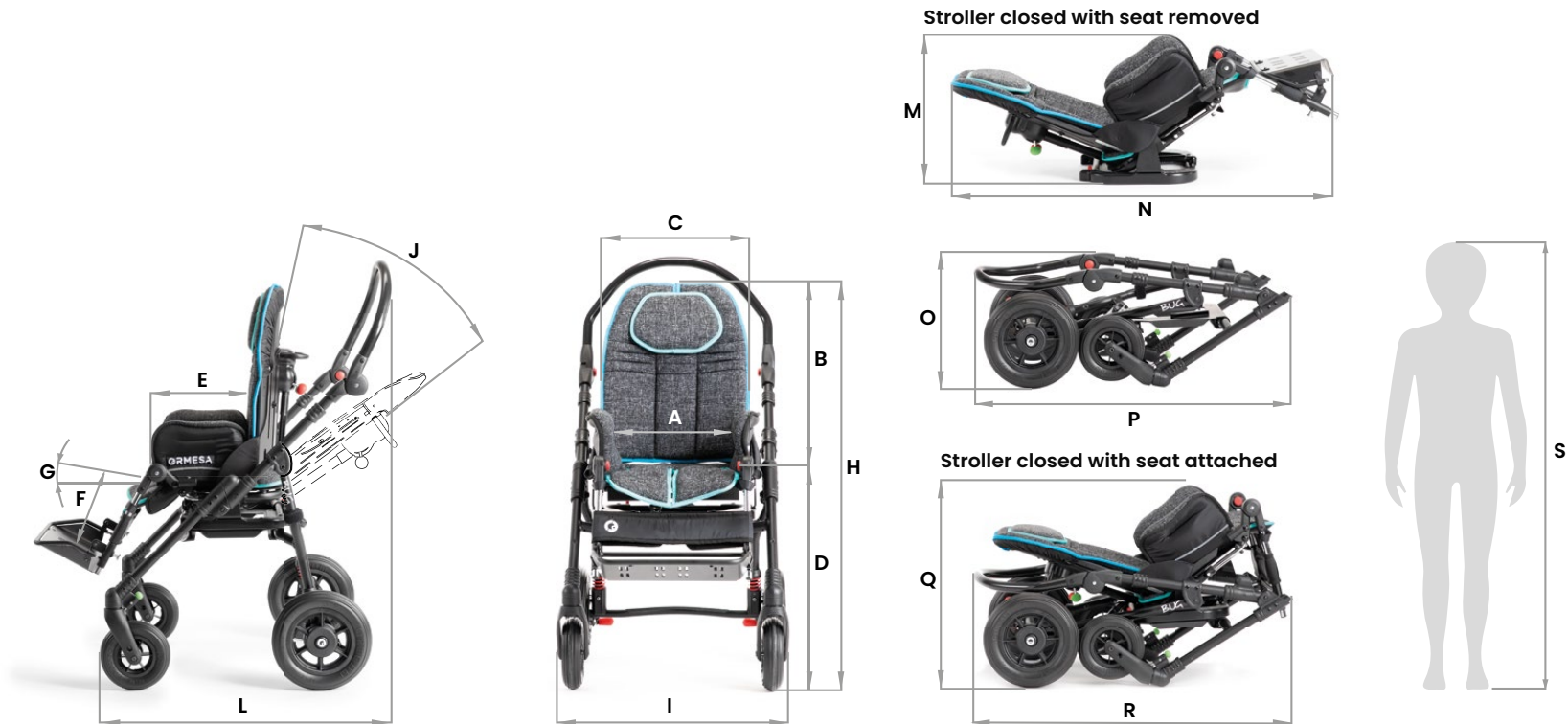

	A (cm)	B (cm)	C (cm)	D (cm)	E (cm)	F (cm)	G (degrees)	J (degrees)	H (cm)	I (cm)	L (cm)	M (cm)	N (cm)	O (cm)	P (cm)	Q (cm)	R (cm)	S (cm)	max. load	frame weight	seat weight	total weight
<b>Small</b>	20-30	52-62	37	52	20-30	14-30	0°-35°	35°	100-110	57	67-96	35	88	33	82	50	81	90-120	35 kg	9,5 kg	8,5 kg	18 kg
<b>Medium</b>	25-36	60-75	46	55	30-40	20-38	5°-30°	40°	112-127	62	80-107	37	112	39	94	52	94	115-140	45 kg	11 kg	12,7 kg	23,7 kg

Hypen-separated numbers indicate the minimum and maximum range of adjustment

# Bug | 856 base sizes and weight

Modular postural system on indoor high-low feeding base.




	A (cm)	B (cm)	C (cm)	D (cm)	E (cm)	F (cm)	G (degrees)	J (degrees)	H (cm)	I (cm)	L (cm)	M (cm)	N (cm)	O (cm)	P (cm)	max. load	frame weight	seat weight	total weight
<b>Small</b>	20-30	52-62	37	47-75	20-30	14-30	0°-35°	35°	95-132	56	34-61	80	40	25	90-120	35 kg	14,5 kg	8,5 kg	23 kg
<b>Medium</b>	25-36	60-75	46	47-75	30-40	20-38	0°-25°	40°	106-148	62	34-61	83	45	32	115-140	45 kg	16,2 kg	12,7 kg	28,9 kg

Hyphen-separated numbers indicate the minimum and maximum range of adjustment

# Bug | 851 base sizes and weight

Modular postural system on indoor high-low compact base.

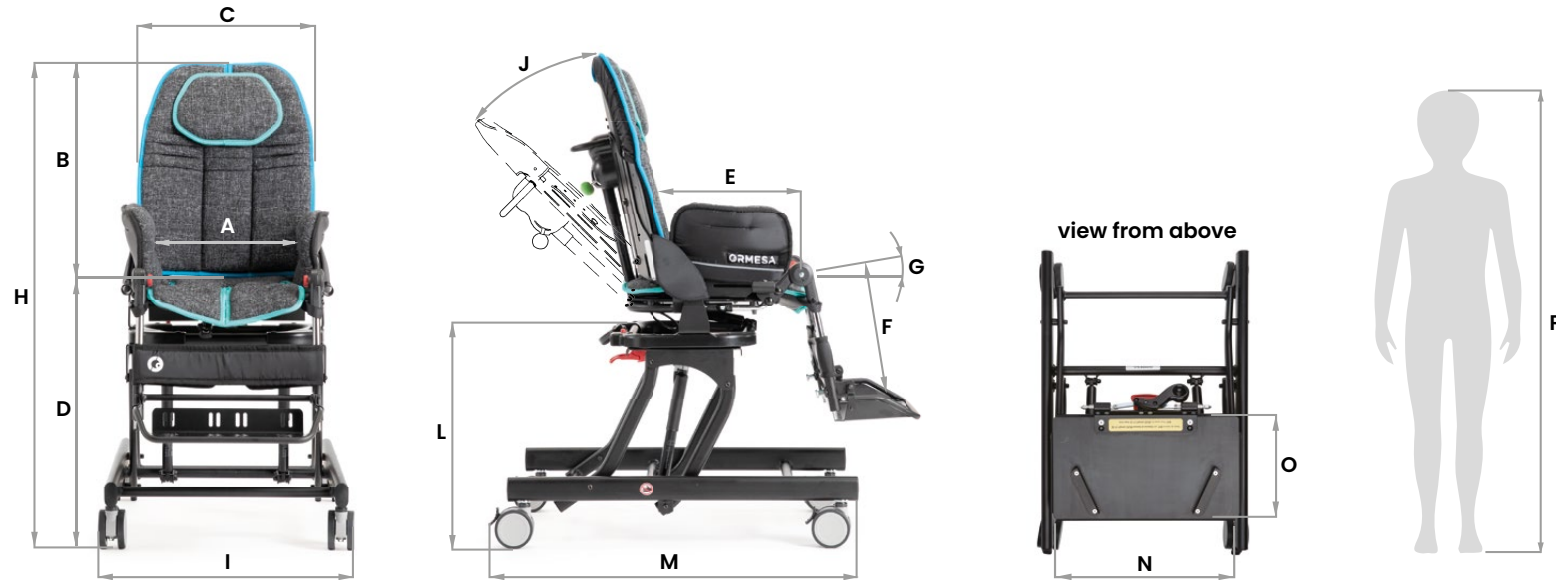


	A (cm)	B (cm)	C (cm)	D (cm)	E (cm)	F (cm)	G (degrees)	J (degrees)	H (cm)	I (cm)	L (cm)	M (cm)	N (cm)	O (cm)	P (cm)	max. load	frame weight	seat weight	total weight
<b>Small</b>	20-30	52-62	37	34-50	20-30	14-30	0°-35°	35°	84-110	50	23-39	72	41	23	90-120	35 kg	8,8 kg	8,5 kg	17,3 kg
<b>Medium</b>	25-36	60-75	46	34-50	30-40	20-38	0°-25°	40°	95-126	55	23-39	76	46	28	115-140	45 kg	10 kg	12,7 kg	22,7 kg

Hyphen-separated numbers indicate the minimum and maximum range of adjustment



# Bug | 850 base sizes and weight

Modular postural system on indoor compact floor base.

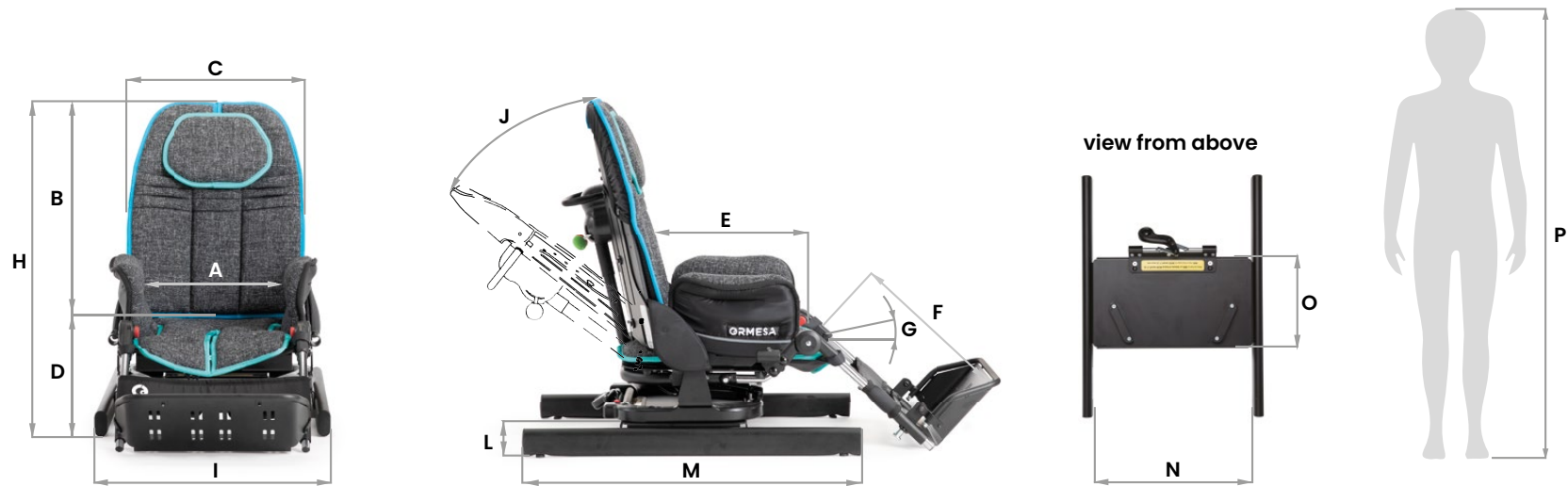

	A (cm)	B (cm)	C (cm)	D (cm)	E (cm)	F (cm)	G (degrees)	J (degrees)	H (cm)	I (cm)	L (cm)	M (cm)	N (cm)	O (cm)	P (cm)	max. load	frame weight	seat weight	total weight
<b>Small</b>	20-30	52-62	37	17	20-30	14-30	0°-35°	35°	68-78	46	6	61	41	23	90-120	35 kg	4,5 kg	8,5 kg	13 kg
<b>Medium</b>	25-36	60-75	46	17	30-40	20-38	0°-25°	40°	77-92	51	6	63	46	28	115-140	45 kg	5,5 kg	12,7 kg	18,2 kg


Hyphen-separated numbers indicate the minimum and maximum range of adjustment





ORMESA

**ORMESA<sup>®</sup>**  
• MADE IN ITALY SINCE 1980 •

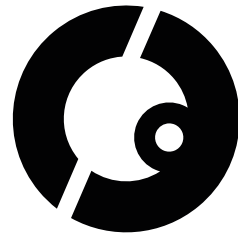

COMPANY WITH  
QUALITY SYSTEM  
CERTIFIED BY DNV GL  
= ISO 13485 =

---

**Ormesa srl** - Via A. Da Sangallo, 1 | 06034 Foligno (PG) - ITALY  
Tel. +39 0742 22927 | FAX +39 0742 22637 | [info@ormesa.com](mailto:info@ormesa.com)

**[www.ormesa.com](http://www.ormesa.com)**