

Bovie®

OPERATING WITH PRECISION



QUALITY SURGICAL PRODUCTS, WITH PROCEDURE-FOCUSED DESIGN

Since its establishment in the late 1970s, Bovie Medical has been striving to continue with the pioneering work of Dr. Bovie, credited as the inventor of modern electrosurgery, by continuing to develop and manufacture advanced electrosurgical products. Manufactured in the USA, Bovie Medical is recognised for their quality design and manufacturing as well as excellent customer service. It has become a leading manufacturer that sets the standard for electrosurgical generators and accessories worldwide.

BOVIE HIGH FREQUENCY DESICCATORS - IT'S YOUR CHOICE

Advantages of choosing Bovie;

Three, low-end electrosurgical devices available
Choice of modes for varying applications
Low-cost models for everyday procedures
Standard accessories across the range
Huge choice of dermal tips for a wide range of procedures
Dermal tips are compatible with (and often cheaper than),
the Conmed Hyfrecator® product
Made in the USA with a 4-year warranty






ELECTROSURGERY - BETTER BY BOVIE!



Superior products, quality and service.
www.lsrhealthcare.com.au

MEDICAL EQUIPMENT SPECIALISTS™

BOVIE'S SUITE OF HIGH FREQUENCY DESICCATORS FOR OFFICE USE

	DERM 102	DERM 942	BANTAM PRO A952
			
Overview	The Derm 102 for Family Practice is an economical, high frequency desiccator, designed to give practitioners the ability to perform minor skin procedures in-office.	Derm 942 - designed and engineered to be the most reliable and durable 40 watt desiccator available today for a wider range of applications.	The Bovie Bantam PRO, A952 is a High Frequency Desiccator with all the capabilities of the Derm 942, plus Cutting ability.
Max. Power	10 watts	40 watts	50 watts (40W Bipolar)
Adjustment	0.1W	0.1W (up to 10W) then 1W	0.1W (up to 10W) then 1W
Monopolar	Yes	Yes	Yes
Bipolar	Yes	Yes	Yes
Warranty	4 years	4 years	4 years

Products and specifications are subject to change without notice. ©2023

Common Office Procedures

DESICCATION	LESION	SUGGESTED POWER RANGE	LESION	SUGGESTED POWER RANGE	HEMOSTASIS		SUGGESTED POWER RANGE
	Angiomas (Cherry)	1.2 - 2.0 watts	Pyogenic Granuloma	5.0 - 10.0 watts		Ellipse excision	5.0 - 10.0 watts
Telangiectasias	1.5 - 2.0 watts	C and D	5.0 - 10.0 watts	Punch excision	5.0 - 10.0 watts		
Seborrheic Keratosis	3.0 - 4.0 watts	Skin Tags (Acrochordons)	1.0 - 3.0 watts	Shave excision	5.0 - 10.0 watts		
Warts	5.0 - 15.0 watts						

Accessories

Bovie Medical provides a complete portfolio of accessory options for your delicate procedural needs. Below are examples of most commonly used accessories. For a more complete list, please visit our website - www.LSRHealthcare.com.au



Sharp Disposable Tips - Sterile/Non Sterile
BV-A805/BV-A804



Blunt Disposable Dermal Tips - Sterile/Non Sterile
BV-A806/BV-A807



Bipolar Forceps
BV-381-250A/BV-381-250B/
BV-381-250C



Reusable Angled Fine Tip Tungsten Epilation Needle
BV-A34T



Disposable Handpiece Sheaths
BV-A910




3-button Power Control Pencil and Lead
BV-A902



Bipolar Cable for DERM 102
BV-A827V



Footswitch
BV-A803

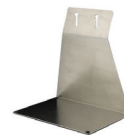

Desk Bracket
BV-A813



Mobile Stand and Bracket
BV-1580ARN



Superior products, quality and service.
[WWW.LSRHEALTHCARE.COM.AU](http://www.LSRHEALTHCARE.COM.AU)

LSR Healthcare Pty Ltd
1300 859 767
www.lsrhealthcare.com.au
enquiries@lsrhealthcare.com.au

