

# Bovie®

## Portable Cautery Unit



**The Bovie portable cautery unit is a lightweight, cordless, thermocautery device used for localized soft tissue trimming and hemostasis.**


Battery-powered. Safe & Simple to use. Convenient & Compact. No wires, foot pedal or grounding plates required. Affordable.

- Provides pinpoint hemostasis
- Removes unhealthy soft tissue and lesions
- Trims gingiva before impressions
- Safer to use than AC-powered cauterizing systems
- DC battery power eliminates the risk of electric shock to patients
- **Can be safely used on patients with pacemakers, defibrators and metal implants as device uses heat not electricity.**

The Bovie unit quickly reaches 2,200°F (1,200°C) for fast and convenient cautery. The unit instantly cools once the activation button is released. Penlight-sized handle weighs only 3 oz. Kit includes cautery unit, four single use disposable tips, and two AA batteries.

### Clinical Uses;

- To control hemostasis prior to final impressions
- For removal of unhealthy tissue in subgingival restorative cases
- Tissue tag removal (Safe around implant placements)
- Removing flaps over partially erupted wisdom teeth
- Fibroma removal

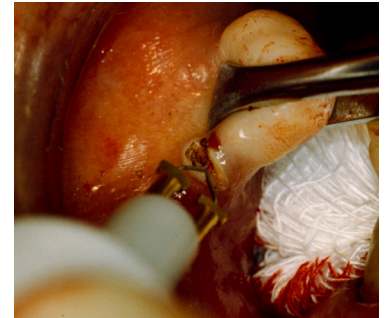

The surgical removal of a fibroma.  
Photography courtesy of Dr. Larry Gaum



H101-Fine tip



H103-Loop tip



DEL1-Portable Cautery Unit

BOVIE PORTABLE CAUTERY UNIT	430102
BOVIE DISPOSABLE FINE TIPS 10-PACK	430101
BOVIE DISPOSABLE LOOP TIPS 10-PACK	430103