


Increased Flexibility, Efficient Cutting, Heat Matrix Technology

Increased safety reduces stress and creates better patient outcomes

Instrument separation is a serious concern in endodontics. Nickel-titanium (NiTi) endodontic instruments were introduced to facilitate instrumentation of curved canals. NiTi instruments are super-elastic and will flex far more than stainless-steel instruments before exceeding their elastic limit reaching their separation point.

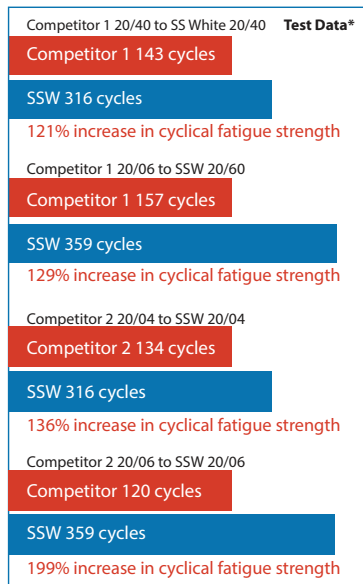


Cyclic fatigue is the leading cause of file separation. Exact Taper H rotary files offer a remarkable leap forward in resistance to cyclic fatigue. The numbers say it all.

This flexibility is an important property that allows preparation of curved canals while minimizing transportation. Despite this increased flexibility, separation is still a concern with NiTi instruments. The clinical concern is that NiTi files have been reported to undergo fatigue failure which occurs unexpectedly without any sign of previous deformation and therefore visual inspection would not seem to be the ideal way of evaluating nickel titanium instruments in order to prevent fracture.

SS White EXACTTaper™ Files offer greater flexibility that is more resistant to cyclic fatigue when tested against both Dentsply Protaper® Gold and Edge Endo EdgeTaper™ Platinum files.



Heat Treatment Provides Better Flexibility to Easily Navigate Tougher Canals

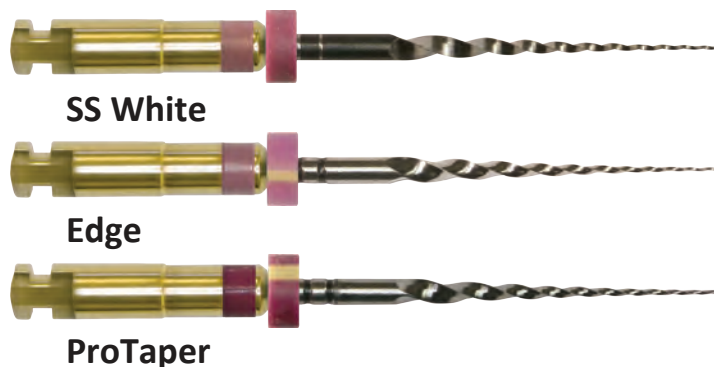

- Same sizes, lengths and tapers enable you to use the same protocol or technique you use with ProTaper® Gold and Edge Taper Platinum™ file systems
- Shaping files cut dentin throughout each canal zone with greater efficiency, while reducing the potential for over-preparation
- Finishing files insure the proper apical fit for the matching EXACTTaper, or your existing Gutta Percha points
- Shaping files have multiple tapers that ensure flexibility and cut dentin in specific canal zones

SKU	DESCRIPTION	PK Size	PROTAPER GOLD	PROTAPER GOLD DESCRIPTION	EDGE SKU	EDGE SEQUEL SAPPHIRE EDGE TAPER PLATINUM
318145	EXACTTaper ^H ™ SX 6 pack 19 mm	6 Pk	PTGRSX19	ProTaper Gold Rotary File Shaper X 19mm	ETSX19HT	SX 19 mm
318121	EXACTTaper ^H ™ SX-F3 ASST 21mm	6 Pk	PTGRA21	ProTaper GOLD Rotary Files Asstd (SX-S2, F1-F3) 21mm	ETP21MM	
318139	EXACTTaper ^H ™ S1 6 pack: 21 mm	6 Pk	PTGRS121	ProTaper Gold Rotary File Shaper 1 Purple 21mm	ETS121HT	S1 21 mm
318142	EXACTTaper ^H ™ S2 6 pack: 21 mm	6 Pk	PTGRS221	ProTaper Gold Rotary File Shaper 2 White 21mm	ETS221HT	S2 21 mm
318124	EXACTTaper ^H ™ F1 6 pack: 21 mm	6 Pk	PTGRF121	ProTaper Gold Rotary File Finisher 1 Yellow 21mm	ETF121HT	F1 21 mm
318127	EXACTTaper ^H ™ F2 6 pack: 21 mm	6 Pk	PTGRF221	ProTaper Gold Rotary File Finisher 2 Red 21mm	ETF221HT	F2 21 mm
318130	EXACTTaper ^H ™ F3 6 pack: 21 mm	6 Pk	PTGRF321	ProTaper Gold Rotary File Finisher 3 Blue 21mm	ETF321HT	F3 21 mm
318133	EXACTTaper ^H ™ F4 6 pack: 21 mm	6 Pk	PTGRF421	ProTaper Gold Rotary File Finisher 4 Double Black 21mm	ETF421HT	F4 21 mm
318136	EXACTTaper ^H ™ F5 6 pack: 21 mm	6 Pk	PTGRF521	ProTaper Gold Rotary File Finisher 5 Double Yellow 21mm	ETF521HT	F5 21 mm
318122	EXACTTaper ^H ™ SX-F3 ASST 25mm	6 Pk	PTGRA25	ProTaper Universal Rotary Files Asstd (SX-S2, F1-F3) 25mm	ETP25MM	
318140	EXACTTaper ^H ™ S1 6 pack: 25 mm	6 Pk	PTGRS125	ProTaper Gold Rotary File Shaper 1 Purple 25mm	ETS125HT	S1 25 mm
318143	EXACTTaper ^H ™ S2 6 pack: 25 mm	6 Pk	PTGRS225	ProTaper Gold Rotary File Shaper 2 White 25mm	ETS225HT	S2 25 mm
318125	EXACTTaper ^H ™ F1 6 pack: 25 mm	6 Pk	PTGRF125	ProTaper Gold Rotary File Finisher 1 Yellow 25mm	ETF125HT	F1 25 mm
318128	EXACTTaper ^H ™ F2 6 pack: 25 mm	6 Pk	PTGRF225	ProTaper Gold Rotary File Finisher 2 Red 25mm	ETF225HT	F2 25 mm
318131	EXACTTaper ^H ™ F3 6 pack: 25 mm	6 Pk	PTGRF325	ProTaper Gold Rotary File Finisher 3 Blue 25mm	ETF325HT	F3 25 mm
318134	EXACTTaper ^H ™ F4 6 pack: 25 mm	6 Pk	PTGRF425	ProTaper Gold Rotary File Finisher 4 Double Black 25mm	ETF425HT	F4 25 mm
318137	EXACTTaper ^H ™ F5 6 pack: 25 mm	6 Pk	PTGRF525	ProTaper Gold Rotary File Finisher 5 Double Yellow 25mm	ETF525HT	F5 25 mm
318123	EXACTTaper ^H ™ SX-F3 ASST 31mm	6 Pk	PTGRA31	ProTaper Universal Rotary Files Asstd (SX-S2, F1-F3) 31mm	ETP31MM	
318141	EXACTTaper ^H ™ S1 6 pack: 31 mm	6 Pk	PTGRS131	ProTaper Gold Rotary File Shaper 1 Purple 31mm	ETS131HT	S1 31 mm
318144	EXACTTaper ^H ™ S2 6 pack: 31 mm	6 Pk	PTGRS231	ProTaper Gold Rotary File Shaper 2 White 31mm	ETS231HT	S2 31 mm
318126	EXACTTaper ^H ™ F1 6 pack: 31 mm	6 Pk	PTGRF131	ProTaper Gold Rotary File Finisher 1 Yellow 31mm	ETF131HT	F1 31 mm
318129	EXACTTaper ^H ™ F2 6 pack: 31 mm	6 Pk	PTGRF231	ProTaper Gold Rotary File Finisher 2 Red 31mm	ETF231HT	F2 31 mm
318132	EXACTTaper ^H ™ F3 6 pack: 31 mm	6 Pk	PTGRF331	ProTaper Gold Rotary File Finisher 3 Blue 31mm	ETF331HT	F3 31 mm
318135	EXACTTaper ^H ™ F4 6 pack: 31 mm	6 Pk	PTGRF431	ProTaper Gold Rotary File Finisher 4 Double Black 31mm	ETF431HT	F4 31 mm
318138	EXACTTaper ^H ™ F5 6 pack: 31 mm	6 Pk	PTGRF531	ProTaper Gold Rotary File Finisher 5 Double Yellow 31mm	ETF531HT	F5 31 mm