


V3 SYSTEM





Get the results you want with the nickel-titanium rings, anatomically shaped matrix bands, and adaptive wedges included in the Triodent V3 Sectional Matrix System.

V3—the perfect system for creating perfect Class IIs

- Superelastic nickel-titanium spring for optimal separation force and tight contacts
- Wide indications for use — glass-fibre tines will not collapse into wide cavities, even when a cusp is missing
- Tines are anatomically shaped for superior matrix retention and adaptation—and minimal finishing
- Tines' V-shape accommodates the wedge
- V3 Tab-Matrix with pronounced marginal ridge and vertical profile for correct anatomical form
- Wave-Wedge is stackable, self-locating, and self-seating for optimal seal on matrix
- Spring tilt allows easy stacking of multiple rings
- Ring is very stable on undercuts and in forceps during placement

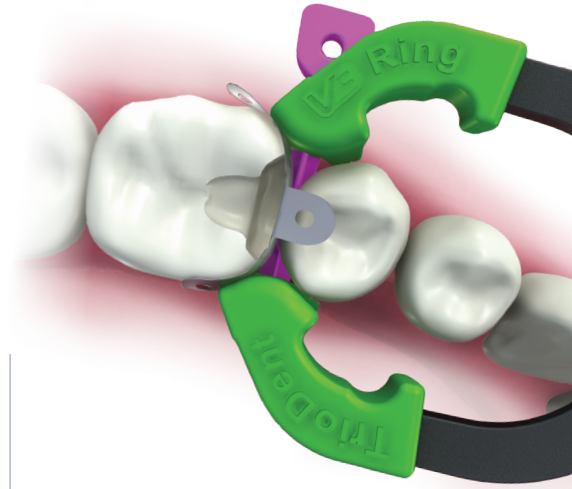

**PERFECT
CONTACTS IDEAL
ANATOMY**

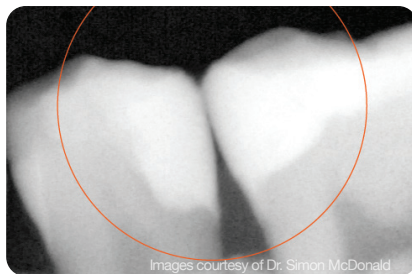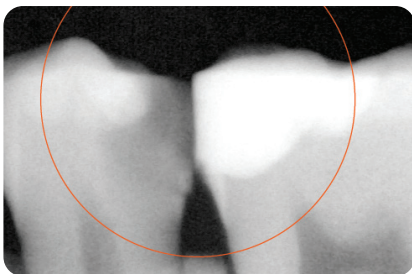

Images courtesy of Dr. Simon McDonald

+ Circumferential Band

- Fails to restore proximal anatomy
- Produces a thin contact at the marginal ridge

The Result: Increased risk of fracture, occlusal interference, caries, and periodontal disease.

+ Triodent V3 Sectional Matrix System

- Creates tight contact points
- Produces natural contours

The Result: Lifelike, durable restorations with anatomically accurate contours and tight contacts.



Starter Kit

V3 SYSTEM COMPONENTS



UNIVERSAL V3 RING - GREEN

- NiTi for ideal spring strength and elasticity
- Long-lasting tines, safe in autoclave
- V-shaped tines accommodate wedge
- Tines shaped for ideal matrix retention and adaptation


NARROW V3 RING - YELLOW

- Tines set closer together for optimal separation of smaller teeth
- Up to 4kg of separating force, even on premolars



SUPERCURVE MATRIX

- Matrix hugs the tooth, avoiding interference with ring placement
- Micro-thin, non-stick finish makes removal easy
- Shaped gingival apron leaves no gaps in the gingival-axial corner
- Rolled ridge on occlusal margin provides beautiful marginal ridges

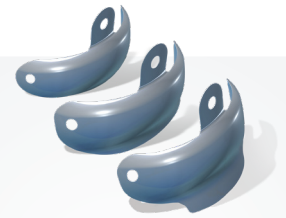

V3 TAB-MATRIX

- Only V3 Tab-Matrix has a pronounced marginal ridge for natural anatomy
- Pronounced curvature on horizontal plan means more curvature remains as matrix wraps around tooth
- Gingival apron is properly shaped to commonly encountered cavity form, avoiding gaps in the gingival-axial corner



WAVE-WEDGE

- Fine wave-shaped wings compress on entry and flare again upon exit for easy placement and ultimate seal
- Hollow underside allows another wedge to be placed from the opposite side of the embrasure



WEDGEGUARD

- WedgeGuard prevents the bur from nicking the adjacent tooth
- After tooth preparation, the guard is detached while the wedge remains in place
- Saves time and protects gingival papillae



PIN-TWEEZERS & FORCEPS

- Pin located in tweezer tip grips holes in Triodont matrices and wedges for absolute control
- Tweezers' naturally closed position means you will never drop a matrix or wedge
- Handy ball burnisher in tweezer tip
- Forceps' locking function and angled-grip arms hold the V3 Ring securely, both mesially and distally

<p>V3 MATRIX KIT</p> <p>1 x V3 Rings Universal (green) 1 x V3 Ring Narrow (yellow) 25 x Wave-Wedges Small (white) 25 x Wave-Wedges Medium (pink) 25 x Wave-Wedges Large (purple) 25 x V3 Tab-Matrices 3.5mm 25 x V3 Tab-Matrices 4.5mm 25 x V3 Tab-Matrices 5.5mm 25 x V3 Tab-Matrices 6.5mm 25 x V3 Tab-Matrices 7.5mm 2 x Pin-Tweezers 1 x Forceps</p>	<p>303301</p>	<p>V3 TAB-MATRIX</p> <p>4.5mm 50pk 4.5mm 100pk 5.5mm 50pk 5.5mm 100pk 6.5mm 50pk 6.5mm 100pk</p>	<p>303372 303375 303373 303376 303374 303377</p>	<p>SUPERCURVE MATRIX</p> <p>Starter Pack 3.5mm (gray) 90pk 4.5mm (green) 90pk 5.5mm (blue) 90pk 6.5mm (purple) 90pk 3.5mm (gray) 50pk 4.5mm (green) 50pk 5.5mm (blue) 50pk 6.5mm (purple) 50pk</p>	<p>303402 303403 303404 303405 303406 303417 303418 303419 303420</p>
<p>V3 RING 2PKS</p> <p>V3 Ring Universal (green) 2pk V3 Ring Narrow (yellow) 2pk V3 Ring Mixed (1 green, 1 yellow) 2pk</p>	<p>303342 303343 303341</p>	<p>V3 SUPERCURVE KIT</p> <p>1 x V3 Rings Universal (green) 1 x V3 Ring Narrow (yellow) 20 x SuperCurve Matrices 3.5mm (gray) 20 x SuperCurve Matrices 4.5mm (green) 20 x SuperCurve Matrices 5.5mm (blue) 20 x SuperCurve Matrices 6.5mm (purple) 20 x SuperCurve Matrices 7.5mm (gray) 25 x Wave-Wedges Small (white) 25 x Wave-Wedges Medium (pink) 25 x Wave-Wedges Large (purple) 2 x Pin-Tweezers 1 x Forceps</p>	<p>303316</p>	<p>WAVE-WEDGE</p> <p>Small (white) 100pk Medium (pink) 100pk Large (purple) 100pk</p>	<p>302066 302065 302064</p>
<p>INSTRUMENTS</p> <p>Pin-Tweezers Forceps</p>	<p>302075 302061</p>			<p>WEDGEGUARD</p> <p>Small (white) 50pk Medium (pink) 50pk Large (purple) 50pk</p>	<p>303383 303384 303385</p>