

# CONTOURS™

COARSE ANATOMY TRIMMERS

## Life-Like Restorations Require Contours



CONTOURS™ from CLINICIAN'S CHOICE® are 90 micron aluminum oxide coarse anatomy trimmers designed to create secondary anatomy on your microfill or hybrid composite restorations. After creating the gross anatomy with your diamond or carbide burs (we recommend EDGE™ Diamonds and Carbides for consistency and efficiency), use CONTOURS for refining the gross anatomy and creating secondary grooves, mamelons, and characterization in the restoration.

Where conventional polishers are inefficient and prone to wear, CONTOURS are unique in that they trim composite quickly, without significant degradation to the trimmer. Moreover, CONTOURS are autoclavable and may be reused several times.


CLINICIAN'S CHOICE has made these long-lasting coarse anatomy trimmers available in the same 6 shapes made popular by the REALITY 5-STAR rated D•FINE™ Polishers for Hybrid Composites, so there's a CONTOURS trimmer perfectly suited for every anterior and posterior restoration you undertake.



CONTOURS Assorted 12-pack  
# 078912



Dentistry and Photography courtesy of Dr. Robert Margoss.

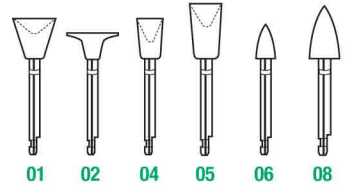

# 078912 CONTOURS 12-pack  
12 trimmers – 2 of each shape

# 0782 CONTOURS 20-pack  
+ 2 digit code of desired shape  
20 trimmers – single shape (specify shape)

CLINICIAN'S  
CHOICE®  
DENTAL PRODUCTS INC.

designed for better dentistry

1-800-265-3444

[www.clinicianschoice.com](http://www.clinicianschoice.com)

BACKED BY OUR TRY-IT GUARANTEE We guarantee you'll love CONTOURS Coarse Anatomy Trimmers, or return them within 30 days for a full refund!