


# Perfect contacts ideal anatomy



FROM  
 triodent

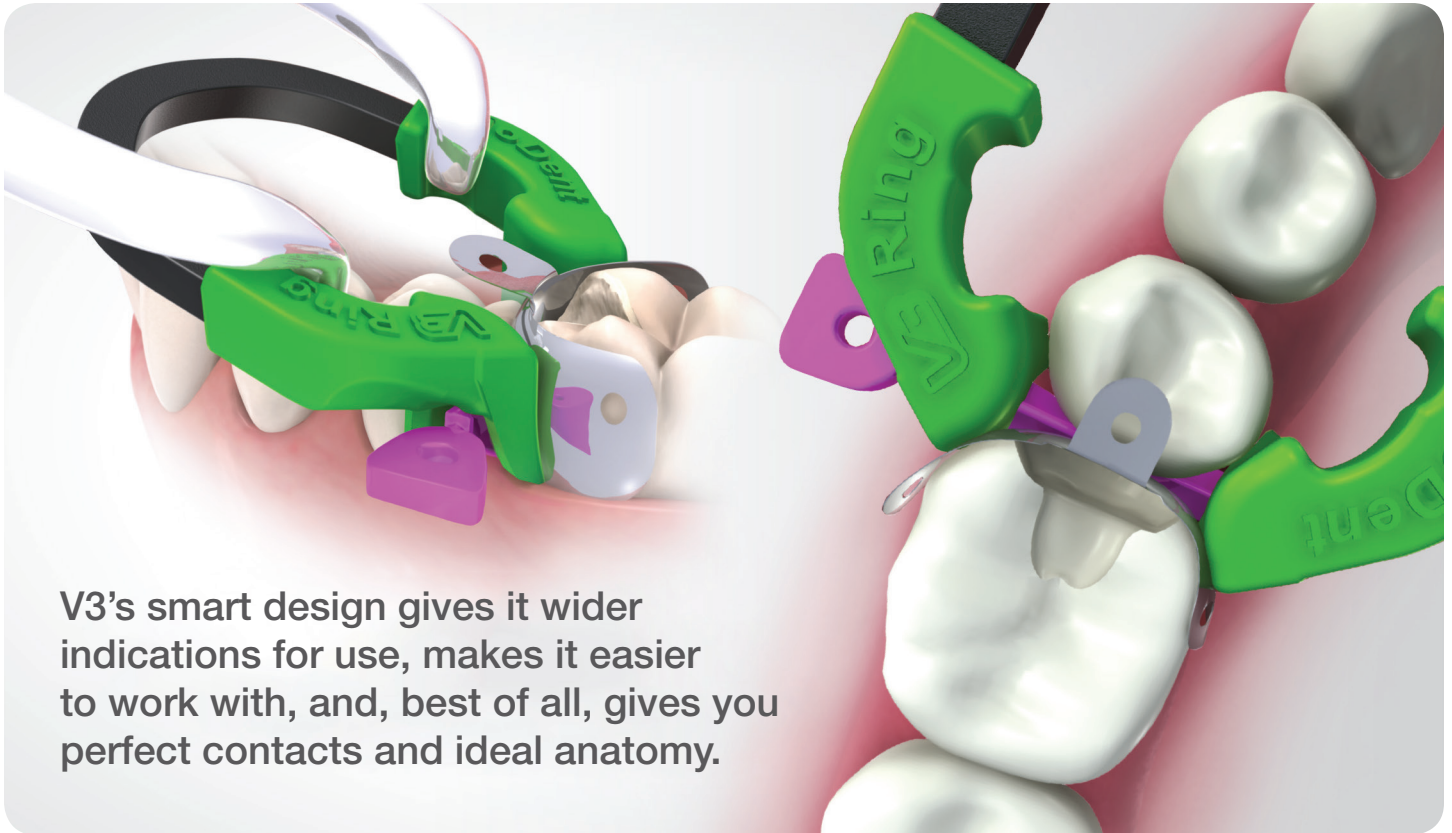

V3's smart design gives it wider indications for use, makes it easier to work with, and, best of all, gives you perfect contacts and ideal anatomy.

## V3—the perfect system for creating perfect Class IIs


- Superelastic nickel-titanium spring for optimal separation force and tight contacts
- Wide indications for use — glass-fibre tines will not collapse into wide cavities, even when a cusp is missing
- Tines are anatomically shaped for superior matrix retention and adaptation—and minimal finishing
- Tines' V-shape accommodates the wedge
- V3 Tab-Matrix with pronounced marginal ridge and vertical profile for correct anatomical form
- Wave-Wedge is stackable, self-locating, and self-seating for optimal seal on matrix
- Spring tilt allows easy stacking of multiple rings
- Ring is very stable on undercuts and in forceps during placement



1.800.265.3444  
clinicalresearchdental.com

Clinical Research Dental

# V3 System Components

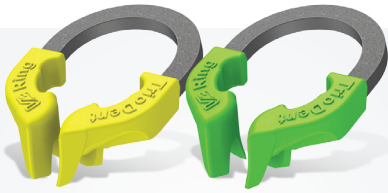

## UNIVERSAL V3 RING - GREEN

- NiTi for ideal spring strength and elasticity
- Long-lasting tines, safe in autoclave
- V-shaped tines accommodate wedge
- Tines shaped for ideal matrix retention and adaptation

## NARROW V3 RING - YELLOW

- Tines set closer together for optimal separation of smaller teeth
- Up to 4kg of separating force, even on premolars



## SUPERCURVE MATRIX

- Matrix hugs the tooth, avoiding interference with ring placement
- Micro-thin, non-stick finish makes removal easy
- Shaped gingival apron leaves no gaps in the gingival-axial corner
- Rolled ridge on occlusal margin provides beautiful marginal ridges



## PIN-TWEEZERS & FORCEPS

- Pin located in tweezer tip grips holes in Trident matrices and wedges for absolute control
- Tweezers' naturally closed position means you will never drop a matrix or wedge
- Handy ball burnisher in tweezer tip
- Forceps' locking function and angled-grip arms hold the V3 Ring securely, both mesially and distally




## WAVE-WEDGE

- Fine wave-shaped wings compress on entry and flare again upon exit for easy placement and ultimate seal
- Hollow underside allows another wedge to be placed from the opposite side of the embrasure



## WEDGEGUARD

- WedgeGuard prevents the bur from nicking the adjacent tooth
- After tooth preparation, the guard is detached while the wedge remains in place
- Saves time and protects gingival papillae

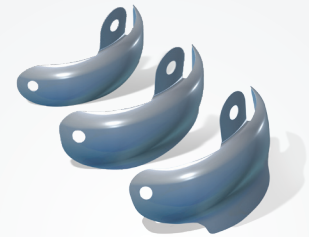

## V3 TAB-MATRIX

- Only V3 Tab-Matrix has a pronounced marginal ridge for natural anatomy
- Pronounced curvature on horizontal plane means more curvature remains as matrix wraps around tooth
- Gingival apron is properly shaped to commonly encountered cavity form, avoiding gaps in the gingival-axial corner

## V3 products are available in the following kit configurations:

### #303301 V3 Matrix Kit

Contains:

- 1 x V3 Rings Universal (green)
- 1 x V3 Ring Narrow (yellow)
- 25 x Wave-Wedges Small (white)
- 25 x Wave-Wedges Medium (pink)
- 25 x Wave-Wedges Large (purple)
- 25 x V3 Tab-Matrices 3.5mm
- 25 x V3 Tab-Matrices 4.5mm
- 25 x V3 Tab-Matrices 5.5mm
- 25 x V3 Tab-Matrices 6.5mm
- 25 x V3 Tab-Matrices 7.5mm
- 2 x Pin-Tweezers
- 1 x Forceps

### V3 Ring 2pk

Available in:

- #303342 - V3 Ring Universal (green) 2pk
- #303343 - V3 Ring Narrow (yellow) 2pk
- #303341 - V3 Ring Mixed (1 green, 1 yellow) 2pk

### Instruments

Available in:

- #302075 - Pin-Tweezers
- #302061 - Forceps

### V3 Tab-Matrix

Available in:

- #303372 - 4.5mm 50pk
- #303375 - 4.5mm 100pk
- #303373 - 5.5mm 50pk
- #303376 - 5.5mm 100pk
- #303374 - 6.5mm 50pk
- #303377 - 6.5mm 100pk

### #303316 V3 SuperCurve Kit

Contains:

- 1 x V3 Rings Universal (green)
- 1 x V3 Ring Narrow (yellow)
- 20 x SuperCurve Matrices 3.5mm (gray)
- 20 x SuperCurve Matrices 4.5mm (green)
- 20 x SuperCurve Matrices 5.5mm (blue)
- 20 x SuperCurve Matrices 6.5mm (purple)
- 20 x SuperCurve Matrices 7.5mm (gray)
- 25 x Wave-Wedges Small (white)
- 25 x Wave-Wedges Medium (pink)
- 25 x Wave-Wedges Large (purple)
- 2 x Pin-Tweezers
- 1 x Forceps

### SuperCurve Matrix

Available in:

- #303402 - Starter Pack
- #303403 - 3.5mm (gray) 90pk
- #303404 - 4.5mm (green) 90pk
- #303405 - 5.5mm (blue) 90pk
- #303406 - 6.5mm (purple) 90pk
- #303417 - 3.5mm (gray) 50pk
- #303418 - 4.5mm (green) 50pk
- #303419 - 5.5mm (blue) 50pk
- #303420 - 6.5mm (purple) 50pk

### Wave-Wedge

Available in:

- #302066 - Small (white) 100pk
- #302065 - Medium (pink) 100pk
- #302064 - Large (purple) 100pk

### WedgeGuard

Available in:

- #303383 - Small (white) 50pk
- #303384 - Medium (pink) 50pk
- #303385 - Large (purple) 50pk