

VALOTM
GRAND

ECLIPSING THE COMPETITION

VALO Grand Cordless Kit

- 1 x VALO Grand LED curing light
- 4 x Batteries (rechargeable)
- 1 x Battery charger
- 1 x Power supply with universal plugs for battery charger
- 1 x Blue light blocking light shield
- 1 x Handpiece bracket holder
- 1 x Sample pack barrier sleeves

- **5972** Black ● **4864** Sapphire
- **4865** Red Rock ● **4866** Midnight

VALO Cordless Kit

- 1 x VALO Cordless LED curing light
- 4 x Batteries (rechargeable)
- 1 x Battery charger
- 1 x Power supply with universal plugs for battery charger
- 1 x Blue light blocking light shield
- 1 x Handpiece bracket holder
- 1 x Sample pack barrier sleeves

- **5941** Black ● **5946** Teal
- **5943** Gold ● **5945** Fuchsia

VALO Grand Corded Kit

- 1 x VALO LED curing light - 7' cord
- 1 x Power supply with universal plugs - 6' cord
- 1 x Handpiece bracket holder
- 1 x Blue light blocking light shield
- 1 x Sample pack barrier sleeves

- **5971** Black

VALO Corded Kit

- 1 x VALO LED curing light - 7' cord
- 1 x Power supply with universal plugs - 6' cord
- 1 x Handpiece bracket holder
- 1 x Blue light blocking light shield
- 1 x Sample pack barrier sleeves

- **5919** Black

BLACK | RED ROCK | SAPPHIRE | MIDNIGHT

ALL THE BENEFITS OF
THE REVOLUTIONARY, AWARD-WINNING
VALO™ CURING LIGHT,
WITH A 50% BIGGER LENS.

AVAILABLE ON ALL VALO CURING LIGHTS



UNPARALLELED POWER

Ultra-high-energy broadband LEDs produce a high-intensity light with an irradiance of 800-2300 mW/cm²



MULTIPLE WAVELENGTHS

Custom LEDs in three wavelengths polymerize all light-cured dental materials



DURABLE UNIBODY CONSTRUCTION

A bar of tempered, high-grade aerospace aluminum is CNC precision milled to create an extremely durable, lightweight handpiece



ACCESSIBILITY

The 11 mm head height allows easy and direct access to all curing sites



COLLIMATED BEAM

Optimally collimated beam delivers consistent, uniform power



MORE UNIFORM CURE

Unique glass lens system forms the collimated blended beam to give an even cure



AWARD-WINNING

The VALO family has won multiple Reality 5 Star Awards, Townie Choice Awards, and Pride Institute "Best of Class" Technology Awards, and has been named one of Dentistry Today's Top 50 Technology Products multiple times

EXCLUSIVE TO THE VALO GRAND LIGHT



LARGER TIP DIAMETER

The VALO Grand light has a 12 mm diameter lens, capable of covering 107 mm²



ADDITIONAL CURING BUTTON


Curing button on the underside facilitates intuitive operation



All VALO LED curing lights use a custom, multiwavelength light-emitting diode (LED) for producing high-intensity light at 385–515 nm, which is capable of polymerizing all light-cured dental materials.

Every member of the VALO curing light family also offers three curing modes to accommodate your preferences.


CURING DEPTH COMPARISON

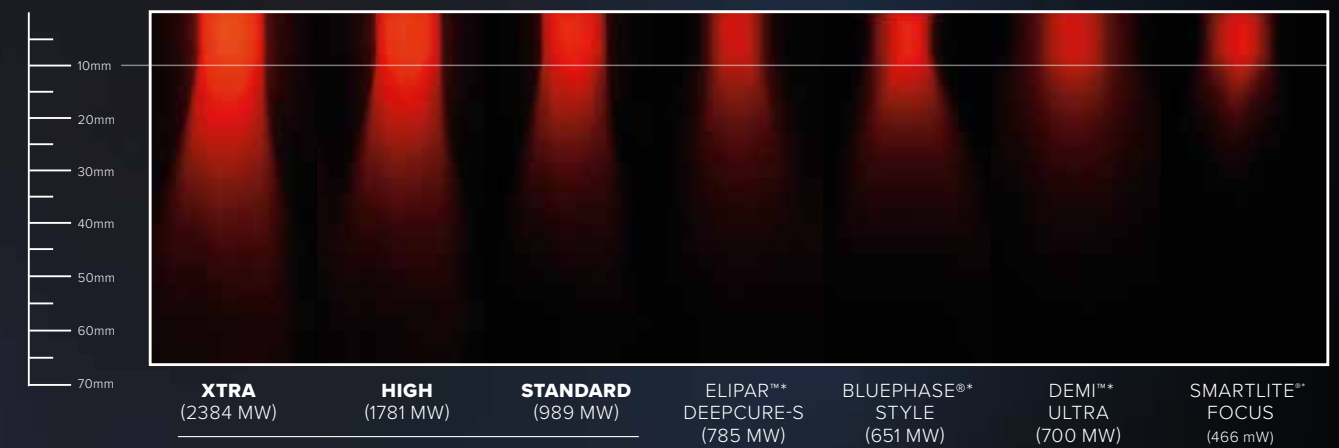

VALO™ GRAND

*Trademark of a company other than Ultradent

AREA COVERAGE COMPARISON

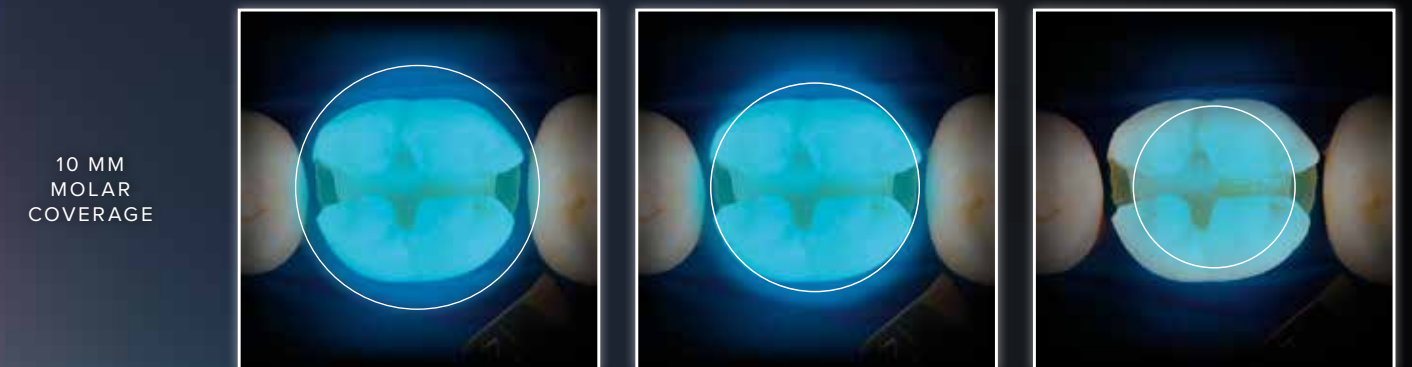

10 MM
MOLAR
COVERAGE

VALO™
GRAND

VALO™

AVERAGE
COMPETITOR

SURFACE
AREA

107
mm²

78
mm²

46
mm²