


Painless Steel® Dental Injection Needles

BETTER COMFORT, NO ANXIETY.



EASE PATIENT ANXIETY & BUILD STRONGER LOYALTY
with these surgically crafted dental injection needles.

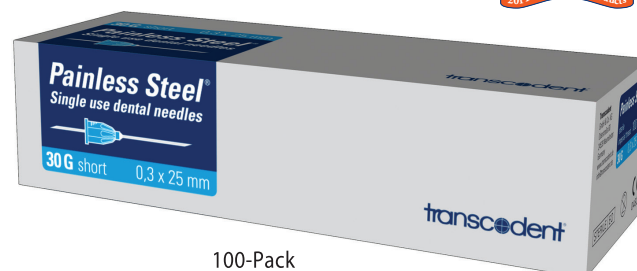

100-Pack



Bevel mark

For an optimal injection position.



Three-edge lancet grinding

The needle penetrates the tissue easily.



Siliconized canula

The needle glides smoothly.



Threaded hub


Ensures secure fit and easy handling



Heat-sealed Needle-container

Using heat sealing eliminates the possibility of contaminating the needle with paper particles of a band

30G

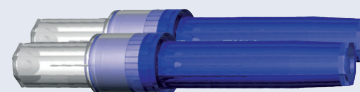


short
x-short

27G



long
short

25G

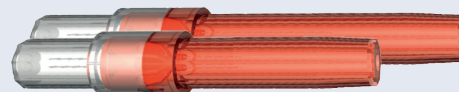

long
short

Size	Diameter	Length		Hub Color
		mm	inch	
30 G short	0.3	25 mm	1"	blue
30 G x-short	0.3	13 mm	1/2"	blue
27 G long	0.4	35 mm	1 3/8"	yellow
27 G short	0.4	25 mm	1"	yellow
25 G long	0.5	35 mm	1 3/8"	red
25 G short	0.5	25 mm	1"	red

Ordering Codes

Blue - Infiltration anesthesia	386861
Blue - Intraligamentary anesthesia	386860
Yellow - Nerve block anesthesia	386863
Yellow - Infiltration anesthesia	386862
Red - Nerve block anesthesia	386865
Red - Infiltration anesthesia	386864



1-800-265-3444 CLINICALRESEARCHDENTAL.COM